Magni Documentation

Version 1.6.0

Christian Schou Oxvig and Patrick Steffen Pedersen in collaboration with Jan Østergaard, Thomas Arildsen, Tobias L. Jensen, and Torben Larsen

Magni Documentation

magni is a Python package developed by Christian Schou Oxvig and Patrick Steffen Pedersen in collaboration with Jan Østergaard, Thomas Arildsen, Tobias L. Jensen, and Torben Larsen. The work was supported by 1) the Danish Council for Independent Research | Technology and Production Sciences - via grant DFF-1335-00278 for the project Enabling Fast Image Acquisition for Atomic Force Microscopy using Compressed Sensing and 2) the Danish e-Infrastructure Cooperation - via a grant for a high performance computing system for the project "High Performance Computing SMP Server for Signal Processing".

This page gives an Introduction to the package, briefly describes How to Read the Documentation, and explains how to actually use The Package.

Introduction

magni [5] is a Python package which provides functionality for increasing the speed of image acquisition using Atomic Force Microscopy (AFM) (see e.g. [1] for an introduction). The image acquisition algorithms of magni are based on the Compressed Sensing (CS) signal acquisition paradigm (see e.g. [2] or [3] for an introduction) and include both sensing and reconstruction. The sensing part of the acquisition generates sensed data from regular images possibly acquired using AFM. This is done by AFM hardware simulation. The reconstruction part of the acquisition reconstructs images from sensed data. This is done by CS reconstruction using well-known CS reconstruction algorithms modified for the purpose. The Python implementation of the above functionality uses the standard library, a number of third-party libraries, and additional utility functionality designed and implemented specifically for magni. The functionality provided by magni can thus be divided into five groups:

- Atomic Force Microscopy (magni.afm): AFM specific functionality including AFM image acquisition, AFM hardware simulation, and AFM data file handling.
- **Compressed Sensing** (magni.cs): General CS functionality including signal reconstruction and phase transition determination.
- **Imaging** (magni.imaging): General imaging functionality including measurement matrix and dictionary construction in addition to visualisation and evaluation.
- Reproducibility (magni.reproducibility): Tools that may aid in increasing the reproducibility of result obtained using magni.
- **Utilities** (magni.utils): General Python utilities including multiprocessing, tracing, and validation.

See Other Resources as well as Notation for further documentation related to the project and the Tests and Examples to draw inspiration from.

References

- [1] B. Bhushan and O. Marti, "Scanning Probe Microscopy Principle of Operation, Instrumentation, and Probes", *in Springer Handbook of Nanotechnology*, pp. 573-617, 2010.
- [2] D.L. Donoho, "Compressed Sensing", *IEEE Transactions on Information Theory*, vol. 52, no. 4, pp. 1289-1306, Apr. 2006.

11/2/16 Page 2 of 209

[3] E.J. Candès, J. Romberg, and T. Tao, "Robust Uncertainty Principles: Exact Signal Reconstruction From Highly Incomplete Frequency Information", *IEEE Transactions on Information Theory*, vol. 52, no.2, pp. 489-509, Feb. 2010.

Footnotes

[5] In Norse mythology, Magni is son of Thor and the god of strength. However, the word MAGNI could as well be an acronym for almost anything including "Making AFM Grind the Normal Impatience".

How to Read the Documentation

The included subpackages, modules, classes and functions are documented through Python docstrings using the same format as the third-party library, numpy, i.e. using the numpydoc standard. A description of any entity can thus be found in the source code of in the docstring of that entity. For readability, the documentation has been compiled using sphinx to produce this HTML page which can be found in the magni folder under '/doc/build/html/index.html'. The entire documentation is also available as a PDF file in the magni folder under '/doc/build/pdf/index.pdf'.

Building the Documentation

The HTML documentation may be built from source using the supplied Makefile in the magni folder under '/doc/'. Make sure the required Dependencies for building the documentation are installed. The build process consists of running three commands:

\$ make sourceclean \$ make docapi \$ make html

Note: In the *make docapi* command it is assumed that the python interpreter is available on the PATH under the name *python*. If the python interpreter is available under another name, the PYTHONINT variable may be set, e.g. "make PYTHONINT=python2 docapi" if the python interpreter is named python2.

Run make clean to remove all builds created by **Sphinx** under '/doc/build'.

The Package

The source code of magni is released under the BSD 2-Clause license, see the License section. To install magni, follow the instructions given under Download and Installation.

magni has been tested with Anaconda (64-bit) on Linux. It may or may not work with other Python distributions and/or operating systems. See also the list of Dependencies for magni.

License

Magni is licensed under the OSI-approved BSD 2-Clause License. See https://opensource.org/licenses/BSD-2-Clause for further information.

11/2/16 Page 3 of 209

Copyright (c) 2014-2016,

Primary developers

Christian Schou Oxvig and Patrick Steffen Pedersen.

Additional developers

Jan Østergaard, Thomas Arildsen, Tobias L. Jensen, and Torben Larsen.

Institution

Aalborg University, Department of Electronic Systems, Signal and Information Processing, Fredrik Bajers Vej 7, DK-9220 Aalborg, Denmark.

All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

- 1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
- 2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT HOLDER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

Download and Installation

All official releases of magni are available for download atdoi:10.5278/VBN/MISC/Magni. The source code is hosted on GitHub at https://github.com/SIP-AAU/Magni.

To use Magni, extract the downloaded archive and include the extracted magni folder in your PYTHONPATH.

Note: The magni package (excluding examples and documentation) is also available on an "as is" basis in source form at PyPi and as a conda package at Anaconda.org.

Dependencies

magni has been designed for use with Python 2 >= 2.7 or Python 3 >= 3.3

Note: The below listed dependency version requirements are the absolute minimum tested versions. Some of these libraries may not work with newer versions of Python unless a newer version of those libraries are used.

11/2/16 Page 4 of 209

Required third party dependencies for magni are:

- Matplotlib (Tested on version >= 1.3)
- Numpy (Tested on version >= 1.8)
- PyTables (Tested on version >= 3.1)
- Scipy (Tested on version >= 0.14)
- Setuptools (Tested on version >= 11.3)

Optional third party dependencies for magni are:

- Bottleneck (Tested on version >=1.0.0) (For speed-up of some algorithms)
- Conda (Tested on version >= 3.7.0) (For automatic metadata capture of a Conda managed Python environment)
- Coverage (Tested on version >= 3.7) (For running the test suite script)
- IPython (Tested on version >= 2.1) or Jupyter [4] (Tested on version >= 1.0) (For running the IPython notebook examples)
- Math Kernel Library (MKL) (Tested on version >= 11.1) (For accelerated vector operations)
- Nose (Tested on version >= 1.3) (For running unittests and doctests)
- PEP8 (Tested on version >= 1.5) (For running style check tests)
- PIL (or Pillow) (Tested on version >= 1.1.7) (For running the IPython notebook examples as tests)
- Pyflakes (Tested on version >= 0.8) (For running style check tests)
- Radon (Tested on version >= 1.2) (For running style check tests)
- scikit-learn (Tested on version >=0.15.1) (For using the scikit-learn logistic regression solver in magni.cs.phase_transition)
- **Sphinx** (Tested on version >= 1.3.1) (For building the documentation from source)
- [4] The IPython project has evolved into Jupyter. When using Jupyter make sure to install the "full" Jupyter metapackage which includes all Jupyter and IPython components.

Note: When using the magni.utils.multiprocessing subpackage, it is generally a good idea to restrict acceleration libraries like MKL or OpenBLAS to use a single thread. If MKL is installed, this is done automatically at runtime in the magni.utils.multiprocessing subpackage. If other libraries than MKL are used, the user has to manually set an appropriate evironmental variable, e.g. OMP_NUM_THREADS.

You may use the *dep_check.py* script found in the Magni folder under '/magni/tests/' to check for missing dependencies for Magni. Simply run the script to print a dependency report, e.g.:

\$ python dep_check.py

Tests

A test suite consisting of unittests, doctests, the IPython notebook examples, and several style checks is included in magni. The tests are organized in python modules found in the Magni folder under '/magni/tests/'. Each module features one or more unittest.TestCase classes containing the tests. Thus, the tests may be invoked using any test runner that supports the unittest.TestCase. E.g. running the wrapper for the doctests usingNose is done by issuing:

11/2/16 Page 5 of 209

\$ nosetests magni/tests/wrap_doctests.py

The entire test suite may be run by executing the convenience scriptrun tests.py:

\$ magni/tests/run_tests.py

Note: This convenience script assumes that magni is available on the PYTHONPATH as explained under Download and Installation.

Bug Reports

Found a bug? Bug report may be submitted using the magniGitHub issue tracker. Please include all relevant details in the bug report, e.g. version of Magni, input/output data, stack traces, etc. If the supplied information does not entail reproducibility of the problem, there is no way we can fix it.

Note: Due to limited funds, we are unfortunately unable make any guarantees, whatsoever, that reported bugs will be fixed.

Other Resources

Papers published in relation to the Enabling Fast Image Acquisition for Atomic Force Microscopy using Compressed Sensing project:

- C. S. Oxvig, T. Arildsen, and T. Larsen, "Structure Assisted Compressed Sensing Reconstruction of Undersampled AFM Images", *Ultramicroscopy*, vol. 172, pp. 1-9, Jan. 2017, doi:10.1016/j.ultramic.2016.09.011.
- C. S. Oxvig, T. Arildsen, and T. Larsen, "Storing Reproducible Results from Computational Experiments using Scientific Python Packages", in *Proceedings of the 15th Python in Science Conference*, pp. 45-50, Austin, Texas, USA, July 11 17, 2016, http://conference.scipy.org/proceedings/scipy2016/christian oxvig.html
- P. S. Pedersen, C. S. Oxvig, J. Østergaard, and T. Larsen, "Validating Function Arguments in Python Signal Processing Applications", in *Proceedings of the 15th Python in Science Conference*, pp. 106-113, Austin, Texas, USA, July 11 - 17, 2016, http://conference.scipy.org/proceedings/scipy2016/patrick_pedersen.html
- P. S. Pedersen, J. Østergaard and T. Larsen, "Modelling reconstruction quality of Lissajous undersampled atomic force microscopy images," in 2016 IEEE 13th International Symposium on Biomedical Imaging (ISBI), pp. 245-248, Prague, Czech Republic, 2016, doi:10.1109/ISBI.2016.7493255.
- T. Arildsen, C. S. Oxvig, P. S. Pedersen, J. Østergaard, and T. Larsen, "Reconstruction Algorithms in Undersampled AFM Imaging", *IEEE Journal of Selected Topics in Signal Processing*, vol. 10, no. 1, pp. 31-46, Feb. 2016, doi:10.1109/JSTSP.2015.2500363.
- P. S. Pedersen, J. Østergaard, and T. Larsen, "Predicting reconstruction quality within compressive sensing for atomic force microscopy," in *2015 IEEE Global Conference on Signal and Information Processing (GlobalSIP)*, pp. 418-422 Orlando, FL, 2015, doi:10.1109/GlobalSIP.2015.7418229.
- C. S. Oxvig, P. S. Pedersen, T. Arildsen, J. Østergaard, and T. Larsen, "Magni: A Python Package for Compressive Sampling and Reconstruction of Atomic Force Microscopy Images", *Journal of Open Research Software*, vol. 2, no. 1, p. e29, Oct.

11/2/16 Page 6 of 209

Page 7 of 209

2014, doi:10.5334/jors.bk.

• T. L. Jensen, T. Arildsen, J. Østergaard, and T. Larsen, "Reconstruction of Undersampled Atomic Force Microscopy Images: Interpolation versus Basis Pursuit", in *International Conference on Signal-Image Technology & Internet-Based Systems* (SITIS), pp. 130-135, Kyoto, Japan, December 2 - 5, 2013,doi:10.1109/SITIS.2013.32.

Notation

To the extent possible, a consistent notation has been used in the documentation and implementation of algorithms that are part of magni. All the details are described in A Note on Notation.

A Note on Notation

As much as possible, a consistent notation is used in the magni package. This implies that variable names are shared between functions that are related. Furthermore a consistent coordinate system is used for the description of related surfaces.

The Compressed Sensing Reconstruction Problem

In the magni.cs subpackage, a consistent naming scheme is used for variables, i.e., vectors and matrices. This section gives an overview of the chosen notation. For the purpose of illustration, consider the Basis Pursuit CS reconstruction problem [1]:

minimize
$$||\alpha||_1$$

subject to $\mathbf{z} = \mathbf{A}\alpha$

Here

 $\mathbf{A} \in \mathbb{C}^{m imes n}$ is the matrix product of a sampling matrix

 $\mathbf{\Phi} \in \mathbb{C}^{m imes p}$ and a dictionary matrix

 $\mathbf{\Psi} \in \mathbb{C}^{p imes n}$. The dictionary coefficients are denoted

 $lpha \in \mathbb{C}^{n imes 1}$ whereas the (noiseless) measurements are denoted by

 $\mathbf{z} \in \mathbb{C}^{m \times 1}$.

Thus, the following relations are used:

$$\mathbf{A} = \mathbf{\Phi} \mathbf{\Psi}$$

$$\mathbf{x}=\mathbf{\Psi}\alpha$$

$$\mathbf{z} = \mathbf{\Phi} \mathbf{x}$$

$$= \mathbf{\Phi} \mathbf{\Psi} \alpha$$

$$= \mathbf{A} \alpha$$

Here the vector

 $\mathbf{x} \in \mathbb{C}^{p \times 1}$ represents the signal of interrest. That is, it is the signal that is assumed to have a sparse representation in the dictionary

 Ψ . The sparsity of the coefficient vector

lpha, that is the size of the support set, is denoted

11/2/16

 $k = |\operatorname{supp}(\alpha)|$. Furthermore, we have the undersampling ratio

 $\delta = m/n$, the sparsity level is

 $\rho = k/m$, and the signal density is

au=k/n.

Oftentimes, one considers noisy measurements, e.g. measurements corrupted by Additive White Gaussian Noise (AWGN)

 $\mathbf{e} \in \mathbb{R}^{m imes 1}$ with

$$\mathbf{y} = \mathbf{z} + \mathbf{e}$$

= $\mathbf{A}\alpha + \mathbf{e}$

Here the (noisy) measurements are denoted

 $\mathbf{y} \in \mathbb{C}^{m \times 1}$.

Note: All vectors in magni are considered to be column vectors, i.e. they explicitly have a second dimension of length 1 when defined as numpy arrays.

Note: Even though the above example involves complex vectors and matrices, the algorithms provided in magni.cs may be restricted to inputs and outputs that are real.

References

[1] S. Chen, D. L. Donoho, and M. A. Saunders, "Atomic Decomposition by Basis Pursuit", *Siam Review*, vol. 43, no. 1, pp. 129-159, Mar. 2001.

Handling Images as Matrices and Vectors

In parts of magni.imaging, an image is considered a matrix

 $\mathbf{M} \in \mathbb{R}^{h imes w}$. That is, the image height is

h whereas the width is


w. In the ${f magni.cs}$ subpackage, the image must be represented as a vector. This is done by stacking the columns of

M to form the vector

x. Thus, the dimension of the image vector representation is

 $n=h\cdot w$. The magni.imaging._util.vec2mat() (available as magni.imaging.vec2mat()) and magni.imaging._util.mat2vec() (available as magni.imaging.mat2vec()) may be use to convert between the matrix and vector notations.

When the matrix representation is used, the following coordinate system is used for its visual representation:

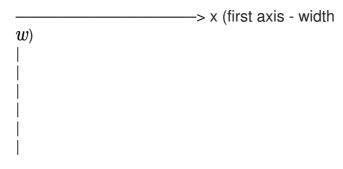

11/2/16 Page 8 of 209

```
v
y (second axis - height h)
```

This way, a position on an AFM sample of size $w \times h$ is specified by a (x,y) coordinate pair.

Examples

The magni package includes a large number of examples showing its capabilities. See the dedicated Examples page for all the details.

Examples

All the examples are available as IPython Notebooks in the magni folder under '/examples/'. For an introduction to getting started with IPython Notebook see the official documentation.

Starting the IPython Notebook

Starting the IPython Notebook basically boils down to running:

ipython notebook

from a shell with the working directory set to the Magni '/examples/' folder. Remember to make sure that magni is available as described in Download and Installation prior to starting the IPython Notebook.

Examples overview

An overview of the available examples is given in the below table:

IPython Notebook Name	Example illustrates	Magni functionality used
afm-io	 Reading data from a mi-file. Handling the resulting buffers and images. 	● magni.afm.io.read_mi_file
cs- phase_transition- config	 Using Magni configuration modules including setting and getting configuration values. 	 magni.cs.phase_transition.config magni.utils.config

11/2/16 Page 9 of 209

IPython .		Ğ
Notebook Name	Example illustrates	Magni functionality used
cs- phase_transition	 Estimating phase transitions using simulations. Plotting phase transitions. Plotting phase transition probability colormaps. 	 magni.cs.phase_transitionutil.determine magni.cs.phase_transition.io magni.cs.phase_transition.plotting
cs-reconstruction	 Reconstruction of compressively sampled 1D signals. 	magni.cs.reconstruction
imaging- dictionaries	 Handling compressed sensing dictionaries using Magni. 	magni.imaging.dictionaries
imaging-domains	 Easy handling of an image in the three domains: image, measurement and sparse (dictionary). 	magni.imaging.domains.MultiDomainImage
imaging- measurements	 Handling sampling/measurement patterns using Magni. Sampling a surface. Sampling an image. Illustrating sampling patterns. 	magni.imaging.measurements
imaging- preprocessing	 Pre-processing an image prior to sampling De-tilting AFM images. 	magni.imaging.preprocessing
magni	 The typical work flow in compressively sampling and reconstructing AFM images using Magni. 	magni.afmmagni.imaging
reporducibility- data	 Obtaining various platform and runtime data for annotations and chases. 	magni.reproducibility.data

11/2/16 Page 10 of 209

IPython Notebook Name	Example illustrates	Magni functionality used
reporducibility-io	 Annotating an HDF5 database to help in improving the reproducibility of the results it contains. 	magni.reproducibility.io
util-matrices	 Using the special Magni Matrix and MatrixCollection classes. 	 magni.utils.matrices.Matrix magni.utils.matrices.MatrixCollection
utils- multiprocessing	 Doing multiprocessing using Magni. 	magni.utils.multiprocessing
utils-plotting	 Using the predefined plotting options in Magni to create clearer and more visually pleasing plots. 	• magni.utils.plotting
utils-validation	 Validation of function parameters. Disabling input validation to reduce computation overhead. 	• magni.utils.validation

API Overview

An overview of the high level magni API is given below:

magni package

Package providing a toolbox for compressed sensing for atomic force microscopy.

Routine listings

afm

Subpackage providing atomic force miscroscopy specific functionality.

CS

Subpackage providing generic compressed sensing functionality.

imaging

Subpackage providing generic imaging functionality.

tests

Subpackage providing unittesting of the other subpackages.

utils

Subpackage providing support functionality for the other subpackages.

11/2/16 Page 11 of 209

Notes

See the README file for additional information.

Subpackages

magni.afm package

Subpackage providing atomic force miscroscopy specific functionality.

The present subpackage includes functionality for handling AFM files and data and functionality for utilizing the other subpackages for such AFM data.

Routine listings

config

Configger providing configuration options for this subpackage.

io

Subpackage providing input/output functionality for .mi files.

reconstruction

Module providing reconstruction and analysis of reconstructed images.

types

Subpackage providing data container classes for .mi files.

Subpackages

magni.afm.io package

Subpackage providing input/output functionality for .mi files.

Currently, only input functionality is provided.

Routine listings

read_mi_file(path)

Read an .mi file given by a path.

Submodules

magni.afm.io. data attachment module

Module providing functionality for attaching data.

The reading of an .mi file is logically separated into four steps of which the functionality provided by this module performs the third step.

Routine listings

attach_data(obj, data)

Attach data to a hierarchical object-structure.

magni.afm.io._data_attachment. attach_data (obj, data)

11/2/16 Page 12 of 209

Attach data to a hierachical object-structure.

In the case of .mi image files, the data should be attached to the buffers of the file. In the case of .mi spectroscopy files, the data should be attached to the chunks of the file.

- **Parameters:** obj (object) The hierarchical object-structure to attach data to.
 - data (numpy.ndarray) The 1-dimensional data.

See also:

magni.afm.io.read_mi_file()

Function using the present function.

magni.afm.io._data_attachment. attach_image_data(obj, data)

Attach data to a hierarchical object-structure representing an image.

The data should be attached to the buffers of the file.

- **Parameters:** obj (object) The hierarchical object-structure to attach data
 - data (numpy.ndarray) The 1-dimensional data.

See also:

attach_data()

Function using the present function.

magni.afm.io. data attachment. **attach spectroscopy data** (obj. data)

Attach data to a hierarchical object-structure representing a spectroscopy.

The data should be attached to the chunks of the file.

- **Parameters:** obj (object) The hierarchical object-structure to attach data to.
 - data (numpy.ndarray) The 1-dimensional data.

See also:

attach_data()

Function using the present function.

magni.afm.io._data_attachment._handle_format_inconsistency (obj, data)

Handle format inconsistency.

The inconsistency is the optional precense of an undocumented thumbnail. If the thumbnail is present, the file header contains a thumbnail parameter specifying the resolution of the thumbnail image.

- **Parameters: obj** (*object*) The hierarchical object-structure to attach data to.
 - data (numpy.ndarray) The 1-dimensional data.

11/2/16 Page 13 of 209

Returns:

- **obj** (*object*) The hierarchical object-structure to attach data to.
- data (numpy.ndarray) The 1-dimensional data excluding any thumbnail data.
- **thumbnail** (*list*) A list of a thumbnail object-structure if any.

See also:

_attach_image_data()

Function using the present function.

Notes

The thumbnail is stored as 8-bit unsigned integers with a red, a green, and a blue channel. If the thumbnail is present, it is represented by three buffer object-structures; one for each channel.

magni.afm.io. object building module

Module providing functionality for building a hierarchical object-structure.

The reading of an .mi file is logically separated into four steps of which the functionality provided by this module performs the second step.

Routine listings

build_object(params)

Build a hierarchical object-structure from header parameters.

magni.afm.io. object building. **build object**(params)

Build a hierarchical object-structure from header parameters.

The hierarchical object-structure mimics that of magni.afm.types.

Parameters: params (*list or tuple*) – The header

parameters.

Returns: obj (*object*) – A hierarchical object-structure.

See also:

magni.afm.io.read_mi_file()

Function using the present function.

Notes

This function splits the header parameters into file-related parameters and buffer-related parameters.

magni.afm.io._object_building._build_buffer(file_type, index, params)

Build a buffer-like object-structure from buffer parameters.

For spectroscopy buffers, this function converts parameters which contain "sub-parameters" to objects with attributes.

11/2/16 Page 14 of 209

Parameters: • file_type (str) – The .mi file type.

• **index** (*int*) – The index of the buffer.

• params (list or tuple) - The buffer

parameters.

Returns: obj (*object*) – The buffer-like object-structure.

See also:

build_object()

Function using the present function.

magni.afm.io._object_building._expand_buffers(buffers)

Expand the buffer-like object-structures.

Grid parameters specify implicit point parameters. Furthermore, some chunk parameters may be implicitly specified. These implicit parameters should be made explicit.

Parameters: buffers (*list or tuple*) — The buffer-like object-structures.

See also:

build_object()

Function using the present function.

_generate_grid_points()

Make implicit point parameters explicit.

_generate_implicit_chunks()

Make implicit chunk parameters explicit.

Notes

This function relies on the two functions, _generate_grid_points and _generate_implicit_chunks to make the implicit parameters explicit.

magni.afm.io._object_building. **generate_grid_points**(attrs, index)

Make implicit grid point parameters explicit.

Parameters: • attrs (list or tuple) - The attributes of the

• **index** (*int*) – The index of the first grid point.

Returns: points (*tuple*) – The grid points.

See also:

expand buffers()

Function using the present function.

magni.afm.io._object_building._generate_implicit_chunks(npoints, explicit_chunks)

11/2/16 Page 15 of 209

Make implicit chunk parameters explicit.

Parameters: • npoints (int) – The number of points.

• **explicit_chunks** (*list*) – The explicit chunks.

Returns:

chunks (list) - The chunks including both explicit and implicit chunks.

See also:

_expand_buffer()

Function using the present function.

magni.afm.io._object_building._handle_format_inconsistency(obj)

Handle format inconsistency.

The inconsistency is the usage of string values for the 'trace' header parameter which is specified to have a boolean value.

Parameters: obj (*object*) – The hierarchical object-structure.

Returns:

obj (*object*) – The updated hierarchical object-structure.

See also:

build_object()

Function using the present function.

magni.afm.io. stream conversion module

Module providing functionality for converting a file stream.

The reading of an .mi file is logically separated into four steps of which the functionality provided by this module performs the first step.

Routine listings

convert stream(stream)

Convert a file stream to flat, basic variables.

magni.afm.io._stream_conversion.convert_stream(stream)

Convert a file stream to flat, basic variables.

The flat, basic variables are separated into header parameters and the data. The header parameters are name, value pairs whereas the data is a 1-dimensional numeric array.

Parameters: stream (*str*) – The file stream.

Returns:

- params (tuple) The header parameters.
- data (numpy.ndarray) The data.

See also:

magni.afm.io.read_mi_file()

11/2/16 Page 16 of 209 Function using the present function.

Notes

This function splits the file stream into header and data, and splits the header into parameter name, value pairs.

magni.afm.io._stream_conversion._convert_image_data(buffer_, data_type)

Convert the file stream image data to a 1-dimensional numeric array.

Parameters: • **buffer** (*str*) – The file stream image data.

• data_type (*str*) – The data type of the file stream image data.

Returns:

data (numpy.ndarray) - The data as a 1-dimensional numeric

array.

See also:

convert_stream()

Function using the present function.

magni.afm.io. stream conversion. **convert parameter value**(string)

Convert a file stream parameter value to a basic python value.

The converted value is either a boolean, string, floating point, or integer value or a list containing a mix of these.

Parameters: string (*str*) – The file stream parameter value.

Returns:

value (None) - The converted, basic python

value.

See also:

convert_stream()

Function using the present function.

magni.afm.io._stream_conversion._convert_spectroscopy_data(buffer_, data_type)

Convert the file stream spectroscopy data to a 1-dimensional numeric array.

- **Parameters: buffer** (*str*) The file stream spectroscopy data.
 - data_type (str) The data type of the file stream spectroscopy data.

Returns:

data (*numpy.ndarray*) – The data as a 1-dimensional numeric array.

See also:

convert_stream()

Function using the present function.

magni.afm.io. util module

11/2/16 Page 17 of 209 Module providing the public functionality of the magni.afm.io subpackage.

```
magni.afm.io._util. read_mi_file(path)
```

Read an .mi file given by a path.

The supported .mi file types are Image and Spectroscopy.

```
Parameters: path (str) – The path of the .mi file to read. file_ (magni.afm.types.File) – The read .mi file.
```

```
See also:

magni.afm.io._stream_conversion()

First step of the reading.

magni.afm.io._object_building()

Second step of the reading.

magni.afm.io._data_attachment()

Third step of the reading.
```

Notes

Depending on the .mi file type, the returned value will be an instance of a subclass of magni.afm.types.File: magni.afm.types.image.Image Or magni.afm.types.spectroscopy.Spectroscopy.

The reading of an .mi file is logically separated into four steps:

- 1. Converting the file stream to flat, basic variables.
- 2. Converting the header parameters to a hierarchical object-structure mimicking that of magni.afm.types.
- 3. Attaching the actual data to the hierarchical object-structure.
- 4. Instantiating the classes in magni.afm.types from the hierarchical object structure.

The functionality needed for each step is grouped in separate modules with the functionality needed for the fourth step being grouped in the this module.

Examples

An example of how to use this function to read the example .mi file provided with the package:

```
>>> import os, magni
>>> from magni.afm.io import read_mi_file
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... mi_file = read_mi_file(path)
```

magni.afm.io. util. instantiate spectroscopy (obi)

Instantiate a Spectroscopy object from a hierarchical object-structure.

Parameters: obj (*object*) – The hierarchical object-structure.

11/2/16 Page 18 of 209

Returns:

spectroscopy (*magni.afm.types.spectroscopy*.*Spectroscopy*) – The instantiated Spectroscopy object.

See also:

read_mi_file()

Function using the present function.

Notes

This function instantiates the sweep buffer and the spectroscopy file.

magni.afm.io._util._instantiate_spectroscopy_buffer(data_buffer, sweep_buffer)

Instantiate a spectroscopy Buffer object from a hierarchical object-structure.

Parameters: obj (*object*) – The hierarchical object-structure.

Returns: buffer_ (*magni.afm.types.spectroscopy.Buffer*) – The instantiated

spectroscopy Buffer.

See also:

_instantiate_spectroscopy()

Function using the present function.

magni.afm.types package

Subpackage providing data container classes for .mi files.

Routine listings

image

Module providing data container classes for .mi image files.

spectroscopy

Module providing data container classes for .mi spectroscopy files.

BaseClass(object)

Base class of every magni.afm.types data class.

File(BaseClass)

Base class of the magni.afm.types file classes.

FileCollection(BaseClass)

Data class for collections of File instances with identical settings.

Submodules

magni.afm.types. util module

Module providing the common functionality of the subpackage.

class magni.afm.types._util. BaseClass(attrs)

Bases: object

11/2/16 Page 19 of 209

Base class of every magni.afm.types data class.

The present class validates the attributes passed to its constructor against the allowed attributes of the class of the given instance and exposes these attributes through a read-only dictionary property, attrs. Furthermore, the present class exposes the allowed attributes of the class of the given instance through a read-only dictionary static property, params.

```
Parameters: attrs (dict) – The desired attributes of the instance.
```

attrs

magni.utils.types.ReadOnlyDict – The attributes of the instance.

params

magni.utils.types.ReadOnlyDict - The allowed attributes of the instance.

Examples

An example could be the subclass, 'Person' which allows only the string attribute, name:

```
>>> from magni.afm.types._util import BaseClass
>>> class Person(BaseClass):
... def __init__(self, attrs):
... BaseClass.__init__(self, attrs)
... _params = {'name': str}
```

This class can then be initiated with a name:

```
>>> person = Person({'name': 'Murphy'})
>>> print('{\lr\frac{\rr}}'.format(person.attrs['name']))
'Murphy'
```

Any other attributes, than 'name', passed to the class are ignored:

```
>>> person = Person({'name': 'Murphy', 'age': 42})
>>> for name in person.attrs.keys():
... print('{\ll r\rrangle'}'.format(name))
'name'
```

attrs

```
params = ReadOnlyDict()
```

class magni.afm.types. util. File (attrs, buffers)

Bases: magni.afm.types._util.BaseClass

Base class of the magni.afm.types file classes.

The present class specifies the allowed file-level attributes and exposes the read-only property, buffers which all .mi files have. Furthermore, the present class provides a method for accessing buffers by their 'bufferLabel' attribute.

```
attrs (list or tuple) – The desired attributes of the file.
buffers (list or tuple) – The buffers of the file.
```

11/2/16 Page 20 of 209

buffers

tuple – The buffers of the file.

See also:

BaseClass

Superclass of the present class.

Examples

A subclass of the present class is implicitly instantiated when loading, for example, the .mi file provided with the package:

```
>>> import os, magni
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... file_ = magni.afm.io.read_mi_file(path)
```

This file has a number of buffers which each has the 'bufferLabel' attribute:

```
>>> if os.path.isfile(path):
... for buffer_ in file_.buffers[::2][:3]:
... label = buffer_.attrs['bufferLabel']
... print("Buffer with 'bufferLabel': {!r}".format(label))
... else:
... for label in ('Topography', 'Deflection', 'Friction'):
... print("Buffer with 'bufferLabel': {!r}".format(label))
Buffer with 'bufferLabel': 'Topography'
Buffer with 'bufferLabel': 'Deflection'
Buffer with 'bufferLabel': 'Friction'
```

If only, for example, buffers with 'bufferLabel' equal to 'Topography' are desired, the method, get_buffer can be called:

```
>>> if os.path.isfile(path):
... buffers = len(file_.get_buffer('Topography'))
... print("Buffers with 'bufferLabel' == 'Topography': {}"
... .format(buffers))
... else:
... print("Buffers with 'bufferLabel' == 'Topography': 2")
Buffers with 'bufferLabel' == 'Topography': 2
```

buffers

get_buffer(label)

Get the buffers that have the specified buffer label.

Parameters: label (*str*) – The desired buffer label of the buffer. **Beturns:** buffers (*list*) – The buffers that have the desired buffer label.

class magni.afm.types._util. FileCollection (files, paths)

Bases: magni.afm.types._util.BaseClass

Data class for collections of File instances with identical settings.

11/2/16 Page 21 of 209

The settings are the following attributes: 'fileType', 'mode', 'xPixels', 'yPixels', 'xOffset', 'yOffset', 'xLength', 'yLength', 'scanSpeed', 'acMac', 'acACMode', 'plotType'.

The present class exposes the files of the collection through a read-only tuple property, files and the paths of these files through a read-only tuple property paths.

- **Parameters:** files (*list or tuple*) The files of the collection.
 - paths (list or tuple) The paths of the files of the collection.

files

tuple - The files of the collection.

paths

tuple – The paths of the files of the collection.

See also:

BaseClass

Superclass of the present class.

Examples

No example .mi file collection is distributed with magni.

files

paths

magni.afm.types.image module

Module providing data container classes for .mi image files.

The classes of this module can be used either directly or indirectly through themagni.afm.io subpackage by loading an .mi image file.

Routine listings

Buffer(magni.afm.types.BaseClass)

Data class of the .mi image file buffers.

Image(magni.afm.types.File)

Data class of the .mi image files.

See also:

magni.afm.io

.mi file loading.

class magni.afm.types.image. **Buffer**(attrs, data)

Bases: magni.afm.types._util.BaseClass

Data class of the .mi image file buffers.

11/2/16 Page 22 of 209

```
attrs (dict) – The attributes of the buffer.
data (numpy.ndarray) – The 2D data of the buffer.
```

apply_clipping

bool – A flag indicating if clipping should be applied to the data.

data

numpy.ndarray - The 2D data of the buffer.

See also:

magni.utils.types.BaseClass

Superclass of the present class.

Examples

A subclass of the present class is implicitly instantiated when loading, for example, the .mi file provided with the package:

```
>>> import os, magni
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... image = magni.afm.io.read_mi_file(path)
... buffer_ = image.buffers[0]
```

This buffer can have a number of attributes including 'bufferLabel':

```
>>> if os.path.isfile(path):
... print('{!r}'.format(buffer_.attrs['bufferLabel']))
... else:
... print("Topography")
'Topography'
```

The primary purpose of this class is, however, to contain the 2D data of a buffer:

```
>>> if os.path.isfile(path):
... print('Buffer data shape: {\llower lr'}'.format(
... tuple(int(s) for s in buffer_.data.shape)))
... else:
... print('Buffer data shape: (256, 256)')
Buffer data shape: (256, 256)
```

apply clipping

Get the apply clipping property of the buffer.

The clipping is specified by the 'DisplayOffset', 'DisplayRange', and 'filter' attributes of the buffer.

Returns: apply_clipping (bool) – A flag indicating if clipping should be applied to the data.

data

Get the data property of the buffer.

11/2/16 Page 23 of 209

If apply_clipping is False, the data is returned as-is. Otherwise, the filtering specified by the 'filter' attribute of the buffer and the clipping specified by the 'DisplayOffset' and 'DisplayRange' attributes of the buffer are applied to the data before it is returned.

```
Returns: data (numpy.ndarray) – The 2D data of the buffer.
```

class magni.afm.types.image. Image (attrs, buffers)

```
Bases: magni.afm.types._util.File
```

Data class of the .mi image files.

```
attrs (dict) – The attributes of the image.
buffers (list or tuple) – The buffers of the image.
```

See also:

```
magni.utils.types.File
```

Superclass of the present class.

Examples

The present class is implicitly instantiated when loading, for example, the .mi file provided with the package:

```
>>> import os, magni
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... image = magni.afm.io.read_mi_file(path)
```

This file has a number of buffers which each has the 'bufferLabel' attribute:

```
>>> if os.path.isfile(path):
... for buffer_ in image.buffers[::2][:3]:
... label = buffer_.attrs['bufferLabel']
... print("Buffer with 'bufferLabel': {!r}".format(label))
... else:
... for label in ('Topography', 'Deflection', 'Friction'):
... print("Buffer with 'bufferLabel': {!r}".format(label))
Buffer with 'bufferLabel': 'Topography'
Buffer with 'bufferLabel': 'Deflection'
Buffer with 'bufferLabel': 'Friction'
```

magni.afm.types.spectroscopy module

Module providing data container classes for .mi spectroscopy files.

The classes of this module can be used either directly or indirectly through their module by loading an .mi spectroscopy file.

Routine listings

Buffer(magni.afmtypes.BaseClass)

11/2/16 Page 24 of 209

Data class for .mi spectroscopy buffer.

Chunk(magni.afmtypes.BaseClass)

Data class for .mi spectroscopy chunk.

Grid(magni.afmtypes.BaseClass)

Data class for .mi spectroscopy grid.

Point(magni.afmtypes.BaseClass)

Data class for .mi spectroscopy point.

Spectroscopy(magni.afm.types.File)

Data class for .mi spectroscopy.

See also:

magni.afm.io

.mi file loading.

class magni.afm.types.spectroscopy. **Buffer** (attrs, data)

Bases: magni.afm.types._util.BaseClass

Data class of the .mi spectroscopy buffers.

- **Parameters:** attrs (*dict*) The attributes of the buffer.
 - data (list or tuple) The grids, points, or chunks of the buffer.

data

numpy.ndarray – The grids, points, or chunks of the buffer.

See also:

magni.utils.types.BaseClass

Superclass of the present class.

Examples

No example .mi spectroscopy file is distributed with magni.

data

class magni.afm.types.spectroscopy. Chunk (attrs, data)

Bases: magni.afm.types._util.BaseClass

Data class of the .mi spectroscopy chunks.

- **Parameters:** attrs (*dict*) The attributes of the chunk.
 - data (numpy.ndarray) The data of the chunk.

data

numpy.ndarray – The data of the chunk.

11/2/16 Page 25 of 209

See also:

magni.utils.types.BaseClass

Superclass of the present class.

Examples

No example .mi spectroscopy file is distributed with magni.

data

sweep

Get the sweep property of the chunk.

The sweep property is the series of the entity which was swept.

Returns: sweep (numpy.ndarray) – The sweep property.

Notes

To reduce the memory footprint of chunks, the series does not exist explicitly, until it is requested.

time

Get the time property of the chunk.

Returns: time (numpy.ndarray) – The time property.

Notes

To reduce the memory footprint of chunks, the series does not exist explicitly, until it is requested.

class magni.afm.types.spectroscopy. **Grid** (attrs, points)

Bases: magni.afm.types._util.BaseClass

Data class of the .mi spectroscopy grids.

attrs (dict) – The attributes of the grid.
points (list or tuple) – The points of the grid.

points

tuple – The points of the grid.

See also:

magni.utils.types.BaseClass

Superclass of the present class.

11/2/16 Page 26 of 209

Notes

The points are input and output as a 2D tuple of Point instances.

Examples

No example .mi spectroscopy file is distributed with magni.

points

class magni.afm.types.spectroscopy. **Point**(attrs, chunks=())

Bases: magni.afm.types._util.BaseClass

Data class of the .mi spectroscopy points.

- **Parameters:** attrs (*dict*) The attributes of the point.
 - **chunks** (*list or tuple, optional*) The chunks of the point. (the default is (), which implies no chunks)

chunks

tuple – The chunks of the point.

See also:

magni.utils.types.BaseClass

Superclass of the present class.

Examples

No example .mi spectroscopy file is distributed with magni.

chunks

class magni.afm.types.spectroscopy. **Spectroscopy** (*attrs*, *buffers*)

Bases: magni.afm.types._util.File

Data class of the .mi spectroscopy files.

- **Parameters:** attrs (*dict*) The attributes of the image.
 - buffers (list or tuple) The buffers of the image.

See also:

magni.utils.types.File

Superclass of the present class.

Examples

11/2/16 Page 27 of 209 No example .mi spectroscopy file is distributed with magni.

Submodules

magni.afm._config module

Module providing configuration options for the magni.afm subpackage.

See also:

magni.utils.config.Configger

The Configger class used

Notes

This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

algorithm

: {'amp', 'gamp', 'it', 'iht', 'sl0'}

The compressed sensing reconstruction algorithm subpackage to use (the default is 'it').

magni.afm.reconstruction module

Module providing AFM image reconstruction and analysis of reconstructed images.

Routine listings

analyse(x, Phi, Psi)

Sample an image, reconstruct it, and analyse the reconstructed image.

reconstruct(y, Phi, Psi)

Reconstruct an image from compressively sensed measurements.

```
magni.afm.reconstruction. analyse (x, Phi, Psi, **kwargs)
```

Sample an image, reconstruct it, and analyse the reconstructed image.

- **Parameters:** x (numpy.ndarray) The original image vector.
 - (magni.utils.matrices.Matrix Phi *numpy.matrix*) The measurement matrix.
 - **Psi** (*magni.utils.matrices.Matrix or numpy.matrix*) The dictionary.

Returns:

x (*numpy.ndarray*) – The reconstructed image vector.

See also:

magni.afm.config()

Configuration options.

magni.imaging.evaluation()

Image reconstruction quality evaluation.

11/2/16 Page 28 of 209

Notes

Additional **kwargs are passed on to the reconstruction algorithm.

Examples

Prior to the actual example, data is loaded and a measurement matrix and a dictionary are defined. First, the example MI file provided with the package is loaded:

```
>>> import os, numpy as np, magni
>>> from magni.afm.reconstruction import analyse
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... mi_file = magni.afm.io.read_mi_file(path)
... mi_buffer = mi_file.get_buffer('Topography')[0]
... mi_data = mi_buffer.data
... x = magni.imaging.mat2vec(mi_data)
```

Next, a measurement matrix is defined. This matrix is equal to the matrix generated by running np.eye(len(x))[::2, :] but for speed, the matrix is instead defined with fast operations:

Next, a dictionary is defined. This dictionary is the DCT basis likewise defined with fast operations:

```
>>> if os.path.isfile(path):
... Psi = magni.imaging.dictionaries.get_DCT(mi_data.shape)
```

Finally, the actual example:

```
>>> if os.path.isfile(path):
... print('MSE: {:.2f}, PSNR: {:.2f}'.format(*analyse(x, Phi, Psi)))
... else:
... print('MSE: 0.24, PSNR: 6.22')
MSE: 0.24, PSNR: 6.22
```

magni.afm.reconstruction. **reconstruct**(y, Phi, Psi, **kwargs)

Reconstruct an image from compressively sensed measurements.

Parameters: • y (numpy.ndarray) – The measurement vector.

Phi (magni.utils.matrices.Matrix or numpy.matrix) – The measurement matrix.

• **Psi** (*magni.utils.matrices.Matrix or numpy.matrix*) – The dictionary.

Returns: x (*numpy.ndarray*) – The reconstructed image vector.

See also:

11/2/16 Page 29 of 209

```
magni.afm.config()
Configuration options.
magni.cs.reconstruction()
Compressed sensing reconstruction algorithms.
```

Notes

Additional **kwargs are passed on to the reconstruction algorithm.

Examples

Prior to the actual example, data is loaded and a measurement matrix and a dictionary are defined. First, the example MI file provided with the package is loaded:

```
>>> import os, numpy as np, magni
>>> from magni.afm.reconstruction import reconstruct
>>> path = magni.utils.split_path(magni.__path__[0])[0]
>>> path = path + 'examples' + os.sep + 'example.mi'
>>> if os.path.isfile(path):
... mi_file = magni.afm.io.read_mi_file(path)
... mi_buffer = mi_file.get_buffer('Topography')[0]
... mi_data = mi_buffer.data
... x = magni.imaging.mat2vec(mi_data)
```

Next, a measurement matrix is defined. This matrix is equal to the matrix generated by running np.eye(len(x))[::2, :] but for speed, the matrix is instead defined with fast operations:

Next, a dictionary is defined. This dictionary is the DCT basis likewise defined with fast operations:

```
>>> if os.path.isfile(path):
... Psi = magni.imaging.dictionaries.get_DCT(mi_data.shape)
```

Finally, the actual example:

```
>>> if os.path.isfile(path):
... y = Phi.dot(x)
... print('Maximum absolute pixel error: {:.3f}'
... .format(np.abs(reconstruct(y, Phi, Psi) - x).max()))
... else:
... print('Maximum absolute pixel error: 0.960')
Maximum absolute pixel error: 0.960
```

magni.cs package

11/2/16 Page 30 of 209

Subpackage providing generic compressed sensing functionality.

Routine listings

indicators

Module prodiving performance indicator determination functionality.

phase transition

Subpackage providing phase transition determination functionality.

reconstruction

Subpackage providing implementations of generic reconstruction algorithms.

Subpackages

magni.cs.phase transition package

Subpackage providing phase transition determination functionality.

Routine listings

config

Configger providing configuration options for this subpackage.

io

Module providing input/output functionality for stored phase transitions.

plotting

Module providing plotting for this subpackage.

determine(algorithm, path, label='default', overwrite=False)

Determine the phase transition of a reconstruction algorithm.

Notes

See <u>util</u> for documentation of determine.

The phase transition of a reconstruction algorithm describes the reconstruction capabilities of that reconstruction algorithm. For a description of the concept of phase transition, see [1].

References

[1] C. S. Oxvig, P. S. Pedersen, T. Arildsen, and T. Larsen, "Surpassing the Theoretical 1-norm Phase Transition in Compressive Sensing by Tuning the Smoothed IO Algorithm", in IEEE International Conference on Acoustics, Speech and Signal Processing (ICASSP), Vancouver, Canada, May 26-31, 2013, pp. 6019-6023.

Submodules

magni.cs.phase transition. analysis module

Module providing functionality for analysing the simulation results.

Routine listings

run(path, label)

Determine the phase transition from the simulation results.

11/2/16 Page 31 of 209

See also:

magni.cs.phase transition.config Configuration options.

Notes

For a description of the concept of phase transition, see[1].

References

[1] C. S. Oxvig, P. S. Pedersen, T. Arildsen, and T. Larsen, "Surpassing the Theoretical 1-norm Phase Transition in Compressive Sensing by Tuning the Smoothed IO Algorithm", in IEEE International Conference on Acoustics, Speech and Signal Processing (ICASSP), Vancouver, Canada, May 26-31, 2013, pp. 6019-6023.

magni.cs.phase transition. analysis. **run** (path, label)

Determine the phase transition from the simulation results.

The simulation results should be present in the HDF5 database specified bypath in the pytables group specified by label in an array named 'dist'. The determined phase transition (50% curve) is stored in the same HDF5 database, in the same HDF group in an array named 'phase_transition'. Additionally, the 10%, 25%, 75%, and 90% percentiles are stored in an array named 'phase_transition_percentiles'.

- **Parameters:** path (*str*) The path of the HDF5 database.
 - label (str) The path of the pytables group in the HDF5 database.

See also:

estimate PT()

The actual phase transition estimation.

Notes

A simulation is considered successful if the simulation result is less than a normalised mean squared error tolerance computed as 10^(-SNR/10) with SNR configured in the configuration module.

magni.cs.phase_transition._analysis._built_in_logit_solver(dist, percentiles)

Fit a logistic regression model using the built-in solver.

- **Parameters: dist** (*ndarray*) The simulated signal "distances" in the phase space.
 - **percentiles** (*list or tuple*) The percentiles to estimate.

Returns:

rho (*ndarray*) – The "len(percentiles)"-by-"len(delta)" array of estimated phase transition rho vectors. The phase transition rho vectors are (in order): 50% (the phase transition esitmate), smaller to larger percentiles.

11/2/16 Page 32 of 209 magni.cs.phase transition. analysis. **estimate PT**(*rho*, *success*, *percentiles*)

Estimate the phase transition location for a given delta.

The phase transition location is estimated using logistic regression. The algorithm used for this is Newton's method.

- **Parameters:** rho (ndarray) The rho values.
 - **success** (*ndarray*) The success indicators.
 - percentiles (list or tuple) The percentiles to estimate.

Returns:

rho (*float*) – The estimated phase transition location.

Notes

The function includes a number of non-standard ways of handling numerical and convergence related issues. This will be changed in a future version of the code.

magni.cs.phase_transition._analysis._**sklearn_logit_solver**(*dist*, *percentiles*)

Fit a logistic regression model using the solver from scikit-learn.

- **Parameters: dist** (*ndarray*) The simulated signal "distances" in the phase space.
 - **percentiles** (*list or tuple*) The percentiles to estimate.

Returns:

rho (*ndarray*) – The "len(percentiles)"-by-"len(delta)" array of estimated phase transition rho vectors. The phase transition rho vectors are (in order): 50% (the phase transition esitmate), smaller to larger percentiles.

magni.cs.phase transition. backup module

Module providing backup capabilities for the monte carlo simulations.

The backup stores the simulation results and the simulation timings pointwise for the points in the delta-rho simulation grid. The set function targets a specific point while the get function targets the entire grid in order to keep the overhead low.

Routine listings

create(path)

Create the HDF5 backup database with the required arrays.

get(path)

Return which of the results have been stored.

set(path, ij tuple, stat time, stat dist)

Store the simulation data of a specified point.

See also:

magni.cs.phase_transition.config Configuration options.

Notes

11/2/16 Page 33 of 209

In practice, the backup database includes an additional array for tracking for which points data has been stored. By first storing the data and then modifying this array, the data is guaranteed to have been stored, when the array is modified.

magni.cs.phase_transition._backup. **create**(path)

Create the HDF5 backup database with the required arrays.

The required arrays are an array for the simulation results, an array for the simulation timings, and an array for tracking the status.

Parameters: path (*str*) – The path of the HDF5 backup

database.

See also:

magni.cs.phase_transition.config()

Configuration options.

magni.cs.phase_transition._backup. **get**(*path*)

Return which of the results have been stored.

The returned value is a copy of the boolean status array indicating which of the results have been stored.

Parameters: path (*str*) – The path of the HDF5 backup

database.

Returns: **status** (*ndarray*) – The boolean status array.

magni.cs.phase_transition._backup. **set** (path, ij_tuple, stat_time, stat_dist, stat_mse, stat_norm)

Store the simulation data of a specified point.

- **Parameters:** path (*str*) The path of the HDF5 backup database.
 - ij tuple (tuple) A tuple (i, j) containing the parameters i, j as listed below.
 - i (int) The delta-index of the point in the delta-rho grid.
 - i (int) The rho-index of the point in the delta-rho grid.
 - **stat dist** (*ndarray*) The simulation results of the specified point.
 - **stat_time** (*ndarray*) The simulation timings of the specified point.
 - stat mse (ndarray) The simulation mean squared error of the specified point.
 - **stat_norm** (*ndarray*) The simulation true vector 2-norm.

magni.cs.phase transition. config module

Module providing configuration options for the phase_transition subpackage.

See also:

magni.utils.config.Configger

The Configger class used

Notes

11/2/16 Page 34 of 209 This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

algorithm_kwargs

: dict

The keyword arguments passed on to the reconstruction algorithm.

coefficients

: {'rademacher', 'gaussian', 'laplace', 'bernoulli'}

The distribution which the non-zero coefficients in the coefficient vector are drawn from.

custom noise factory

: callable

The callable that generates a custom m-by-1 vector of noise from the four arguments "m", "n", "k", "noise power".

custom system matrix factory

: callable

The callable that generates a custom m-by-n system matrix from the two arguments "m", "n".

delta

: list or tuple

The delta values of the delta-rho grid whose points are used for the monte carlo simulations (the default is [0., 1.]).

logit solver

: {'built-in', 'sklearn'}

The solver to use in the logistic regression fit of the phase transition.

maxpoints

: int

The maximum number of phase space grid points to be handled by a process before it is replaced by a new process to free ressources (the default is None, which implies that processes are not replaced).

monte carlo

: int

The number of monte carlo simulations to run in each point of the delta-rho grid (the default is 1).

noise

: {'AWGN', 'AWLN', 'custom'}

The type of noise to use (the default is None, which implies that noiseless measurements are used).

problem_size

: int

The length of the coefficient vector (the default is 800).

rho

: list or tuple

The rho values of the delta-rho grid whose points are used for the monte carlo simulations (the default is [0., 1.]).

seed

: int

The seed used when picking seeds for generating data for the monte carlo simulations

11/2/16 Page 35 of 209

or using the scikit-learn logistic regression solver (the default is None, which implies an arbitrary seed).

SNR

: int or float

The signal-to-noise ratio in decibel to use (the default is 40 dB).

support_distribution

: ndarray

The n x 1 support distribution array to use in the generation of test vectors (the default is None, which implies that the active entries in test vectors are the low index entries).

system matrix

: {'USE', 'RandomDCT2D', 'custom'}

The system matrix to use in the simulation (the default is USE, which implies that the system matrix is drawn from the Uniform Spherical Ensemble).

magni.cs.phase transition. data module

Module providing problem suite instance and noise generation functionality.

The problem suite instances consist of a matrix, A, and a coefficient vector, alpha, with which the measurement vector, y, can be generated (with or without noise from the noise vector e)

Routine listings

generate matrix(m, n)

Generate a matrix belonging to a specific problem suite.

generate_noise(m, n, k)

Generate a noise vector of a specific type.

generate vector(n, k)

Generate a vector belonging to a specific problem suite.

See also:

magni.cs.phase_transition._config

Configuration options.

Notes

The matrices and vectors generated in this module use the numpy.random submodule. Consequently, the calling script or function should control the seed to ensure reproducibility.

The choice of non-zero indices in the coefficient vector is controlled by the configuration option 'support_structure' whereas the distribution of the non-zero coefficients is controlled by the configuration option 'coefficient'.

Examples

For example generate a sample from the USE/Rademacher problem suite:

11/2/16 Page 36 of 209

```
>>> import numpy as np, magni
>>> from magni.cs.phase_transition import _data
>>> m, n, k = 400, 800, 100
>>> A = _data.generate_matrix(m, n)
>>> alpha = _data.generate_vector(n, k)
>>> y = np.dot(A, alpha)
```

Or generate a problem suite instance with "linear" support distribution.

```
>>> support_distrib = np.reshape(np.arange(n, dtype=np.float) + 1, (n, 1))
>>> support_distrib /= np.sum(support_distrib)
>>> magni.cs.phase_transition.config['support_distribution'] = support_distrib
>>> A = _data.generate_matrix(m, n)
>>> alpha = _data.generate_vector(n, k)
>>> y = np.dot(A, alpha)
```

Or generate an AWGN noise vector based on a 40 dB SNR

```
>>> magni.cs.phase_transition.config['noise'] = 'AWGN'
>>> e = _data.generate_noise(m, n, k)
```

magni.cs.phase_transition._data. **generate_matrix**(m, n)

Generate a matrix belonging to a specific problem suite.

The available matrices are

- A random matrix drawn from the Uniform Spherical Ensemble (USE).
- A fixed uniformly row sub-sampled DCT matrix ensemble (RandomDCT2D).
- An option to use custom matrix factory (see notes below).

See Notes for a description of these matrices. Which of the available matrices is used, is specified as a configuration option.

```
Parameters: • m (int) – The number of rows.
```

• **n** (*int*) – The number of columns.

Returns: A (*ndarray*) – The generated matrix.

See also:

magni.cs.phase_transition.config()

Configuration options.

magni.utils.matrices.MatrixCollection()

Fast transform implementation.

Notes

The Uniform Spherical Ensemble:

The matrices of this ensemble have i.i.d. Gaussian entries of mean zero and variance one. Its columns are then normalised to have unit length.

Fixed uniformly row sub-sampled DCT ensemble:

The matrices of this ensemble correspond to the combination of a 2D array subsampled using a uniform point pattern and a 2D Discrete Cosine Transform (DCT) matrix. The combined matrix is implemented as a fast transform with a DCT based on an FFT routine.

11/2/16 Page 37 of 209

Custom matrix factory:

The matrix generation is delegated to the configured custom system matrix factory callable which is expected to take the arguments *m*, *n* and return *A*.

magni.cs.phase transition. data. **generate noise** (m, n, k)

Generate a noise vector of a specific type.

The available types are:

- AWGN: Additive White Gaussian Noise
- AWLN : Additive White Laplacian Noise
- custom : The noise generation delegated to the is configured custom noise factory callable which is expected to take the arguments m, n, k, noise_power.

Which of the available types is used, is specified as a configuration option.

- **Parameters:** m (*int*) The number of rows.
 - **n** (*int*) The number of columns.
 - **k** (int) The number of non-zero coefficients.

Returns:

e (*ndarray*) – The generated noise vector.

See also:

magni.cs.phase_transition.config()

Configuration options.

Notes

The noise power is computed from the configured SNR and the theoretical ensemble variance of the coefficients generated by *generate cofficients*.

magni.cs.phase_transition._data. **generate_vector**(n, k)

Generate a vector belonging to a specific problem suite.

The available ensembles are:

- Gaussian
- Rademacher
- Laplace
- Bernoulli

See Notes for a description of the ensembles. Which of the available ensembles is used, is specified as a configuration option. Note, that the non-zero k non-zero coefficients are the *k* first entries if no support structure specified in the configuration.

- **Parameters: n** (*int*) The length of the vector.
 - **k** (int) The number of non-zero

coefficients.

Returns:

alpha (*ndarray*) – The generated vector.

See also:

Page 38 of 209 11/2/16

magni.cs.phase_transition.config()

Configuration options.

Notes

The Gaussian ensemble:

The non-zero coefficients are drawn from the normal Gaussian distribution.

The Rademacher ensemble:

The non-zero coefficients are drawn from the constant amplitude with random signs ensemble.

The Laplace ensemble:

The non-zero coefficients are drawn from the zero-mean, unit scale Laplace distribution (variance = 2).

The Bernoulli ensemble:

The non-zero coefficients are all equal to one.

magni.cs.phase_transition._simulation module

Module providing the actual simulation functionality.

Routine listings

run(algorithm, path, label)

Simulate a reconstruction algorithm.

See also:

magni.cs.phase_transition._config

Configuration options.

Notes

The results of the simulation are backed up throughout the simulation. In case the simulation is interrupted during execution, the simulation will resume from the last backup point when run again.

magni.cs.phase_transition._simulation. **run** (algorithm, path, label, pre_simulation_hook=None) Simulate a reconstruction algorithm.

The simulation results are stored in a HDF5 database rather than returned by the function.

11/2/16 Page 39 of 209

- Parameters: algorithm (function) A function handle to the reconstruction algorithm.
 - path (str) The path of the HDF5 database where the results should be stored.
 - **label** (*str*) The label assigned to the simulation results
 - pre_simumlation_hook (callable) A handle to a callable which should be run just before the call to the reconstruction algorithm (the default is None, which implies that no pre hook is run).

magni.cs.phase_transition._simulation._simulate(algorithm, ij_tuple, seed, path, pre simulation hook=None)

Run a number of monte carlo simulations in a single delta-rho point.

The result of a simulation is the simulation error distance, i.e., the ratio between the energy of the coefficient residual and the energy of the coefficient vector. The time of the simulation is the execution time of the reconstruction attempt.

- Parameters: algorithm (function) A function handle to the reconstruction algorithm.
 - **ij_tuple** (*tuple*) A tuple (i, j) containing the parameters i, j as listed below.
 - i (int) The delta-index of the point in the delta-rho grid.
 - **j** (*int*) The rho-index of the point in the delta-rho grid.
 - **seed** (*int*) The seed to use in the random number generator when generating the problem suite instances.
 - path (str) The path of the HDF5 backup database.
 - pre simulation hook (callable) A handle to a callable which should be run just before the call to the reconstruction algorithm (the default is None, which implies that no pre hook is run).

See also:

magni.cs.phase_transition._config()

Configuration options.

magni.cs.phase transition. data.generate matrix()

Matrix generation.

magni.cs.phase_transition._data.generate_vector()

Coefficient vector generation.

Notes

The pre simulation hook may be used to setup the simulation to match the specfic simulation parameters, e.g. if an oracle estimator is used in the reconstruction algorithm. The pre simulation hook takes one argument which is the locals() dict.

The following reconstruction statistics are computed:

- time: Measured algorithm run time in seconds.
- dist: Normalised mean squared error (NMSE) (||alpha_hat alpha|| / ||alpha||)^2
- mse: Mean squared error (MSE) 1/n * ||alpha_hat alpha||^2

11/2/16 Page 40 of 209 norm: True vector norm - ||alpha||

magni.cs.phase_transition._util module

Module providing the public function of the magni.cs.phase_transition subpackage.

magni.cs.phase_transition._util. **determine**(algorithm, path, label='default', overwrite=False, pre simulation hook=None)

Determine the phase transition of a reconstruction algorithm.

The phase transition is determined from a number of monte carlo simulations on a deltarho grid for a given problem suite.

- Parameters: algorithm (function) A function handle to the reconstruction algorithm.
 - path (str) The path of the HDF5 database where the results should be stored.
 - label (str) The label assigned to the phase transition (the default is
 - overwrite (bool) A flag indicating if an existing phase transition should be overwritten if it has the same path and label (the default is False).
 - pre simulation hook (callable) A handle to a callable which should be run just before the call to the reconstruction algorithm (the default is None, which implies that no pre hook is run).

See also:

magni.cs.phase_transition._simulation.run()

The actual simulation.

magni.cs.phase_transition._analysis.run()

The actual phase determination.

Notes

The pre simulation hook may be used to setup the simulation to match the specific simulation parameters, e.g. if an oracle estimator is used in the reconstruction algorithm. The *pre_simulation_hook* takes one argument which is the locals() dict.

Examples

An example of how to use this function is provided in the examples folder in the csphase_transition.ipynb ipython notebook file.

magni.cs.phase transition.io module

Module providing input/output functionality for stored phase transitions.

Routine listings

load_phase_transition(path, label='default')

Load the coordinates of a phase transition from a HDF5 file.

11/2/16 Page 41 of 209 magni.cs.phase transition.io.load phase transition(path, label='default', delta=None)

Load the coordinates of a phase transition from a HDF5 file.

This function is used to load the phase transition from the output file generated by magni.cs.phase_transition.determine.

- **Parameters:** path (str) The path of the HDF5 file where the phase transition is stored.
 - **label** (str) The label assigned to the phase transition (the default is 'default').
 - delta (ndarray) The undersampling ratio values to use (the default is None which implies that a uniform spacing of values in the interval (0,1] is assumed.)

Returns:

- **delta** (*ndarray*) The delta values of the phase transition points.
- **rho** (*ndarray*) The rho values of the phase transition points.

See also:

magni.cs.phase_transition.determine()

Phase transition determination.

magni.cs.phase_transition.plotting()

Phase transition plotting.

Examples

An example of how to use this function is provided in the examples folder in the csphase transition.ipynb ipython notebook file.

magni.cs.phase transition.plotting module

Module providing plotting for the magni.cs.phase_transition subpackage.

Routine listings

plot phase transitions(curves, plot I1=True, output path=None)

Function for plotting phase transition boundary curves.

plot phase transition colormap(dist, delta, rho, plot I1=True,

output path=None) Function for plotting reconstruction probabilities in the phase space.

plot phase transition computation times(time, delta, rho, output path=None)

Function for plotting average computation times in the phase space.

magni.cs.phase_transition.plotting. **plot_phase_transitions**(curves, plot_l1=True,

output path=None, legend loc='upper left', errorevery=None, reference curves=None)

Plot of a set of phase transition boundary curves.

The set of phase transition boundary curves are plotted an saved under the output path, if specified. The curves must be a list of dictionaries each having keys delta, rho, and label. delta must be an ndarray of

 δ values in the phase space. *rho* must be an *ndarray* of the corresponding

 ρ values in the phase space. *label* must be a str describing the curve. Optionally, a

11/2/16 Page 42 of 209

curves dictionary may have an yerr key with a corresponding 2-by-"len(delta)" array as value. The first row in this array specifices the location of the upper error bars whereas the second row specifies the location of the lower error bars.

- **Parameters:** curves (*list*) A list of dicts describing the curves to plot.
 - **plot_l1** (bool, optional) Whether or not to plot the theoretical ℓ_1 curve (the default is True).
 - **output_path** (*str, optional*) Path (including file type extension) under which the plot is saved (the default value is None, which implies that the plot is not saved).
 - **legend_loc** (*str*) Location of legend as a matplotlib legend location string (the default is 'upper left', which implies that the legend is placed in the upper left corner of the plot.)
 - errorevery (int) The subsampling of error bars if used (the default value is None, which implies that the default Matplotlib value for errorevery is used.)
 - reference_curves (list, optional) The list of dicts describing the reference curves to plot.

Notes

The plotting is done using matplotlib, which implies that an open figure containing the plot will result from using this function.

Tabulated values of the theoretical

 ℓ_1 phase transition boundary is available at http://people.maths.ox.ac.uk/tanner/polytopes.shtml

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.phase transition.plotting import plot phase transitions
>>> delta = np.array([0.1, 0.2, 0.9])
>>> rho = np.array([0.1, 0.3, 0.8])
>>> curves = [{'delta': delta, 'rho': rho, 'label': 'data1'}]
>>> output path = 'phase transitions.pdf'
>>> plot phase transitions(curves, output path=output path)
```

magni.cs.phase transition.plotting. plot phase transition colormap (dist, delta, rho, plot 11=True, output path=None)

Create a colormap of the phase space reconstruction probabilities.

The delta and rho values span a 2D grid in the phase space. Reconstruction probabilities are then calculated from the *dist* 3D array of reconstruction error distances. The resulting 2D grid of reconstruction probabilities is visualised over the square centers in this 2D grid using a colormap. Values in this grid at lower indicies correspond to lower values of

 δ and

ρ. If plot 11 is True, then the theoretical 11 curve is overlayed the colormap. The colormap is saved under the *output path*, if specified.

11/2/16 Page 43 of 209

- **Parameters: dist** (*ndarray*) A 3D array of reconstruction error distances.
 - delta (ndarray) - δ values used in the 2D grid.
 - rho (ndarray) ρ values used in the 2D grid.
 - plot_l1 (bool) Whether or not to plot the theoretical ℓ_1 curve. (the default is True)
 - **output_path** (*str, optional*) Path (including file type extension) under which the plot is saved (the default value is None which implies, that the plot is not saved).

See also:

```
magni.cs.phase_transition.plotting.plot_phase_transition_computation_times()
```

Plot average computation times.

magni.cs.phase transition.io.load phase transition()

Loading phase transitions from an HDF database.

Notes

The plotting is done using matplotlib, which implies that an open figure containing the plot will result from using this function.

The values in *delta* and *rho* are assumed to be equally spaced.

Due to the centering of the color coded rectangles, they are not necessarily square towards the ends of the intervals defined by delta and rho.

Tabulated values of the theoretical

 ℓ_1 boundary available phase transition is at http://people.maths.ox.ac.uk/tanner/polytopes.shtml

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.phase_transition.plotting import (
   plot_phase_transition_colormap)
\rightarrow delta = np.array([0.2, 0.5, 0.8])
>>> rho = np.array([0.3, 0.6])
>>> dist = np.array([[[1.35e-08, 1.80e-08], [1.08, 1.11]],
... [[1.40e-12, 8.32e-12], [8.57e-01, 7.28e-01]], [[1.92e-13, 1.17e-13],
... [2.10e-10, 1.12e-10]])
>>> out_path = 'phase_transition_cmap.pdf'
>>> plot phase_transition_colormap(dist, delta, rho, output_path=out_path)
```

magni.cs.phase transition.plotting. plot phase transition computation times (time, delta, rho, output path=None)

Create a colormap of the phase space average algorithm computation time.

The delta and rho values span a 2D grid in the phase space. Average computation times are then calculated from the time 3D array of reconstruction times. The resulting

11/2/16 Page 44 of 209

2D grid of average reconstruction times is visualised over the square centers in this 2D grid using a colormap. Values in this grid at lower indicies correspond to lower values of

ρ. The colormap is saved under the *output path*, if specified.

- **Parameters:** time (ndarray) A 3D array of algorithm computation times.
 - delta (ndarray) δ values used in the 2D grid.
 - rho (ndarray) - ρ values used in the 2D grid.
 - **output_path** (*str, optional*) Path (including file type extension) under which the plot is saved (the default value is None which implies, that the plot is not saved).

See also:

```
magni.cs.phase_transition.plotting.plot_phase_transition_colormap()
```

Plot phase transition probabilities.

```
magni.cs.phase_transition.io.load_phase_transition()
```

Loading phase transitions from an HDF database.

Notes

The plotting is done using matplotlib, which implies that an open figure containing the plot will result from using this function.

The values in *delta* and *rho* are assumed to be equally spaced.

Due to the centering of the color coded rectangles, they are not necessarily square towards the ends of the intervals defined by delta and rho.

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.phase transition.plotting import (
  plot phase transition computation times)
>>> delta = np.array([0.2, 0.5, 0.8])
>>> rho = np.array([0.3, 0.6])
>>> times = np.array([[[1.35e-08, 1.80e-08], [1.08, 1.11]],
... [[1.40e-12, 8.32e-12], [8.57e-01, 7.28e-01]], [[1.92e-13, 1.17e-13],
... [2.10e-10, 1.12e-10]]])
>>> out_path = 'computation_times_cmap.pdf'
>>> plot_phase_transition_computation_times(
... times, delta, rho, output_path=out_path)
```

magni.cs.phase_transition.plotting. **plot_extra_curves**(axes, plot 11, reference curves)

Plot any extra curves that might be needed in a phase space plot.

11/2/16 Page 45 of 209

- **Parameters:** axes (matplotlib.axes.Axes) The matplotlib Axes instance on which the theoretical I1 phase transition should be plotted.
 - **plot I1** (bool, optional) Whether or not to plot the theoretical ℓ_1 curve (the default is True).
 - reference curves (list) The list of dicts describing the reference curves to plot.

magni.cs.phase_transition.plotting. **plot_theoretical_l1** (axes)

Plot the theoretical I1 phase transition on the axes.

Parameters: axes (matplotlib.axes.Axes) – The matplotlib Axes instance on which the theoretical I1 phase transition should be plotted.

Notes

The plotted theoretical

 $\ell 1$ phase transition based tabulated values of available is on at http://people.maths.ox.ac.uk/tanner/polytopes.shtml

magni.cs.reconstruction package

Subpackage providing implementations of generic reconstruction algorithms.

Each subpackage provides a family of generic reconstruction algorithms. Thus each subpackage has a config module and a run function which provide the interface of the given family of reconstruction algorithms.

Routine listings

amp

Subpackage providing implementations of Approximate Message Passing (AMP).

gamp

Subpackage providing implementations of Generalised Approximate Message Passing (GAMP).

it

Subpackage providing implementations of Iterative Thresholding (IT).

iht

Subpackage providing implementations of Iterative Hard Thresholding (IHT). (Deprecated)

sl0

Subpackage providing implementations of Smoothed IO Norm (SLO).

Subpackages

magni.cs.reconstruction.amp package

Subpackage providing implementations of Approximate Message Passing (AMP).

Routine listings

config

Configger providing configuration options for this subpackage.

11/2/16 Page 46 of 209 threshold operator

Module providing AMP threshold operators.

run(y, A)

Run the AMP reconstruction algorithm.

stop criterion

Module providing AMP stop criteria.

util

Module providing AMP utilities.

Notes

An implementation of the Donoho, Maleki, Montanari (DMM) AMP using soft thresholding from [1] is provided with threshold options as specified in [2].

The implementation allows for using custom thresholding functions and stop criteria.

References

- [1] D. L. Donoho, A. Maleki, and A. Montanari, "Message-passing algorithms for compressed sensing", *Proceedings of the National Academy of Sciences of the United States of America*, vol. 106, no. 45, p. 18914-18919, Nov. 2009.
- [2] A. Montanari, "Graphical models concepts in compressed sensing" in Compressed Sensing: Theory and Applications, Y. C. Eldar and G. Kutyniok (Ed.), Cambridge University Press, ch. 9, pp. 394-438, 2012.

Submodules

magni.cs.reconstruction.amp._algorithm module

Module providing the core Approximate Message Passing (AMP) algorithm.

Routine listings

run(y, A)

Run the AMP reconstruction algorithm.

See also:

magni.cs.reconstruction.amp._config

Configuration options.

magni.cs.reconstruction.amp.threshold_operator

Threshold operators.

magni.cs.reconstruction.amp.stop_criterion

Stop criteria.

magni.cs.reconstruction.amp.util

Utilities.

Notes

The default configuration of the AMP algorithm provides the Donoho, Maleki, Montanari (DMM) AMP from [1] using the soft threshold with the "residual" thresholding strategy from

11/2/16 Page 47 of 209

[2].

References

- [1] D. L. Donoho, A. Maleki, and A. Montanari, "Message-passing algorithms for compressed sensing", Proceedings of the National Academy of Sciences of the United States of America, vol. 106, no. 45, p. 18914-18919, Nov. 2009.
- [2] A. Montanari, "Graphical models concepts in compressed sensing" in Compressed Sensing: Theory and Applications, Y. C. Eldar and G. Kutyniok (Ed.), Cambridge University Press, ch. 9, pp. 394-438, 2012.

magni.cs.reconstruction.amp._algorithm. $\mathbf{run}(y, A)$

Run the AMP reconstruction algorithm.

- **Parameters:** y (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix, MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.

Returns:

- **alpha** (*ndarray*) The n x 1 reconstructed coefficient vector.
- **history** (*dict, optional*) The dictionary of various measures tracked in the GAMP iterations.

See also:

magni.cs.reconstruction.amp._config()

Configuration options.

magni.cs.reconstruction.amp.threshold_operator()

Threshold operators.

magni.cs.reconstruction.amp.stop_criterion()

Stop criteria.

magni.cs.reconstruction.amp.util()

Utilities.

Notes

Optionally, the algorithm may be configured to save and return the iteration history along with the reconstruction result. The returned history contains the following:

- alpha bar: Coefficient estimates (the reconstruction coefficients).
- MSE: solution mean squared error (if the true solution is known).
- threshold parameters: The state of the threshold parameters.
- stop_criterion : The currently used stop criterion.
- stop_criterion_value : The value of the stop criterion.
- stop_iteration: The iteration at which the algorithm stopped.
- stop reason: The reason for termination of the algorithm.

Examples

For example, recovering a vector from noiseless measurements using AMP with soft

11/2/16 Page 48 of 209 thresholding

```
>>> import numpy as np
>>> from magni.cs.reconstruction.amp import run, config
>>> from magni.cs.reconstruction.amp.util import theta mm
>>> np.random.seed(seed=6028)
>>> k, m, n = 10, 200, 400
>>> A = 1 / np.sqrt(m) * np.random.randn(m, n)
>>> alpha = np.zeros((n, 1))
>>> alpha[:k] = np.random.normal(scale=2, size=(k, 1))
>>> np_printoptions = np.get_printoptions()
>>> np.set_printoptions(suppress=True, threshold=k+2)
>>> alpha[:k + 2]
array([[ 1.92709461],
    [ 0.74378508],
    [-3.2418159],
    [-1.32277347],
    [0.90118].
    [-0.19157262],
    [0.82855712],
    [0.24817994],
    [-1.43034777],
    [-0.21232344],
    [ 0.
            ],
    [ 0.
           ]])
>>> y = A.dot(alpha)
>>> threshold_params = {'threshold_level_update_method': 'residual',
... 'theta': theta_mm(float(m) / n), 'tau_hat_sq': 1.0}
>>> config['threshold_parameters'] = threshold_params
>>> alpha_hat = run(y, A)
>>> alpha_hat[:k + 2]
array([[ 1.92612707],
    [ 0.74185284],
    [-3.24179313],
    [-1.32200929],
    [0.90089208],
    [-0.19097308],
    [ 0.82658038],
    [0.24515825],
    [-1.42980997],
    [-0.2111469].
    [0.00002815],
    [0.00047293]])
>>> np.sum(np.abs(alpha - alpha hat) > 1e-2)
```

or recover the same vector returning a history comparing the pr. iteration solution to the true vector

11/2/16 Page 49 of 209

```
>>> config['report_history'] = True
>>> config['true_solution'] = alpha
>>> alpha_hat, history = run(y, A)
>>> alpha hat[:k + 2]
array([[ 1.92612707],
    [0.74185284],
    [-3.24179313],
    [-1.32200929],
    [0.90089208],
    [-0.19097308],
    [ 0.82658038],
    [ 0.24515825],
    [-1.42980997],
    [-0.2111469],
    [ 0.00002815].
    [ 0.00047293]])
>>> np.array(history['MSE']).reshape(-1, 1)[1:11]
array([[ 0.01273323],
    [0.0053925],
    [ 0.00276334],
    [ 0.00093472],
    [ 0.0004473 ],
    [0.00021142],
    [0.00009618],
    [0.00004371],
    [0.00002275],
    [ 0.00001053]])
```

or recover the same vector using the median based update method

```
>>> config['report history'] = False
>>> threshold params['threshold level update method'] = 'median'
>>> config['threshold_parameters'] = threshold_params
>>> alpha_hat = run(y, A)
>>> alpha_hat[:k + 2]
array([[ 1.92577079],
    [0.7413161],
    [-3.24131125],
    [-1.32044146],
    [ 0.89998407],
    [-0.18981235],
    [0.82584016],
    [0.24445332],
    [-1.42934708],
    [-0.21043221],
    [ 0.
    [ 0.
>>> np.sum(np.abs(alpha - alpha_hat) > 1e-2)
>>> np.set_printoptions(**np_printoptions)
```

magni.cs.reconstruction.amp._algorithm._get_amp_initialisation(y, A)

Return an initialisation of the AMP algorithm.

- **Parameters:** y (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix, MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.

Returns:

init (*dict*) – The initialisation of the the AMP algorithm.

magni.cs.reconstruction.amp._config module

Module providing configuration options for the magni.cs.reconstruction.amp subpackage.

11/2/16 Page 50 of 209

See also:

magni.cs.reconstruction._config.Configger

The Configger class used.

Notes

This module instantiates the *Configger* class provided by magni.cs.reconstruction._config.Configger. The configuration options are the following:

iterations

: int

The maximum number of iterations to do (the default is 300).

precision float

: {np.float, np.float16, np.float32, np.float64, np.float128,

np.complex64, np.complex128, np.complex256} The floating point precision used for the computations (the default is np.float64).

report history

: bool

The indicator of whether or not to return the progress history along with the result (the default is False).

stop criterion

: magni.utils.validation.types.StopCriterion

The stop criterion to use in the iterations (the default is magni.cs.reconstruction.gamp.stop criterion.MSEConvergence).

threshold

: magni.utils.validation.types.ThresholdOperator

The threshold operator to use (the default is magni.cs.reconstruction.amp.threshold_operator.SoftThreshold).

threshold parameters

: dict

The parameters used in the threshold operator (no default is provided, which implies that this must be specified by the user).

tolerance

: float

The least acceptable stop criterion tolerance to break the interations (the default is 1e-6).

true solution

: ndarray or None

The true solution to allow for tracking the convergence of the algorithm in the artificial setup where the true solution is known a-priori (the default is None, which implies that no true solution tracking is used).

warm start

: ndarray

The initial guess of the solution vector (alpha_bar) (the default is None, which implies that alpha bar is taken to be a vector of zeros).

magni.cs.reconstruction.amp.stop criterion module

11/2/16 Page 51 of 209

Module providing stop criteria for the Approximate Message Passing (AMP) algorithm.

Routine listings

ValidatedStopCriterion(magni.utils.validation.types.StopCriterion)

A base class for validated magni.cs.reconstruction.amp stop criteria.

MSEConvergence(ValidatedStopCriterion)

A mean square error (MSE) convergence stop criterion.

NormalisedMSEConvergence(ValidatedStopCriterion)

A normalised mean squaure error (NMSE) convergence stop criterion.

Residual(ValidatedStopCriterion)

A residual based stop criterion.

ResidualMeasurementsRatio(ValidatedStopCriterion)

A residual-measurements-ratio based stop criterion.

class magni.cs.reconstruction.amp.stop_criterion. ValidatedStopCriterion(var)

Bases: magni.utils.validation.types.StopCriterion

A base class for validated magni.cs.reconstruction.amp stop criteria.

Parameters: var (*dict*) – The stop criterion state variables.

compute (var)

Compute the stop criterion value.

Parameters: var (*dict*) – The variables used in computing of the stop criterion

value.

• **stop** (*bool*) – The indicator of whether or not the stop criterion is

satisfied.

• **value** (*float*) – The stop criterion value.

Notes

This method honors magni.utils.validation.enable allow validate once.

class magni.cs.reconstruction.amp.stop criterion. **MSEConvergence** (*var*)

Bases: magni.cs.reconstruction.amp.stop_criterion.ValidatedStopCriterion

A mean square error (MSE) convergence stop criterion.

Parameters: var (*dict*) – The stop criterion state

variables.

Notes

The following state variables are used in this stop criterion:

- A
- tolerance

11/2/16 Page 52 of 209

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.reconstruction.amp.stop_criterion import MSEConvergence
>>> state = {'tolerance': 1e-3, 'A': np.ones((10, 10))}
>>> variables = {'alpha_bar_prev': np.ones((10, 1)),
... 'alpha_bar': np.arange(10).reshape(10, 1)}
>>> MSE = MSEConvergence(state)
>>> stop, val = MSE.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute(var)

Compute the MSE convergence stop criterion value.

The AMP algorithm should converge to a fixed point. This criterion is based on the mean squared error of the difference between the proposed solution in this iteration and the proposed solution in the previous solution.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

Returns:

- stop (bool) The indicator of whether or not the stop criterion is satisfied.
- value (float) The stop criterion value.

class magni.cs.reconstruction.amp.stop criterion. NormalisedMSEConvergence (var)

Bases: magni.cs.reconstruction.amp.stop_criterion.ValidatedStopCriterion

A normalised mean squaure error (NMSE) convergence stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

tolerance

Examples

For example,

11/2/16 Page 53 of 209

```
>>> import numpy as np
>>> from magni.cs.reconstruction.amp.stop_criterion import (
... NormalisedMSEConvergence)
>>> state = {'tolerance': 1e-3}
>>> variables = {'alpha_bar_prev': np.ones((10, 1)),
... 'alpha_bar': np.arange(10).reshape(10, 1)}
>>> NMSE = NormalisedMSEConvergence(state)
>>> stop, val = NMSE.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute (var)

Compute the normalised MSE convergence stop criterion value.

The AMP algorithm should converge to a fixed point. This criterion is based on the mean squared error of the difference between the proposed solution in this iteration and the proposed solution in the previous solution normalised by the mean squared error of the proposed solution in the previous iteration.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

• **stop** (*bool*) – The indicator of whether or not the stop criterion is satisfied.

• **value** (*float*) – The stop criterion value.

class magni.cs.reconstruction.amp.stop criterion. Residual (var)

Bases: magni.cs.reconstruction.amp.stop_criterion.ValidatedStopCriterion

A residual based stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

- y
- tolerance
- A

Examples

For example,

11/2/16 Page 54 of 209

```
>>> import numpy as np
>>> from magni.cs.reconstruction.amp.stop_criterion import Residual
>>> state = {'tolerance': 1e-3, 'y': np.ones((10, 1)),
... 'A': np.ones((10,10))}
>>> variables = {'A_dot_alpha_bar': np.arange(10).reshape(10, 1)}
>>> Res = Residual(state)
>>> stop, val = Res.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute(var)

Compute the residual stop criterion value.

If the noise level is (approximately) known, the AMP iterations may be stopped once the residual is on the order of the noise level. This stopping criterion is based on the mean squared error of the residual.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

Returns:

- stop (bool) The indicator of whether or not the stop criterion is satisfied.
- value (float) The stop criterion value.

class magni.cs.reconstruction.amp.stop_criterion. ResidualMeasurementsRatio (var)

 $Bases: magni.cs. reconstruction. amp. stop_criterion. Validated Stop Criterion\\$

A residual-measurements-ratio based stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

- y
- tolerance

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.reconstruction.amp.stop_criterion import (
... ResidualMeasurementsRatio)
>>> state = {'tolerance': 1e-3, 'y': np.ones((10, 1))}
>>> variables = {'A_dot_alpha_bar': np.arange(10).reshape(10, 1)}
>>> ResMeasRat = ResidualMeasurementsRatio(state)
>>> stop, val = ResMeasRat.compute(variables)
>>> stop
False
>>> np.round(val, 2)
14.32
```

compute (var)

11/2/16 Page 55 of 209

Compute the residual-measurements-ratio stop criterion value.

If the noise level is (approximately) known, the AMP iterations may be stopped once the residual is on the order of the noise level. This stopping criterion is based on ratio of the mean squared error of the residual to the mean squared error of the measurements.

Parameters: var (dict) – Dictionary of variables used in calculating of the stop

criterion.

Returns:

 stop (bool) – The indicator of whether or not the stop criterion is satisfied.

• value (float) – The stop criterion value.

magni.cs.reconstruction.amp.threshold_operator module

Module providing threshold functions for the Approximate Message Passing (AMP) algorithm.

Routine listings

ValidatedThresholdOperator(magni.utils.validation.types.ThresholdOperator)

A base class for validated magni.cs.reconstruction.amp threshold operator

SoftThreshold(ValidatedThresholdOperator)

A soft threshold operator.

class magni.cs.reconstruction.amp.threshold operator. ValidatedThresholdOperator(var)

Bases: magni.utils.validation.types.ThresholdOperator

A base class for validated magni.cs.reconstruction.amp threshold operator

Parameters: var (*dict*) – The threshold operator state variables.

compute_deriv_threshold(var)

Compute the entrywise derivative threshold.

Parameters: var (*dict*) – The variables used in computing the derivative

threshold.

Returns: eta_deriv (ndarray) – The computed entrywise derivative threshold.

Notes

This method honors *magni.utils.validation.enable_allow_validate_once*.

compute_threshold(var)

Compute the entrywise threshold.

Parameters: var (*dict*) – The variables used in computing the

threshold.

Returns: eta (*ndarray*) – The computed entrywise threshold.

Notes

11/2/16 Page 56 of 209

This method honors magni.utils.validation.enable allow validate once.

update threshold level (var)

Update the threshold level state.

Parameters: var (*dict*) – The variables used in computing the threshold level update.

Notes

This method honors magni.utils.validation.enable allow validate once.

class magni.cs.reconstruction.amp.threshold_operator. **SoftThreshold**(*var*)

Bases: magni.cs.reconstruction.amp.threshold_operator.ValidatedThresholdOperator

A soft threshold operator.

This soft threshold operator is based on the description of it and its use in AMP as given in [1] with corrections from [2].

- Parameters: threshold level update method ({'residual', 'median'}) The method to use for updating the threshold level.
 - theta (float) The tunable regularisation parameter in the threshold level.
 - **tau hat sq** (*float*) The mean squared error of the (approximated) un-thresholded estimate used to determine the threshold level.

Notes

The above Parameters are the threshold parameters that must be passed in a var dict to the threshold constructor.

References

- [1] A. Montanari, "Graphical models concepts in compressed sensing" in Compressed Sensing: Theory and Applications, Y. C. Eldar and G. Kutyniok (Ed.), Cambridge University Press, ch. 9, pp. 394-438, 2012.
- [2] J. T. Parker, "Approximate Message Passing Algorithms for Generalized Bilinear Inference", PhD Thesis, Graduate School of The Ohio State University, 2014

compute deriv threshold(var)

Compute the entrywise derivative soft threshold.

Parameters: var (*dict*) – The variables used in computing the derivative

threshold.

Returns: **eta_deriv** (*ndarray*) – The computed entrywise derivative threshold.

compute_threshold(var)

Compute the entrywise soft threshold.

Parameters: var (*dict*) – The variables used in computing the threshold.

11/2/16 Page 57 of 209

```
Returns: eta (ndarray) – The computed entrywise threshold.
```

update_threshold_level(var)

Update the threshold level state.

```
Parameters: var (dict) – The variables used in computing the threshold level update.
```

magni.cs.reconstruction.amp.util module

Module providing utilities for the Approximate Message Passing (AMP) algorithm.

Routine listings

```
theta_mm(theta)
```

Return the minimax optimal value of the AMP theta tunable parameter.

```
magni.cs.reconstruction.amp.util.theta_mm (delta)
```

Return the minimax optimal value of the AMP theta tunable parameter.

A simple lookup table is used to determine the minimax optimal value of the AMP theta tunable parameter as specified in [1]. Specifically, the formula in the top left column on page 6 is used without the 1/sqrt(delta) factor (as specified in the SI Appendix) and with kappa=2.

```
Parameters: delta (float) – The undersampling ratio.

Returns: theta_mm (float) – The minimax optimal value of the theta parameter.
```

References

[1] D. L. Donoho, A. Maleki, and A. Montanari, "Message-passing algorithms for compressed sensing", *Proceedings of the National Academy of Sciences of the United States of America*, vol. 106, no. 45, p. 18914-18919, Nov. 2009.

Examples

For example, get the optimal parameters for a few undersampling ratios

```
>>> from magni.cs.reconstruction.amp.util import theta_mm
>>> theta_mm(0.1)
1.736
>>> theta_mm(0.2121)
1.378
>>> theta_mm(0.001)
3.312
>>> theta_mm(0.999)
0.032
>>> theta_mm(0.0001)
3.312
>>> theta_mm(0.0001)
3.312
>>> theta_mm(0.09999)
0.032
```

magni.cs.reconstruction.gamp package

Subpackage providing implementations of Generalised Approximate Message Passing

11/2/16 Page 58 of 209

(GAMP).

Routine listings

channel initialisation

Module providing functionality for initialisations of GAMP channels.

config

Configger providing configuration options for this subpackage.

input_channel

Module providing GAMP input channels.

output channel

Module providing GAMP output channels.

run(y, A, A asq=None)

Run the GAMP reconstruction algorithm.

stop_criterion

Module providing GAMP stop criteria.

Notes

Implementations of Mimimum Mean Squared Error (MMSE) Generalised Approximate Message Passing (GAMP) from [1], [2] based on description of it in [3] are available. The GAMP is a generalisation of the Approximate Message Passing (AMP) algorithm derived independently by Donoho et al. [4] and Krzakala et al. [5], [6].

This implementation allows custom input- and output channels as well as the use of sum approximations of the squared system matrix as detailed in [2], [5]. Furthermore, a simple damping option is available based on the description in [7] (see also [8] for more details on damping in GAMP).

References

- [1] S. Rangan, "Generalized Approximate Message Passing for Estimation with Random Linear Mixing", *in IEEE International Symposium on Information Theory (ISIT)*, St. Petersburg, Russia, Jul. 31 Aug. 5, 2011, pp. 2168-2172.
- [2] (1, 2) S. Rangan, "Generalized Approximate Message Passing for Estimation with Random Linear Mixing", arXiv:1010.5141v2, pp. 1-22, Aug. 2012.
- [3] J. T. Parker, "Approximate Message Passing Algorithms for Generalized Bilinear Inference", PhD Thesis, Graduate School of The Ohio State University, 2014
- [4] D.L. Donoho, A. Maleki, and A. Montanari, "Message-passing algorithms for compressed sensing", *Proceedings of the National Academy of Sciences of the United States of America*, vol. 106, no. 45, pp. 18914-18919, Nov. 2009.
- [5] (1, 2) F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Probabilistic reconstruction in compressed sensing: algorithms, phase diagrams, and threshold achieving matrices", *Journal of Statistical Mechanics: Theory and Experiment*, vol. P08009, pp. 1-57, Aug. 2012.
- [6] F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Statistical-Physics-Based Reconstruction in Compressed Sensing", *Physics Review X*, vol. 2, no. 2, pp. (021005-1)-(021005-18), May 2012.

11/2/16 Page 59 of 209

- [7] S. Rangan, P. Schniter, and A. Fletcher. "On the Convergence of Approximate Message Passing with Arbitrary Matrices", in *IEEE International Symposium on Information Theory (ISIT)*, pp. 236-240, Honolulu, Hawaii, USA, Jun. 29 Jul. 4, 2014.
- [8] J. Vila, P. Schniter, S. Rangan, F. Krzakala, L. Zdeborova, "Adaptive Damping and Mean Removal for the Generalized Approximate Message Passing Algorithm", in IEEE International Conference on Acoustics, Speech, and Signal Processing (ICASSP), South Brisbane, Queensland, Australia, Apr. 19-24, 2015, pp. 2021-2025.

Submodules

magni.cs.reconstruction.gamp._algorithm module

Module providing the core Generalised Approximate Message Passing (GAMP) algorithm.

Routine listings

run(y, A, A asq=None)

Run the GAMP reconstruction algorithm.

See also:

magni.cs.reconstruction.gamp._config

Configuration options.

magni.cs.reconstruction.gamp.input_channel

Available input channels.

magni.cs.reconstruction.gamp.output_channel

Available output channels.

magni.cs.reconstruction.gamp.stop_criterion

Available stop critria.

Notes

The default configuration of the GAMP algorithm provides the s-GB AMP algorithm from[1] using an MSE convergence based stop criterion. Both the input channel, the output channel, and the stop criterion may be changed.

This implementation allows for the use of sum approximations of the squared system matrix as detailed in [1] and [2]. Furthermore, a simple damping option is available based on the description in [3] (see also [4] for more details on damping in GAMP).

References

- [1] (1, 2) F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Probabilistic reconstruction in compressed sensing: algorithms, phase diagrams, and threshold achieving matrices", *Journal of Statistical Mechanics: Theory and Experiment*, vol. P08009, pp. 1-57, Aug. 2012.
- [2] S. Rangan, "Generalized Approximate Message Passing for Estimation with Random Linear Mixing", arXiv:1010.5141v2, pp. 1-22, Aug. 2012.
- [3] S. Rangan, P. Schniter, and A. Fletcher. "On the Convergence of Approximate Message Passing with Arbitrary Matrices", *in IEEE International Symposium on Information Theory (ISIT)*, pp. 236-240, Honolulu, Hawaii, USA, Jun. 29 Jul. 4, 2014.

11/2/16 Page 60 of 209

[4] J. Vila, P. Schniter, S. Rangan, F. Krzakala, L. Zdeborova, "Adaptive Damping and Mean Removal for the Generalized Approximate Message Passing Algorithm", in IEEE International Conference on Acoustics, Speech, and Signal Processing (ICASSP), South Brisbane, Queensland, Australia, Apr. 19-24, 2015, pp. 2021-2025.

magni.cs.reconstruction.gamp._algorithm. **run** (*y*, *A*, *A* asq=None)

Run the GAMP reconstruction algorithm.

- **Parameters: y** (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix, MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.
 - A_asq (ndarray or magni.utils.matrices.{Matrix, MatrixCollection} or None) - The m x n matrix which is the entrywise absolute value squared product of the measurement matrix and the dictionary matrix (the default is None, which implies that a sum approximation is used).

Returns:

- **alpha** (*ndarray*) The n x 1 reconstructed coefficient vector.
- history (dict, optional) The dictionary of various measures tracked in the GAMP iterations.

See also:

magni.cs.reconstruction.gamp._config()

Configuration options.

magni.cs.reconstruction.gamp.input_channel()

Input channels.

magni.cs.reconstruction.gamp.output_channel()

Output channels.

magni.cs.reconstruction.gamp.stop_criterion()

Stop criteria.

Notes

Optionally, the algorithm may be configured to save and return the iteration history along with the reconstruction result. The returned history contains the following:

- alpha_bar : Mean coefficient estimates (the reconstruction coefficients).
- alpha_tilde : Variance coefficient estimates.
- MSE: solution Mean squared error (if the true solution is known).
- input channel parameters : The state of the input channel.
- output_channel_parameters : The state of the output channel.
- stop_criterion : The currently used stop criterion.
- stop criterion value: The value of the stop criterion.
- stop iteration: The iteration at which the algorithm stopped.
- stop reason: The reason for termination of the algorithm.

Examples

For example, recovering a vector from AWGN noisy measurements using GAMP

11/2/16 Page 61 of 209

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp import run, config
>>> np.random.seed(seed=6028)
>>> k, m, n = 10, 200, 400
>>> tau = float(k) / n
>>> A = 1 / np.sqrt(m) * np.random.randn(m, n)
>>> A_asq = np.abs(A)**2
>>> alpha = np.zeros((n, 1))
>>> alpha[:k] = np.random.normal(scale=2, size=(k, 1))
>>> np_printoptions = np.get_printoptions()
>>> np.set_printoptions(suppress=True, threshold=k+2)
>>> alpha[:k + 2]
array([[ 1.92709461],
    [ 0.74378508],
    [-3.2418159],
    [-1.32277347],
    [0.90118],
    [-0.19157262],
    [0.82855712],
    [ 0.24817994],
    [-1.43034777],
    [-0.21232344],
    [ 0.
            ]])
    [0.
>>> sigma = 0.15
>>> y = A.dot(alpha) + np.random.normal(scale=sigma, size=(A.shape[0], 1))
>>> input_channel_params = {'tau': tau, 'theta_bar': 0, 'theta_tilde': 4,
... 'use_em': False}
>>> config['input_channel_parameters'] = input_channel_params
>>> output_channel_params = {'sigma_sq': sigma**2,
... 'noise_level_estimation': 'fixed'}
>>> config['output_channel_parameters'] = output_channel_params
>>> alpha_hat = run(y, A, A_asq)
>>> alpha_hat[:k + 2]
array([[ 1.93810961],
    [ 0.6955502 ],
    [-3.39759349],
    [-1.35533562],
    [1.10524227],
    [-0.00594848],
    [0.79274671],
    [0.04895264],
    [-1.08726071],
    [-0.00142911],
    [ 0.00022861],
    [-0.00004272]])
>>> np.sum(np.abs(alpha - alpha_hat) > sigma * 3)
```

or recover the same vector returning a history comparing the pr. iteration solution to the true vector and printing the A_asq details

11/2/16 Page 62 of 209

```
>>> config['report_A_asq_setup'] = True
>>> config['report_history'] = True
>>> config['true_solution'] = alpha
>>> alpha_hat, history = run(y, A, A_asq)
GAMP is using the A_asq: [[ 0.024 ..., 0.002]
[0...., 0.014]]
The sum approximation method is: None
>>> alpha_hat[:k + 2]
array([[ 1.93810961],
    [ 0.6955502 ],
    [-3.39759349],
    [-1.35533562],
    [1.10524227],
    [-0.00594848],
    [0.79274671],
    [ 0.04895264],
    [-1.08726071],
    [-0.00142911],
    [ 0.00022861],
    [-0.00004272]])
>>> np.array(history['MSE']).reshape(-1, 1)[1:11]
array([[ 0.04562729],
    [ 0.01328304],
    [ 0.00112098],
    [ 0.00074968],
    [0.00080175],
    [ 0.00076615],
    [ 0.00077043]])
```

or recover the same vector using sample variance AWGN noise level estimation

```
>>> config['report_A_asq_setup'] = False
>>> config['report_history'] = False
>>> output_channel_params['noise_level_estimation'] = 'sample_variance'
>>> config['output_channel_parameters'] = output_channel_params
>>> alpha_hat = run(y, A, A_asq)
>>> alpha_hat[:k + 2]
array([[ 1.94820622],
    [ 0.72162206],
    [-3.39978431],
    [-1.35357001],
    [1.10701779],
    [-0.00834467],
    [0.79790879],
    [0.08441384],
    [-1.08946306],
    [-0.0015894],
    [ 0.00020561].
    [-0.00003623]])
>>> np.sum(np.abs(alpha - alpha_hat) > sigma * 3)
```

or recover the same vector using median AWGN noise level estimation

11/2/16 Page 63 of 209

```
>>> output channel params['noise level estimation'] = 'median'
>>> config['output_channel_parameters'] = output_channel_params
>>> alpha_hat = run(y, A, A_asq)
>>> alpha hat[:k + 2]
array([[ 1.93356483],
    [0.65232347],
    [-3.39440429],
    [-1.35437724],
    [1.10312573],
    [-0.0050555],
    [0.78743162],
    [ 0.03616397],
    [-1.08589927],
    [-0.00136802],
    [ 0.00024121].
    [-0.00004498]])
>>> np.sum(np.abs(alpha - alpha_hat) > sigma * 3)
```

or recover the same vector learning the AWGN noise level using expectation maximization (EM)

```
>>> output channel params['noise level estimation'] = 'em'
>>> config['output_channel_parameters'] = output_channel_params
>>> alpha_hat = run(y, A, A_asq)
>>> alpha_hat[:k + 2]
array([[ 1.94118089],
    [0.71553983],
    [-3.40076165],
    [-1.35662005],
    [1.1099417].
    [-0.00688125],
     0.79442879],
    [0.06258856],
    [-1.08792606],
    [-0.00148811],
    [0.00022266]
    [-0.00003785]])
>>> np.sum(np.abs(alpha - alpha_hat) > sigma * 3)
```

```
>>> np.set printoptions(**np printoptions)
```

magni.cs.reconstruction.gamp._algorithm. _get_gamp_initialisation (y, A, A_asq)

Return an initialisation of the GAMP algorithm.

- **Parameters:** y (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix, MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.
 - A_asq (ndarray or magni.utils.matrices.{Matrix, MatrixCollection} or None) – The m x n matrix which is the entrywise absolute value squared product of the measurement matrix and the dictionary matrix (the default is None, which implies that a sum approximation is used).

Returns:

init (*dict*) – The initialisation of the the GAMP algorithm.

magni.cs.reconstruction.gamp. config module

Module providing configuration options for the magni.cs.reconstruction.gamp subpackage.

```
See also:
```

11/2/16 Page 64 of 209 magni.cs.reconstruction._config.Configger

The Configger class used

Notes

This module instantiates the *Configger* class provided by magni.cs.reconstruction._config.Configger. The configuration options are the following:

damping

: float

The damping applied to the variable side updates (the default is 0.0).

input channel

: magni.utils.validation.types.MMSEInputChannel

The input channel to use (the default is magni.cs.reconstruction.gamp.input channel.IIDBG).

input channel parameters

: dict

The parameters used in the input channel (no default is provided, which implies that this must be specified by the user).

iterations

: int

The maximum number of iterations to do (the default is 300).

output channel

: magni.utils.validation.types.MMSEOutputChannel

The output channel to use (the default is magni.cs.reconstruction.gamp.output_channel.AWGN).

output channel parameters

: dict

The parameters used in the output channel (no default is provided, which implies that this must be specified by the user).

precision float

: {np.float, np.float16, np.float32, np.float64, np.float128,

np.complex64, np.complex128, np.complex256} The floating point precision used for the computations (the default is np.float64).

report_A_asq_setup

: bool

The indicator of whether or not to print the A asq details (the default is False).

report history

: bool

The indicator of whether or not to return the progress history along with the result (the default is False).

stop criterion

: magni.utils.validation.types.StopCriterion

The stop criterion to use in the iterations (the default is magni.cs.reconstruction.gamp.stop_criterion.MSEConvergence).

sum_approximation_constant

11/2/16 Page 65 of 209

: dict

The method and constant used in a sum approximation of the squared system transform (the default is {'rangan': 1.0}, which implies that Rangan's uniform variance methods with the system matrix adapted $||A||_F^2/(m^*n)$ constant is used).

tolerance

: float

The least acceptable stop criterion tolerance to break the interations (the default is 1e-

true solution

: ndarray or None

The true solution to allow for tracking the convergence of the algorithm in the artificial setup where the true solution is known a-priori (the default is None, which implies that no true solution tracking is used).

warm start

: list or tuple

The collection containing the initial guess of the solution vector (alpha bar) and the solution variance vector (alpha tilde) (the default is None, which implies that alpha bar is taken to be a vector of zeros and alpha tilde is taken to be a vector of ones).

magni.cs.reconstruction.gamp.channel initialisation module

Module providing utility functions for initilisation of in- and output channels in the Generalised Approximate Message Passing (GAMP) algorithm.

Routine listings

get em bg amp initialisation(problem params, method='vila')

Get initial parameters for EM Bernoulli-Guassian AMP.

rho_se(delta, zeta, resolution=1000)

Return the theoretical noiseless LASSO phase transition.

magni.cs.reconstruction.gamp.channel initialisation.get em bg amp initialisation(problem params, method='vila')

Get initial parameters for EM Bernoulli-Guassian AMP.

If the initialisation *method* is *vila* then the scheme from [1] is used. If it is *krzakala* then the scheme from [2] is used.

- **Parameters:** problem_params (dict) The problem parameters used to compute the initialisation.
 - **method** (*str*) The initialisation method to use.

Returns:

- tau (float) The initial sparsity level.
- **sigma sq** (*float*) The initial AWGN noise level.
- theta_tilde (float) The initial Gaussian prior variance.

See also:

magni.cs.reconstruction.gamp._output_channel.wrap_calculate_using_AWGN() Related output channel.

magni.cs.reconstruction.gamp._input_channels.wrap_calculate_using_iidsGB() Related input channel.

11/2/16 Page 66 of 209

Notes

Independently of the choice of method, the problem params are:

- y: the measurements
- A: the system matrix

If method is vila, one must also specify:

• SNR: the signal-to-noise ratio

References

- [1] J. P. Vila and P. Schniter, "Expectation-Maximization Gaussian-Mixture Approximate Message Passing", *IEEE Transactions on Signal Processing*, 2013, vol. 61, no. 19, pp. 4658-4672, Oct. 2013.
- [2] F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Probabilistic reconstruction in compressed sensing: algorithms, phase diagrams, and threshold achieving matrices", *Journal of Statistical Mechanics: Theory and Experiment*, vol. P08009, pp. 1-57, Aug. 2012.

Examples

For example, get the "vila" initialisation for a SNR of 100

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp.channel_initialisation import ... get_em_bg_amp_initialisation
>>> np.random.seed(6012)
>>> A = np.random.randn(20, 40)
>>> y = np.random.randn(20, 1)
>>> problem_params = {'A': A, 'y': y, 'SNR': 100}
>>> init_tup = get_init(problem_params, method='vila')
>>> [round(float(elem), 3) for elem in init_tup]
[0.193, 0.01, 0.123]
```

or get the corresponding "krzakala" initialisation

```
>>> del problem_params['SNR']
>>> init_tup = get_init(problem_params, method='krzakala')
>>> [round(float(elem), 3) for elem in init_tup]
[0.05, 1.0, 0.479]
```

magni.cs.reconstruction.gamp.channel_initialisation.rho_se(delta, zeta, resolution=1000)

Return the theoretical noiseless LASSO phase transition.

```
    delta (float) – The under sampling ratio.
    zeta ({1, 2}) – The "problem" to get the phase transition for.
    resolution (int) – The resolution used in the brute force optimisation.
    Returns:
```

Notes

11/2/16 Page 67 of 209

The theoretical noiseless LASSO phase transition is computed based on eq. 5 in[3]. A simple brute force optimisation with the specified *resolution* of that expression is used to find the phase transition. The "problems", for which the phase transition may be computed, are:

- 1. Sparse nonnegative vectors
- 2. Sparse signed vectors

References

[3] D.L. Donoho, A. Maleki, and A. Montanari, "Message-passing algorithms for compressed sensing", *Proceedings of the National Academy of Sciences of the United States of America*, vol. 106, no. 45, pp. 18914-18919, Nov. 2009.

Examples

For example, find a phase transition value for a sparse signed vector:

```
>>> from magni.cs.reconstruction.gamp.channel_initialisation import rho_se
>>> round(float(rho_se(0.19, 2)), 3)
0.238
```

or find the corresponding value for a sparse nonnegative vector

```
>>> round(float(rho_se(0.19, 1)), 3) 0.318
```

and a few more examples

```
>>> round(float(rho_se(0.0, 1)), 3)
0.0
>>> round(float(rho_se(0.0, 2)), 3)
0.0
>>> round(float(rho_se(0.5, 1)), 3)
0.558
>>> round(float(rho_se(0.5, 2)), 3)
0.386
>>> round(float(rho_se(1.0, 1)), 3)
0.95
>>> round(float(rho_se(1.0, 2)), 3)
0.95
```

magni.cs.reconstruction.gamp.input channel module

Module providing input channel functions for the Generalised Approximate Message Passing (GAMP) algorithm.

Routine listings

Validated MMSEInput Channel (magni.utils.validation.types. MMSEInput Channel)

A base class for validated magni.cs.reconstruction.gamp input channels.

IIDBG(ValidatedMMSEInputChannel)

An i.i.d. Bernoulli Gauss MMSE input channel.

IIDsGB(ValidatedMMSEInputChannel)

An i.i.d. sparse Gauss Bernoulli MMSE input channel.

11/2/16 Page 68 of 209

class magni.cs.reconstruction.gamp.input_channel. ValidatedMMSEInputChannel(var)

Bases: magni.utils.validation.types.MMSEInputChannel

A base class for validated magni.cs.reconstruction.gamp input channels.

Parameters: var (*dict*) – The input channel state variables.

compute (var)

Compute the input channel value.

Parameters: var (*dict*) – The variables used in computing of the input channel value.

Returns: • **mean** (*ndarray*) – The computed input channel mean.

• **variance** (*ndarray*) – The computed input channel variance.

Notes

This method honors magni.utils.validation.enable allow validate once.

class magni.cs.reconstruction.gamp.input_channel.**IIDBG**(*var*)

Bases: magni.cs.reconstruction.gamp.input_channel.ValidatedMMSEInputChannel

An i.i.d. Bernoulli Gauss MMSE input channel.

This channel is based on equations (6), (7) in [3] and allows for using Expectation Maximization (EM) for learning the channel parameters as detailed in equations (19), (25), and (32) in [3] (see also [4] for an introduction to EM for GAMP).

- Parameters: tau (float or int) The prior signal "density" (fraction of Gauss to Bernouilli).
 - theta bar (float or int) The prior Gaussian mean.
 - theta tilde (float or int) The prior Gaussian variance.
 - **use_em** (*bool*) The indicator of whether or not to use Expectation Maximazion (EM) to learn the prior parameters.

Notes

The above Parameters are the input channel parameters that must be passed in a var dict to the channel constructor.

If use em is True, the values given for tau, theta bar, and theta tilde are used for initialialisation.

This channel is theoretically equivalent to the IIDsGB channel. However, due to numerical subtleties, it may give different results.

References

11/2/16 Page 69 of 209

- [3] (1, 2) J. Vila, J. and P. Schniter, "Expectation-Maximization Bernoulli-Gaussian Approximate Message Passing", in Forty Fifth Asilomar Conference on Signals, Systems and Computers (ASILOMAR), pp. 799-803, Pacific Grove, California, USA, Nov. 6-9, 2011
- [4] J. P. Vila and P. Schniter, "Expectation-Maximization Gaussian-Mixture Approximate Message Passing", IEEE Transactions on Signal Processing 2013, vol. 61, no. 19, pp. 4658-4672, Oct. 2013.

compute (var)

Compute the IIDBG input channel value.

Parameters: var (*dict*) – The variables used in computing of the input channel value.

Returns: • **mean** (*ndarray*) – The computed input channel mean.

• variance (ndarray) – The computed input channel variance.

class magni.cs.reconstruction.gamp.input_channel. IIDsGB(var)

Bases: magni.cs.reconstruction.gamp.input channel.ValidatedMMSEInputChannel

An i.i.d. sparse Gauss Bernoulli MMSE input channel.

This channel is based on equations (68), (69) in[1] and allows for using Expectation Maximization (EM) for learning the channel parameters as detailed in equations (74), (78), and (79) in [1] (see also [2] for an introduction to EM for GAMP).

- Parameters: tau (float or int) The prior signal "density" (fraction of Gauss to Bernouilli).
 - theta bar (float or int) The prior Gaussian mean.
 - theta_tilde (float or int) The prior Gaussian variance.
 - use_em (bool) The indicator of whether or not to use Expectation Maximazion (EM) to learn the prior parameters.
 - em_damping (float or int) The damping of the EM updates (if using EM).

Notes

The above Parameters are the input channel parameters that must be passed in a var dict to the channel constructor.

If use em is True, the values given for tau, theta bar, and theta tilde are used for initialialisation. The *em damping* must be in [0, 1) with 0 being no damping.

This channel is theoretically equivalent to the **IIDBG** channel. However, due to numerical subtleties, it may give different results.

References

[1] (1, 2) F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Probabilistic reconstruction in compressed sensing: algorithms, phase diagrams, and threshold achieving matrices", Journal of Statistical Mechanics: Theory and Experiment, vol. P08009, pp. 1-57, Aug. 2012.

11/2/16 Page 70 of 209 [2] J. P. Vila and P. Schniter, "Expectation-Maximization Gaussian-Mixture Approximate Message Passing", *IEEE Transactions on Signal Processing*, 2013, vol. 61, no. 19, pp. 4658-4672, Oct. 2013.

compute (var)

Compute the IIDsGB input channel value.

Parameters: var (*dict*) – The variables used in computing of the input channel

value.

Returns: • mean (ndarray) – The computed input channel mean.

• **variance** (*ndarray*) – The computed input channel variance.

magni.cs.reconstruction.gamp.output channel module

Module providing output channel functions for the Generalised Approximate Message Passing (GAMP) algorithm.

Routine listings

ValidatedMMSEOutputChannel(magni.utils.validation.types.MMSEOutputChannel)
A base class for validated magni.cs.reconstruction.gamp output channels.

AWGN(ValidatedMMSEOutputChannel)

An Additive White Gaussian Noise (AWGN) MMSE output channel.

class magni.cs.reconstruction.gamp.output_channel. ValidatedMMSEOutputChannel(var)

Bases: magni.utils.validation.types.MMSEOutputChannel

A base class for validated magni.cs.reconstruction.gamp output channels.

Parameters: var (*dict*) – The output channel state variables.

compute (var)

Compute the output channel value.

Parameters: var (*dict*) – The variables used in computing of the output channel

value.

Returns: • mean (ndarray) – The computed output channel mean.

• **variance** (*ndarray*) – The computed output channel variance.

Notes

This method honors magni.utils.validation.enable allow validate once.

class magni.cs.reconstruction.gamp.output_channel. AWGN (var)

Bases: magni.cs.reconstruction.gamp.output_channel.ValidatedMMSEOutputChannel

An Additive White Gaussian Noise (AWGN) MMSE output channel.

This channel is based on equations (41), (42), and (43) in[1] and allows for using Expectation Maximization (EM) for learning the channel parameter as detailed in equation (77) in [2] (see also [2] for an introduction to EM for GAMP).

11/2/16 Page 71 of 209

- Parameters: sigma sq (float or int) The noise level variance (initial noise level when the noise level is estimated).
 - noise level estimation (str) The method used for estimating (learning) the noise level in each iteration.

Notes

The above Parameters are the output channel parameters that must be passed in avar dict to the channel constructor.

Possible values for *noise level estimation* are:

- 'sample variance' Estimate noise level using the sample variance.
- 'median' Estimate noise level from the median.
- 'em' Estimate noise using Expectation Maximization (EM).
- 'fixed' Use a fixed noise level in all iterations.

References

- [1] S. Rangan, "Generalized Approximate Message Passing for Estimation with Random Linear Mixing", arXiv:1010.5141v2, pp. 1-22, Aug. 2012.
- [2] (1, 2) F. Krzakala, M. Mezard, F. Sausset, Y. Sun, and L. Zdeborova, "Probabilistic reconstruction in compressed sensing: algorithms, phase diagrams, and threshold achieving matrices", Journal of Statistical Mechanics: Theory and Experiment, vol. P08009, pp. 1-57, Aug. 2012.
- [3] J. P. Vila and P. Schniter, "Expectation-Maximization Gaussian-Mixture Approximate Message Passing", IEEE Transactions on Signal Processing, 2013, vol. 61, no. 19, pp. 4658-4672, Oct. 2013.

compute (var)

Compute the AWGN output channel value.

Parameters: var (*dict*) – The variables used in computing of the output channel value.

Returns: • **mean** (*ndarray*) – The computed output channel mean.

• **variance** (*ndarray*) – The computed output channel variance.

magni.cs.reconstruction.gamp.stop criterion module

Module providing stop criteria for the Generalised Approximate Message Passing, GAMP, algorithm.

Routine listings

ValidatedStopCriterion(magni.utils.validation.types.StopCriterion)

A base class for validated magni.cs.reconstruction.gamp stop criteria.

MSEConvergence(ValidatedStopCriterion)

A mean square error (MSE) convergence stop criterion.

NormalisedMSEConvergence(ValidatedStopCriterion)

A normalised mean squaure error (NMSE) convergence stop criterion.

11/2/16 Page 72 of 209 Residual(ValidatedStopCriterion)

A residual based stop criterion.

ResidualMeasurementsRatio(ValidatedStopCriterion)

A residual-measurements-ratio based stop criterion.

class magni.cs.reconstruction.gamp.stop_criterion. ValidatedStopCriterion(var)

Bases: magni.utils.validation.types.StopCriterion

A base class for validated magni.cs.reconstruction.gamp stop criteria.

Parameters: var (*dict*) – The stop criterion state variables.

compute(var)

Compute the stop criterion value.

Parameters: var (*dict*) – The variables used in computing of the stop criterion value.

Returns: • sto

 stop (bool) – The indicator of whether or not the stop criterion is satisfied.

• **value** (*float*) – The stop criterion value.

Notes

This method honors magni.utils.validation.enable allow validate once.

class magni.cs.reconstruction.gamp.stop_criterion. MSEConvergence(var)

Bases: magni.cs.reconstruction.gamp.stop_criterion.ValidatedStopCriterion

A mean square error (MSE) convergence stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

- A
- tolerance

Examples

For example,

11/2/16 Page 73 of 209

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp.stop_criterion import MSEConvergence
>>> state = {'tolerance': 1e-3, 'A': np.ones((10,10))}
>>> variables = {'alpha_bar_prev': np.ones((10, 1)),
... 'alpha_bar': np.arange(10).reshape(10, 1)}
>>> MSE = MSEConvergence(state)
>>> stop, val = MSE.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute(var)

Compute the MSE convergence stop criterion value.

The GAMP algorithm should converge to a fixed point. This criterion is based on the mean squared error of the difference between the proposed solution in this iteration and the proposed solution in the previous solution.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

Returns:

- stop (bool) The indicator of whether or not the stop criterion is satisfied.
- value (float) The stop criterion value.

class magni.cs.reconstruction.gamp.stop_criterion. **NormalisedMSEConvergence**(var)

Bases: magni.cs.reconstruction.gamp.stop_criterion.ValidatedStopCriterion

A normalised mean squaure error (NMSE) convergence stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

tolerance

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp.stop_criterion import (
... NormalisedMSEConvergence)
>>> state = {'tolerance': 1e-3}
>>> variables = {'alpha_bar_prev': np.ones((10, 1)),
... 'alpha_bar': np.arange(10).reshape(10, 1)}
>>> NMSE = NormalisedMSEConvergence(state)
>>> stop, val = NMSE.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute (var)

Compute the normalised MSE convergence stop criterion value.

11/2/16 Page 74 of 209

The GAMP algorithm should converge to a fixed point. This criterion is based on the mean squared error of the difference between the proposed solution in this iteration and the proposed solution in the previous solution normalised by the mean squared error of the proposed solution in the previous iteration.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

Returns: • sto

 stop (bool) – The indicator of whether or not the stop criterion is satisfied.

• value (float) – The stop criterion value.

class magni.cs.reconstruction.gamp.stop criterion. **Residual** (var)

Bases: magni.cs.reconstruction.gamp.stop_criterion.ValidatedStopCriterion

A residual based stop criterion.

Parameters: var (*dict*) – The stop criterion state

variables.

Notes

The following state variables are used in this stop criterion:

- y
- tolerance
- A

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp.stop_criterion import Residual
>>> state = {'tolerance': 1e-3, 'y': np.ones((10, 1)),
... 'A': np.ones((10,10))}
>>> variables = {'A_dot_alpha_bar': np.arange(10).reshape(10, 1)}
>>> Res = Residual(state)
>>> stop, val = Res.compute(variables)
>>> stop
False
>>> np.round(val, 2)
20.5
```

compute(var)

Compute the residual stop criterion value.

If the noise level is (approximately) known, the GAMP iterations may be stopped once the residual is on the order of the noise level. This stopping criterion is based on the mean sqaured error of the residual.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

11/2/16 Page 75 of 209

Returns:

- stop (bool) The indicator of whether or not the stop criterion is satisfied.
- value (float) The stop criterion value.

class magni.cs.reconstruction.gamp.stop_criterion. ResidualMeasurementsRatio (var)

Bases: magni.cs.reconstruction.gamp.stop_criterion.ValidatedStopCriterion

A residual-measurements-ratio based stop criterion.

Parameters: var (*dict*) – The stop criterion state variables.

Notes

The following state variables are used in this stop criterion:

- y
- tolerance

Examples

For example,

```
>>> import numpy as np
>>> from magni.cs.reconstruction.gamp.stop_criterion import (
... ResidualMeasurementsRatio)
>>> state = {'tolerance': 1e-3, 'y': np.ones((10, 1))}
>>> variables = {'A_dot_alpha_bar': np.arange(10).reshape(10, 1)}
>>> ResMeasRat = ResidualMeasurementsRatio(state)
>>> stop, val = ResMeasRat.compute(variables)
>>> stop
False
>>> np.round(val, 2)
14.32
```

compute (var)

Compute the residual-measurements-ratio stop criterion value.

If the noise level is (approximately) known, the GAMP iterations may be stopped once the residual is on the order of the noise level. This stopping criterion is based on ratio of the mean squared error of the residual to the mean squared error of the measurements.

Parameters: var (*dict*) – Dictionary of variables used in calculating of the stop criterion.

• **stop** (*bool*) – The indicator of whether or not the stop criterion is satisfied.

• value (float) - The stop criterion value.

magni.cs.reconstruction.iht package

Subpackage providing an implementation of Iterative Hard Thresholding (IHT).

Note: Deprecated in Magni 1.3.0. magni.cs.reconstruction.iht will be removed in a future version. Use the more general magni.cs.reconstruction.it instead.

11/2/16 Page 76 of 209

Warning: Change of variable interpretation. In magni.cs.reconstruction.it the config variable threshold fixed has a different interpretation than in magni.cs.reconstruction.iht.

Routine listings

config

Configger providing configuration options for this subpackage.

run(y, A)

Run the IHT reconstruction algorithm.

Notes

The IHT reconstruction algorithm is described in [1]. The default configuration uses the False Alarm Rate heuristic described in [2].

References

- [1] T. Blumensath and M.E. Davies, "Iterative Thresholding for Sparse Approximations", Journal of Fourier Analysis and Applications, vol. 14, pp. 629-654, Sep. 2008.
- [2] A. Maleki and D.L. Donoho, "Optimally Tuned Iterative Reconstruction Algorithms for Compressed Sensing", IEEE Journal Selected Topics in Signal Processing vol. 3, no. 2, pp. 330-341, Apr. 2010.

magni.cs.reconstruction.iht. **run** (y, A)

Run the IHT reconstruction algorithm.

Note: Deprecated in Magni 1.3.0 magni.cs.reconstruction.iht will be removed in a future version. Use the more general magni.cs.reconstruction.it instead.

- **Parameters: y** (*ndarray*) The m x 1 measurement vector.
 - A (ndarray) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.

Returns:

alpha (*ndarray*) – The n x 1 reconstructed coefficient vector.

Examples

For example, recovering a vector from random measurements

11/2/16 Page 77 of 209

```
>>> import warnings
>>> import numpy as np
>>> from magni.cs.reconstruction.iht import run, config
>>> np.random.seed(seed=6021)
>>> A = 1 / np.sqrt(80) * np.random.randn(80, 200)
>>> alpha = np.zeros((200, 1))
>>> alpha[:10] = 1
>>> y = A.dot(alpha)
>>> with warnings.catch_warnings():
    warnings.simplefilter('ignore')
     alpha_hat = run(y, A)
>>> np.set_printoptions(suppress=True)
>>> alpha_hat[:12]
array([[ 0.99836297],
    [ 1.00029086],
    [0.99760224],
    [0.99927175],
    [0.99899124],
    [0.99899434],
    [ 0.9987368 ],
    [ 0.99801849],
    [ 1.00059408],
    [ 0.9983772 ],
    [ 0.
    [ 0.
            ]])
>>> (np.abs(alpha_hat) > 1e-2).sum()
```

Or recovering the same only using a fixed threshold level:

```
>>> config.update({'threshold': 'oracle', 'threshold fixed': 10./80})
>>> with warnings.catch_warnings():
   warnings.simplefilter('ignore')
    alpha_hat_2 = run(y, A)
>>> alpha_hat_2[:12]
array([[ 0.99877706],
    [0.99931441],
     0.9978366 1.
    [0.99944973],
    [1.00052762],
    [1.00033436],
    [0.99943286],
    [0.99952526],
    [0.99941578],
    [0.99942908],
    [ 0.
    [ 0.
            ]])
>>> (np.abs(alpha_hat_2) > 1e-2).sum()
10
```

Submodules

magni.cs.reconstruction.iht._config module

Module providing configuration options for the magni.cs.reconstruction.iht subpackage.

```
See also:

magni.cs.reconstruction._config.Configger

The Configger class used
```

11/2/16 Page 78 of 209

Notes

This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

iterations

: int

The maximum number of iterations to do (the default is 300).

kappa fixed

: float

The relaxation parameter used in the algorithm (the default is 0.65).

precision float

: {np.float, np.float16, np.float32, np.float64, np.float128}

The floating point precision used for the computations (the default is np.float64).

threshold

: {'far', 'oracle'}

The method for selecting the threshold value.

threshold fixed

: float

The assumed rho value used for selecting the threshold value if using the oracle method.

tolerance

: float

The least acceptable ratio of residual to measurements (in 2-norm) to break the interations (the default is 0.001).

magni.cs.reconstruction.it package

Subpackage providing implementations of Iterative Thresholding (IT).

Routine listings

config

Configger providing configuration options for this subpackage.

run(y, A)

Run the IT reconstruction algorithm.

Notes

Implementations of Iterative Hard Thresholding (IHT)[1], [2] as well as implementations of Iterative Soft Thresholding (IST) [3], [4] are available. It is also possible to configure the subpackage to use a model based approach as described in [5].

References

- [1] T. Blumensath and M.E. Davies, "Iterative Thresholding for Sparse Approximations", Journal of Fourier Analysis and Applications, vol. 14, pp. 629-654, Sep. 2008.
- [2] T. Blumensath and M.E. Davies, "Normalized Iterative Hard Thresholding: Guaranteed Stability and Performance", *IEEE Journal Selected Topics in Signal Processing* vol. 4, no. 2, pp. 298-309, Apr. 2010.

11/2/16 Page 79 of 209

- [3] I. Daubechies, M. Defrise, and C. D. Mol, "An Iterative Thresholding Algorithm for Linear Inverse Problems with a Sparsity Constraint", *Communications on Pure and Applied Mathematics*, vol. 57, no. 11, pp. 1413-1457, Nov. 2004.
- [4] A. Maleki and D.L. Donoho, "Optimally Tuned Iterative Reconstruction Algorithms for Compressed Sensing", *IEEE Journal Selected Topics in Signal Processing* vol. 3, no. 2, pp. 330-341, Apr. 2010.
- [5] R.G. Baraniuk, V. Cevher, M.F. Duarte, and C. Hedge, "Model-Based Compressive Sensing", *IEEE Transactions on Information Theory*, vol. 56, no. 4, pp. 1982-2001, Apr. 2010.

Submodules

magni.cs.reconstruction.it. algorithm module

Module providing the core Iterative Thresholding (IT) algorithm.

Routine listings

run(y, A)

Run the IT reconstruction algorithm.

See also:

magni.cs.reconstruction.it._config

Configuration options.

magni.cs.reconstruction.it._step_size

Available step-size strategies.

magni.cs.reconstruction.it._stop_criterion

Available stop criteria.

magni.cs.reconstruction.it._threshold

Available threshold estimators.

magni.cs.reconstruction.it._threshold_operators

Available threshold operators

Notes

The default configuration of the IT algorithm provides the Iterative Hard Thresholding (IHT) algorithm [1] using the False Alarm Rate (FAR) heuristic from [2]. By default a residual measurements ratio step criterion is used. The algorithm may also be configured to act as the Iterative Soft Thresholding (IST) [3] algorithm with the soft threshold from [4].

References

- [1] T. Blumensath and M.E. Davies, "Iterative Thresholding for Sparse Approximations", Journal of Fourier Analysis and Applications, vol. 14, pp. 629-654, Sep. 2008.
- [2] A. Maleki and D.L. Donoho, "Optimally Tuned Iterative Reconstruction Algorithms for Compressed Sensing", *IEEE Journal Selected Topics in Signal Processing* vol. 3, no. 2, pp. 330-341, Apr. 2010.
- [3] I. Daubechies, M. Defrise, and C. D. Mol, "An Iterative Thresholding Algorithm for Linear Inverse Problems with a Sparsity Constraint", *Communications on Pure and Applied Mathematics*, vol. 57, no. 11, pp. 1413-1457, Nov. 2004.
- [4] D.L. Donoho, "De-Noising by Soft-Thresholding", *IEEE Transactions on Information Theory*, vol. 41, no. 3, pp. 613-627, May. 1995.

11/2/16 Page 80 of 209

magni.cs.reconstruction.it. algorithm. run(y, A)

Run the IT reconstruction algorithm.

- **Parameters:** y (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix. MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.

Returns:

- **alpha** (*ndarray*) The n x 1 reconstructed coefficient vector.
- **history** (*dict*, *optional*) The dictionary of various measures tracked in the IT iterations.

See also:

magni.cs.reconstruction.it._config()

Configuration options.

magni.cs.reconstruction.it._step_size()

Step-size calculation.

magni.cs.reconstruction.it._stop_criterion()

Stop criterion calculation.

magni.cs.reconstruction.it._threshold()

Threshold calculation.

magni.cs.reconstruction.it. threshold operators()

Thresholding operators.

Notes

Optionally, the algorithm may be configured to save and return the iteration history along with the reconstruction result. The returned history contains the following:

- alpha : Reconstruction coefficient estimates
- MSE: Solution mean squared error (if the true solution is known).
- stop criterion: The currently used stop criterion.
- stop criterion value: The value of the stop criterion.
- stop iteration: The iteration at which the algorithm stopped.
- stop_reason : The reason for termination of the algorithm.

Examples

For example, recovering a vector from random measurements using IHT with the FAR heuristic

11/2/16 Page 81 of 209

```
>>> import numpy as np, magni
>>> from magni.cs.reconstruction.it import run
>>> np.random.seed(seed=6021)
>>> A = 1 / np.sqrt(90) * np.random.randn(90, 200)
>>> alpha = np.zeros((200, 1))
>>> alpha[:10] = 1
>>> y = A.dot(alpha)
>>> alpha_hat = run(y, A)
>>> alpha_hat[:12]
array([[ 0.99748945],
    [ 1.00082074],
    [ 0.99726507],
    [0.99987834],
    [ 1.00025857],
    [ 1.00003266],
    [ 1.00021666],
    [0.99838884],
    [ 1.00018356],
    [0.99859105],
    [ 0.
            ]])
    [ 0.
>>> (np.abs(alpha_hat) > 1e-2).sum()
```

Or recover the same vector as above using IST with a fixed number of non-thresholded coefficients

```
>>> it_config = {'threshold_operator': 'soft', 'threshold': 'fixed',
... 'threshold_fixed': 10}
>>> magni.cs.reconstruction.it.config.update(it_config)
>>> alpha hat = run(y, A)
>>> alpha_hat[:12]
array([[ 0.99822443],
     [0.99888724],
     [ 0.9982493 ],
     [ 0.99928642],
    [0.99964131],
    [ 0.99947346],
    [ 0.9992809 ],
    [ 0.99872624],
    [0.99916842],
    [0.99903033],
    [0.
    [ 0.
            ]])
>>> (np.abs(alpha_hat) > 1e-2).sum()
```

Or use a weighted IHT method with a fixed number of non-thresholded coefficietns to recover the above signal

11/2/16 Page 82 of 209

```
>>> weights = np.linspace(1.0, 0.9, 200).reshape(200, 1)
>>> it_config = {'threshold_operator': 'weighted_hard',
... 'threshold weights': weights}
>>> magni.cs.reconstruction.it.config.update(it config)
>>> alpha_hat = run(y, A)
>>> alpha_hat[:12]
array([[ 0.99853888],
    [0.99886104],
    [0.99843149],
    [ 1.0000333 ],
    [ 1.00060273],
    [ 1.00035539],
    [0.99966219],
    [0.99912209],
    [0.99961541],
    [0.99952029],
    [ 0.
>>> (np.abs(alpha_hat) > 1e-2).sum()
>>> magni.cs.reconstruction.it.config.reset()
```

magni.cs.reconstruction.it._config module

Module providing configuration options for the magni.cs.reconstruction.it subpackage.

```
See also:

magni.cs.reconstruction._config.Configger

The Configger class used
```

Notes

This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

iterations

: int

The maximum number of iterations to do (the default is 300).

```
kappa
```

```
: {'fixed', 'adaptive'}
```

The method used to calculate the step-size (relaxation) parameter kappa.

```
kappa_fixed
```

: float

The step-size (relaxation parameter) used in the algorithm when a fixed step-size is used (the default is 0.65).

```
precision_float
```

```
: {np.float, np.float16, np.float32, np.float64, np.float128}
```

The floating point precision used for the computations (the default is np.float64).

```
report_history
```

: bool

The indicator of whether or not to return the progress history along with the result (the default is False).

stop_criterion

11/2/16 Page 83 of 209

: str

The stop criterion to use in the iterations (the default is residual_measurements_ratio).

threshold

: {'far', 'fixed'}

The method used for calculating the threshold value.

threshold fixed

: int

The number of non-zero coefficients in the signal vector when this number is assumed fixed

threshold operator

: {'hard', 'soft', 'weighted_hard', 'weighted_soft', 'none'}

The threshold operator used in the backprojection step.

threshold weights

: ndarray

Array of weights to be used in one of the weighted threshold operators (the default is array([[1]]), which implies that all coefficients are weighted equally).

tolerance

: float

The least acceptable ratio of residual to measurements (in 2-norm) to break the interations (the default is 0.001).

true solution

: ndarray or None

The true solution to allow for tracking the convergence of the algorithm in the artificial setup where the true solution is known a-priori (the default is None, which implies that no true solution tracking is used).

warm start

: ndarray

The initial guess of the solution vector (the default is None, which implies that a vector of zeros is used).

magni.cs.reconstruction.it. step size module

Module providing functions for calculating the step-size (relaxation parameter) used in the Iterative Thresholding algorithms.

Routine listings

```
calculate using adaptive(var)
```

Calculate a step-size value in an 'adaptive' way.

calculate using fixed(var)

Calculate the fixed step-size value.

get_function_handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.it._step_size. wrap_calculate_using_adaptive(var)

Arguments wrapper for culculate_using_adaptive.

Calculate a step-size value in an 'adaptive' way.

11/2/16 Page 84 of 209

magni.cs.reconstruction.it. step size. wrap calculate using fixed(var)

Arguments wrapper for calculate using fixed.

Calculate a fixed step-size value.

magni.cs.reconstruction.it._step_size. **get_function_handle**(*method*, *var*)

Return a function handle to a given calculation method.

Parameters: • method (str) – Identifier of the calculation method to return a handle

• var (dict) – Local variables needed in the step-size method.

f handle (function) - Handle to calculation method defined in this Returns:

globals scope.

magni.cs.reconstruction.it. stop criterion module

Module providing functions for calculating stop criteria used in Iterative Threholding (IT) algorithms.

Routine listings

calculate using mse convergence(var)

Calculate stop criterion based on mse convergence.

calculate using residual(var)

Calculate stop criterion based on residual.

calculate using residual measurments ratio(var)

Calculate stop criterion based on the residual to measurements ratio.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.it._stop_criterion. wrap_calculate_using_mse_convergence(var)

Arguments wrapper for *calculate_using_mse_convergence*.

Calculate stop criterion based on mse convergence.

magni.cs.reconstruction.it._stop_criterion.wrap_calculate_using_residual(var)

Arguments wrapper for calculate using residual.

Calculate stop criterion based on residual.

magni.cs.reconstruction.it. stop criterion. wrap calculate using residual measurements ratio (var)

Arguments wrapper for calculate using residual measurements ratio.

Calculate stop criterion based on the residual to measurements ratio.

magni.cs.reconstruction.it._stop_criterion. **get_function_handle**(*method*, *var*)

Return a function handle to a given calculation method.

Parameters: • method (str) – Identifier of the calculation method to return a handle to.

• **var** (*dict*) – Local variables needed in the calculation method.

11/2/16 Page 85 of 209

f handle (function) - Handle to calculation method defined in this Returns: globals scope.

magni.cs.reconstruction.it. threshold module

Module providing functions for calculating a threshold (level) used in Iterative Threshold algorithms.

Routine listings

calculate far(delta)

Calculate the optimal False Acceptance Rate for a given indeterminacy.

calculate using far(var)

Calculate a threshold level using the FAR heuristic.

calculate using fixed(var)

Calculate a threshold level using a given fixed support size.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.it._threshold. **calculate_far**(*delta*, *it_algorithm*)

Calculate the optimal False Acceptance Rate for a given indeterminacy.

- **Parameters: delta** (*float*) The indeterminacy, m / n, of a system of equations of size m x n.
 - it_algorithm ({IHT, ITS}) The iterative thresholding algorithm to calculate the FAR for.

Returns:

FAR (float) - The optimal False Acceptance Rate for the given indeterminacy.

Notes

The optimal False Acceptance Rate to be used in connection with the interference heuristic presented in the paper "Optimally Tuned Iterative Reconstruction Algorithms for Compressed Sensing" [1] is calculated from a set of optimal values presented in the same paper. The calculated value is found from a linear interpolation or extrapolation on the known set of optimal values.

References

[1] A. Maleki and D.L. Donoho, "Optimally Tuned Iterative Reconstruction Algorithms for Compressed Sensing", IEEE Journal Selected Topics in Signal Processing vol. 3, no. 2, pp. 330-341, Apr. 2010.

magni.cs.reconstruction.it._threshold. wrap_calculate_using_far(var)

Arguments wrapper for calculate using far.

Calculate a threshold level using the FAR heuristic.

magni.cs.reconstruction.it. threshold. wrap calculate using fixed(var)

Arguments wrapper for calculate using fixed.

11/2/16 Page 86 of 209 Calculate a threshold level using a given fixed support size.

magni.cs.reconstruction.it._threshold. **get_function_handle**(*method*, *var*)

Return a function handle to a given calculation method.

• method (str) – Identifier of the calculation method to return a handle

• var (dict) - Local variables needed in the threshold method.

Returns: f_handle (function) - Handle to calculation method defined in this

globals scope.

magni.cs.reconstruction.it. threshold operators module

Module providing thresholding operators used in Iterative Thresholding algorithms.

Routine listings

get function handle(method)

Return a function handle to a given threshold operator.

threshold_hard(var)

The hard threshold operator.

threshold_none(var)

The "no" threshold operator.

threshold_soft(var)

The soft threshold operator.

threshold_weighted_hard(var)

The weighted hard threshold operator.

threshold weighted soft(var)

The weighted soft threshold operator.

magni.cs.reconstruction.it._threshold_operators. **get_function_handle**(*method*)

Return a function handle to a given threshold operator method.

Parameters: method (*str*) – Identifier of the threshold operator to return a handle to.

Returns: f handle (*function*) – Handle to threshold method defined in this globals

scope.

magni.cs.reconstruction.it. threshold operators. **threshold hard**(*var*)

Threshold the entries of a vector using the hard threshold.

Parameters: var (*dict*) – Local variables used in the threshold

operation.

Notes

This threshold operation works "in-line" on the variables in *var*. Hence, this function does not return anything.

Examples

11/2/16 Page 87 of 209

For example, thresholding a vector of values between -1 and 1

```
>>> import copy, numpy as np, magni
>>> from magni.cs.reconstruction.it._threshold_operators import (
... threshold_hard)
>>> var = {'alpha': np.linspace(-1, 1, 10), 'threshold': 0.4}
>>> threshold_hard(copy.copy(var))
>>> var['alpha']
array([-1. , -0.77777778, -0.55555556, 0. , 0. ,
0. , 0. , 0.555555556, 0.77777778, 1. ])
```

magni.cs.reconstruction.it._threshold_operators. threshold_none(var)

Do not threshold the entries of a vector.

```
Parameters: var (dict) – Local variables used in the threshold operation.
```

Notes

This is a dummy threshold operation that does nothing.

```
magni.cs.reconstruction.it._threshold_operators. threshold_soft(var)
```

Threshold the entries of a vector using the soft threshold.

```
Parameters: var (dict) – Local variables used in the threshold operation.
```

Notes

This threshold operation works "in-line" on the variables in *var*. Hence, this function does not return anything.

Examples

For example, thresholding a vector of values between -1 and 1

```
>>> import copy, numpy as np, magni
>>> from magni.cs.reconstruction.it._threshold_operators import (
... threshold_soft)
>>> var = {'alpha': np.linspace(-1, 1, 10), 'threshold': 0.4}
>>> threshold_soft(copy.copy(var))
>>> var['alpha']
array([-0.6 , -0.37777778, -0.15555556, 0. , 0. ,
0. , 0. , 0.15555556, 0.37777778, 0.6 ])
```

magni.cs.reconstruction.it._threshold_operators. threshold_weighted_hard (var)

Threshold the entries of a vector using a weighted hard threshold.

```
Parameters: var (dict) – Local variables used in the threshold operation.
```

Notes

This threshold operation works "in-line" on the variables in *var*. Hence, this function does not return anything.

11/2/16 Page 88 of 209

Examples

For example, thresholding a vector of values between -1 and 1

```
>>> import copy, numpy as np, magni
>>> from magni.cs.reconstruction.it._threshold_operators import (
... threshold_weighted_hard)
>>> var = {'alpha': np.linspace(-1, 1, 10), 'threshold': 0.4,
... 'threshold_weights': 0.7 * np.ones(10)}
>>> threshold_weighted_hard(copy.copy(var))
>>> var['alpha']
array([-1. , -0.777777778, 0. , 0. , 0. ,
0. , 0. , 0. , 0. , 0.77777778, 1. ])
```

magni.cs.reconstruction.it._threshold_operators.threshold_weighted_soft(var)

Threshold the entries of a vector using a weighted soft threshold.

Parameters: var (*dict*) – Local variables used in the threshold operation.

Notes

This threshold operation works "in-line" on the variables in *var*. Hence, this function does not return anything.

Examples

For example, thresholding a vector of values between -1 and 1

```
>>> import copy, numpy as np, magni
>>> from magni.cs.reconstruction.it._threshold_operators import (
... threshold_weighted_soft)
>>> var = {'alpha': np.linspace(-1, 1, 10), 'threshold': 0.4,
... 'threshold_weights': 0.7 * np.ones(10)}
>>> threshold_weighted_soft(copy.copy(var))
>>> var['alpha']
array([-0.42857143, -0.20634921, 0. , 0. , 0. ,
0. , 0. , 0. , 0.20634921, 0.42857143])
```

magni.cs.reconstruction.it. util module

Module providing miscellaneous utility functions for the different implementations of Iterative Thresholding (IT).

Routine listings

```
get_methods(module)
```

Extract relevant methods from module.

```
def get operators(module)
```

Extract relevant operators from module.

```
magni.cs.reconstruction.it. util. get methods(module)
```

Extract relevant methods from module.

Parameters: module (*str*) – The name of the module from which the methods should be extracted.

11/2/16 Page 89 of 209

Returns: methods (tuple) – The names of the relevant methods from module.

Notes

Looks for functions with names starting with 'wrap_calculate_using_' in the module: "module".

magni.cs.reconstruction.it._util. _get_operators(module)

Extract relevant operators from *module*.

Parameters: module (*str*) – The name of the module from which the operators

should be extracted.

Returns: methods (*tuple*) – The names of the relevant operators from *module*.

Notes

Looks for functions with names starting with 'threshold' in the module "module".

magni.cs.reconstruction.sl0 package

Subpackage providing implementations of Smoothed IO Norm (SL0).

The implementations provided are the original SL0 reconstruction algorithm and a modified SL0 reconstruction algorithm. The algorithm used depends on the 'algorithm' configuration option: 'std' refers to the original algorithm while 'mod' (default) refers to the modified algorithm.

Routine listings

config

Configger providing configuration options for this subpackage.

run(y, A)

Run the specified SL0 reconstruction algorithm.

Notes

See <u>algorithm</u> for documentation of run.

Implementations of the original SL0 reconstruction algorithm[1] and a modified Sl0 reconstruction algorithm [3] are available. It is also possible to configure the subpackage to provide customised versions of the SL0 reconstruction algorithm. The projection algorithm [1] is used for small delta (< 0.55) whereas the contraint elimination algorithm[2] is used for large delta (>= 0.55) which merely affects the computation time.

References

- [1] (1, 2) H. Mohimani, M. Babaie-Zadeh, and C. Jutten, "A Fast Approach for Overcomplete Sparse Decomposition Based on Smoothed IO Norm", *IEEE Transactions on Signal Processing*, vol. 57, no. 1, pp. 289-301, Jan. 2009.
- [2] Z. Cui, H. Zhang, and W. Lu, "An Improved Smoothed I0-norm Algorithm Based on Multiparameter Approximation Function", in 12th IEEE International Conference on Communication Technology (ICCT), Nanjing, China, Nov. 11-14, 2011, pp. 942-945.

11/2/16 Page 90 of 209

[3] C. S. Oxvig, P. S. Pedersen, T. Arildsen, and T. Larsen, "Surpassing the Theoretical 1-norm Phase Transition in Compressive Sensing by Tuning the Smoothed IO Algorithm", in IEEE International Conference on Acoustics, Speech and Signal Processing (ICASSP), Vancouver, Canada, May 26-31, 2013, pp. 6019-6023.

Submodules

magni.cs.reconstruction.sl0. L start module

Module providing functions for calculating the starting value of L (the number of gradient descent iterations for each sigma) used in the Smoothed IO algorithms.

Routine listings

calculate using fixed(var)

Calculate the fixed L value.

calculate using geometric(var)

Calculate an L value in a 'geometric' way.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.sl0._L_start. wrap_calculate_using_fixed(var)

Arguments wrapper for calculate using fixed.

magni.cs.reconstruction.sl0._L_start.wrap_calculate_using_geometric(var)

Arguments wrapper for calculate using geometric.

magni.cs.reconstruction.sl0._L_start.get_function_handle(method, var)

Return a function handle to a given calculation method.

- **Parameters:** method (str) Identifier of the calculation method to return a handle
 - **var** (*dict*) Local variables needed in the L method.

Returns:

f handle (function) - Handle to the calculation method defined in this globals scope.

magni.cs.reconstruction.sl0._L_update module

Module providing functions for calculating the updated value of L (the number of gradient descent iterations for each sigma) used in the Smoothed IO algorithms.

Routine listings

calculate using fixed(var)

Calculate the updated, fixed L value.

calculate using geometric(var)

Calculate an updated L value in a 'geometric' way.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.sl0._L_update. wrap_calculate_using_fixed(var)

Arguments wrapper for calculate using fixed.

magni.cs.reconstruction.sl0. L update. wrap calculate using geometric(var)

11/2/16 Page 91 of 209 Arguments wrapper for calculate using geometric.

magni.cs.reconstruction.sl0. L update. **get function handle** (method, var)

Return a function handle to a given calculation method.

Parameters: • method (str) – Identifier of the calculation method to return a handle

• var (dict) – Local variables needed in the L update method.

Returns: f_handle (*function*) – Handle to the calculation method defined in this

globals scope.

magni.cs.reconstruction.sl0. algorithm module

Module providing the core Smoothed I0 (SL0) algorithm.

Routine listings

run(y, A)

Run the SL0 reconstruction algoritm.

See also:

magni.cs.reconstruction.sl0._config
Configuration options.

Notes

Implementations of the original SL0 reconstruction algorithm[1] and a modified Sl0 reconstruction algorithm [3] are available. It is also possible to configure the subpackage to provide customised versions of the SL0 reconstruction algorithm. The projection algorithm [1] is used for small delta (< 0.55) whereas the contraint elimination algorithm[2] is used for large delta (>= 0.55) which merely affects the computation time.

References

- [1] (1, 2) H. Mohimani, M. Babaie-Zadeh, and C. Jutten, "A Fast Approach for Overcomplete Sparse Decomposition Based on Smoothed IO Norm", *IEEE Transactions on Signal Processing*, vol. 57, no. 1, pp. 289-301, Jan. 2009.
- [2] Z. Cui, H. Zhang, and W. Lu, "An Improved Smoothed I0-norm Algorithm Based on Multiparameter Approximation Function", *in 12th IEEE International Conference on Communication Technology (ICCT)*, Nanjing, China, Nov. 11-14, 2011, pp. 942-945.
- [3] C. S. Oxvig, P. S. Pedersen, T. Arildsen, and T. Larsen, "Surpassing the Theoretical 1-norm Phase Transition in Compressive Sensing by Tuning the Smoothed IO Algorithm", in IEEE International Conference on Acoustics, Speech and Signal Processing (ICASSP), Vancouver, Canada, May 26-31, 2013, pp. 6019-6023.

magni.cs.reconstruction.sl0._algorithm. run(y, A)

Run the SL0 reconstruction algorithm.

11/2/16 Page 92 of 209

- **Parameters:** y (*ndarray*) The m x 1 measurement vector.
 - A (ndarray or magni.utils.matrices.{Matrix, MatrixCollection}) The m x n matrix which is the product of the measurement matrix and the dictionary matrix.

Returns:

alpha (*ndarray*) – The n x 1 reconstructed coefficient vector.

See also:

```
magni.cs.reconstruction.sl0.config()
    Configuration options.
```

Notes

The algorithm terminates after a fixed number of iterations or if the ratio between the 2norm of the residual and the 2-norm of the measurements falls below the specified tolerance.

Examples

For example, recovering a vector from random measurements using the original SL0 reconstruction algorithm

```
>>> import numpy as np. magni
>>> from magni.cs.reconstruction.sl0 import run
>>> np.set printoptions(suppress=True)
>>> magni.cs.reconstruction.sl0.config['L'] = 'fixed'
>>> magni.cs.reconstruction.sl0.config['mu'] = 'fixed'
>>> magni.cs.reconstruction.sl0.config['sigma_start'] = 'fixed'
>>> np.random.seed(seed=6021)
\rightarrow > A = 1 / np.sqrt(80) * np.random.randn(80, 200)
>>> alpha = np.zeros((200, 1))
>>> alpha[:10] = 1
>>> y = A.dot(alpha)
>>> alpha hat = run(y, A)
>>> alpha hat[:12]
array([[ 0.99993202],
    [0.99992793],
     [ 0.99998107],
    [0.99998105],
    [1.00005882],
    [ 1.00000843],
    [0.99999138],
    [ 1.00009479],
    [0.99995889].
    [0.99992509],
    [-0.00001509].
    [ 0.00000275]])
>>> (np.abs(alpha hat) > 1e-2).sum()
```

Or recover the same vector as above using the modified SL0 reconstruction algorithm

11/2/16 Page 93 of 209

```
>>> magni.cs.reconstruction.sl0.config['L'] = 'geometric'
>>> magni.cs.reconstruction.sl0.config['mu'] = 'step'
>>> magni.cs.reconstruction.sl0.config['sigma_start'] = 'reciprocal'
>>> alpha hat = run(y, A)
>>> alpha_hat[:12]
array([[ 0.9999963 ],
    [1.00000119],
    [ 1.00000293],
    [0.99999661],
    [1.00000021],
    [ 0.999951 ],
    [ 1.00000103],
    [ 1.00000662],
    [1.00000404],
    [0.99998937],
    [-0.00000075],
    [ 0.0000037]])
>>> (np.abs(alpha_hat) > 1e-2).sum()
```

magni.cs.reconstruction.sl0._config module

Module providing configuration options for the magni.cs.reconstruction.sl0 subpackage.

```
See also:

magni.cs.reconstruction._config.Configger

The Configger class used
```

Notes

This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

```
epsilon
```

: float

The precision parameter used in centering (the default is 0.01).

L

: {'geometric', 'fixed'}

The method for selecting the maximum number of gradient descent iterations for each sigma.

L_fixed

: int

The value used for L if using the 'fixed' method (the default is 2.0).

L_geometric_ratio

: float

The common ratio used for the L update if using the 'geometric' method (the default is 2.0).

L_geometric_start

: float

The starting value used for the L if using the 'geometric' method (the default is 2.0).

mu

: {'step', 'fixed'}

The method for selecting the relative step-size used in gradient descent iteration.

11/2/16 Page 94 of 209

```
mu_fixed
```

: float

The value used for mu if using the 'fixed' method (the default is 1.0).

```
mu step end
```

: float

The value used for mu for the last iterations if using the 'step' method (the default is 1.5).

mu_step_iteration

: int

The iteration where the value used for mu changes if using the 'step' method

mu step start

: float

The value used for mu for the first iterations if using the 'step' method (the default is 0.001).

precision float

: {np.float, np.float16, np.float32, np.float64, np.float128}

The floating point precision used for the computations (the default is np.float64).

sigma_geometric

: float

The common ratio used for the sigma update (the default is 0.7).

sigma start

: {'reciprocal', 'fixed'}

The method for selecting the starting value of sigma.

sigma start fixed

: float

The fixed factor multiplied onto the maximum coefficient of a least squares solution to obtain the value if using the 'fixed' method (the default is 2.0).

sigma start reciprocal

: float

The constant in the factor

 $\frac{1}{constant \ \delta}$ multiplied onto the maximum coefficient of a least squares solution to obtain the value if using the 'reciprocal' method (the default is 2.75).

sigma stop fixed

: float

The minimum value of std. dev. in Gaussian IO approx (the default 0.01).

magni.cs.reconstruction.sl0._mu_start module

Module providing functions for calculating the starting value of mu (the relative step-size used in the gradient descent iteration) used in the Smoothed IO algorithms.

Routine listings

```
calculate using fixed(var)
```

Calculate the fixed mu value.

calculate using step(var)

Calculate an mu value in a 'step' way.

```
get_function_handle(method)
```

11/2/16 Page 95 of 209

Return a function handle to a given calculation method.

magni.cs.reconstruction.sl0. mu start.wrap calculate using fixed(var)

Arguments wrapper for calculate using fixed.

magni.cs.reconstruction.sl0._mu_start.wrap_calculate_using_step(var)

Arguments wrapper for calculate_using_step.

magni.cs.reconstruction.sl0._mu_start.get_function_handle(method, var)

Return a function handle to a given calculation method.

• method (str) – Identifier of the calculation method to return a handle

to.

• var (dict) – Local variables needed in the mu method.

Returns: f_handle (function) – Handle to the calculation method defined in this

globals scope.

magni.cs.reconstruction.sl0._mu_update module

Module providing functions for calculating the updated value of mu (the relative step-size used in the gradient descent iteration) used in the Smoothed IO algorithms.

Routine listings

calculate_using_fixed(var)

Calculate the updated, fixed mu value.

calculate_using_step(var)

Calculate an updated mu value in a 'step' way.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.sl0._mu_update. wrap_calculate_using_fixed(var)

Arguments wrapper for calculate using fixed.

magni.cs.reconstruction.sl0._mu_update. wrap_calculate_using_step(var)

Arguments wrapper for calculate using step.

magni.cs.reconstruction.sl0._mu_update. **get_function_handle**(*method*, *var*)

Return a function handle to a given calculation method.

• method (str) – Identifier of the calculation method to return a handle

to.

• **var** (*dict*) – Local variables needed in the mu update method.

Returns: f_handle (function) – Handle to the calculation method defined in this

globals scope.

magni.cs.reconstruction.sl0._sigma_start module

Module providing functions for calculating the starting value of sigma used in the Smoothed I0 algorithms.

Routine listings

11/2/16 Page 96 of 209

calculate using fixed(var)

Calculate the fixed sigma value.

calculate_using_reciprocal(var)

Calculate a sigma value in a 'reciprocal' way.

get function handle(method)

Return a function handle to a given calculation method.

magni.cs.reconstruction.sl0. sigma start. wrap calculate using fixed(var)

Arguments wrapper for calculate using fixed.

magni.cs.reconstruction.sl0._sigma_start. wrap_calculate_using_reciprocal(var)

Arguments wrapper for calculate using reciprocal.

magni.cs.reconstruction.sl0._sigma_start. **get_function_handle**(*method*, *var*)

Return a function handle to a given calculation method.

Parameters: • method (str) – Identifier of the calculation method to return a handle

• **var** (*dict*) – Local variables needed in the sigma method.

Returns:

f handle (function) - Handle to the calculation method defined in this globals scope.

Submodules

magni.cs.reconstruction. config module

Module providing a CS reconstruction algorithm adapted configger subclass.

Routine listings

Configger(magni.utils.config.Configger)

Provide functionality to access a set of configuration options.

Notes

This module does not itself contain any configuration options and thus has no access to any configuration options unlike the other config modules of magni.

class magni.cs.reconstruction._config. Configger (params, valids)

Bases: magni.utils.config.Configger

Provide functionality to access a set of configuration options.

The present class redefines the methods for retrieving configuration parameters in order to ensure the desired precision of the floating point parameter values.

- **Parameters:** params (*dict*) The configuration options and their default values.
 - valids (dict) The validation schemes of the configuration options.

property

See also:

Page 97 of 209 11/2/16

magni.utils.config.Configger

Superclass of the present class.

__getitem___(name)

Get the value of a configuration parameter.

Parameters: name (*str*) – The name of the parameter.

Returns: value (*None*) – The value of the

parameter.

Notes

If the value is a floating point value then that value is typecast to the desired precision.

Submodules

magni.cs.indicators module

Module prodiving performance indicator determination functionality.

Routine listings

calculate_coherence(Phi, Psi, norm=None)

Calculate the coherence of the Phi Psi matrix product.

calculate mutual coherence(Phi, Psi, norm=None)

Calculate the mutual coherence of Phi and Psi.

calculate relative energy(Phi, Psi, method=None)

Calculate the relative energy of Phi Psi matrix product atoms.

Notes

For examples of uses of the performance indicators, see the related papers on predicting/modelling reconstruction quality [1], [2].

References

- [1] P. S. Pedersen, J. Ostergaard, and T. Larsen, "Predicting Reconstruction Quality within Compressive Sensing for Atomic Force Microscopy," *2015 IEEE Global Conference on Signal and Information Processing (GlobalSIP)*, Orlando, FL, 2015, pp. 418-422. doi: 10.1109/GlobalSIP.2015.7418229
- [2] P. S. Pedersen, J. Ostergaard, and T. Larsen, "Modelling Reconstruction Quality of Lissajous Undersampled Atomic Force Microscopy Images," 2016 IEEE 13th International Symposium on Biomedical Imaging (ISBI), Prague, Czech Republic, 2016, pp. 245-248. doi: 10.1109/ISBI.2016.7493255

magni.cs.indicators. **calculate_coherence**(*Phi*, *Psi*, *norm=None*)

Calculate the coherence of the Phi Psi matrix product.

In the context of Compressive Sensing, coherence usually refers to the maximum

11/2/16 Page 98 of 209

absolute correlation between two columns of the Phi Psi matrix product. This function allows the usage of a different normalised norm where the infinity-norm yields the usual case.

Parameters: • Phi

- Phi (magni.utils.matrices.Matrix or numpy.ndarray) The measurement matrix.
- **Psi** (*magni.utils.matrices.Matrix or numpy.ndarray*) The dictionary matrix.
- **norm** (*int or float*) The normalised norm used for the calculation (the default value is None which implies that the 0-, 1-, 2-, and infinity-norms are returned).

Returns:

coherence (*float or dict*) – The coherence value(s).

Notes

If *norm* is None, the function returns a dict containing the coherence using the 0-, 1-, 2-, and infinity-norms. Otherwise, the function returns the coherence using the specified norm.

The coherence is calculated as:

$$rac{1}{n^2-n} \sum_{i=1}^n \sum_{\substack{j=1 \ j=i}}^n rac{| \ \ _i \ \ \ _j|}{|| \ \ \ _i||_2|| \ \ \ \ _j||_2}$$

where n is the number of columns in Psi. In the case of the 0-norm, the coherence is calculated as:

$$\frac{1}{n^2 - n} \sum_{i=1}^{n} \sum_{\substack{j=1 \ i \ i}}^{n} \frac{|i \quad j|}{|| \quad i||_2 || \quad j||_2}$$

where

() is 1 if a is non-zero and 0 otherwise. In the case of the infinity-norm, the coherence is calculated as:

$$\lim_{\substack{i \ j \ 1 \ i \ j}} n \quad \frac{| \quad i \quad j|}{|| \quad i||_2|| \quad j||_2}$$

Examples

For example,

11/2/16

```
>>> import numpy as np
>>> import magni
>>> from magni.cs.indicators import calculate_coherence
>>> Phi = np.zeros((5, 9))
>>> Phi[0, 0] = Phi[1, 2] = Phi[2, 4] = Phi[3, 6] = Phi[4, 8] = 1
>>> Psi = magni.imaging.dictionaries.get_DCT((3, 3))
>>> for item in sorted(calculate_coherence(Phi, Psi).items()):
... print('{}-norm: {:.3f}'.format(*item))
0-norm: 0.222
1-norm: 0.335
inf-norm: 1.000
```

The above values can be calculated individually by specifying a norm:

```
>>> for norm in (0, 1, 2, np.inf):
... value = calculate_coherence(Phi, Psi, norm=norm)
... print('{}-norm: {:.3f}'.format(norm, value))
0-norm: 0.222
1-norm: 0.141
2-norm: 0.335
inf-norm: 1.000
```

magni.cs.indicators. calculate_mutual_coherence(Phi, Psi, norm=None)

Calculate the mutual coherence of Phi and Psi.

In the context of Compressive Sensing, mutual coherence usually refers to the maximum absolute correlation between two columns of Phi and Psi. This function allows the usage of a different normalised norm where the infinity-norm yields the usual case.

Parameters: • Phi

- **Phi** (magni.utils.matrices.Matrix or numpy.ndarray) The measurement matrix.
- **Psi** (*magni.utils.matrices.Matrix or numpy.ndarray*) The dictionary matrix.
- norm (int or float) The normalised norm used for the calculation (the default value is None which implies that the 0-, 1-, 2-, and infinity-norms are returned).

Returns:

mutual_coherence (*float or dict*) – The mutual_coherence value(s).

Notes

If *norm* is None, the function returns a dict containing the mutual coherence using the 0-, 1-, 2-, and infinity-norms. Otherwise, the function returns the mutual coherence using the specified norm.

The mutual coherence is calculated as:

$$rac{1}{mn}\sum_{i=1}^{m}\sum_{j=1}^{n}|_{i}$$

where m is the number of rows in Phi and n is the number of columns in Psi. In the case of the 0-norm, the mutual coherence is calculated as:

$$\frac{1}{mn} \sum_{i=1}^{m} \sum_{j=1}^{n} (| i \quad j |)$$

11/2/16 Page 100 of 209

where

() is 1 if a is non-zero and 0 otherwise. In the case of the infinity-norm, the mutual coherence is calculated as:

Examples

For example,

```
>>> import numpy as np
>>> import magni
>>> from magni.cs.indicators import calculate_mutual_coherence
>>> Phi = np.zeros((5, 9))
>>> Phi[0, 0] = Phi[1, 2] = Phi[2, 4] = Phi[3, 6] = Phi[4, 8] = 1
>>> Psi = magni.imaging.dictionaries.get_DCT((3, 3))
>>> for item in sorted(calculate_mutual_coherence(Phi, Psi).items()):
... print('{}-norm: {:.3f}'.format(*item))
0-norm: 0.889
1-norm: 0.298
2-norm: 0.333
inf-norm: 0.667
```

The above values can be calculated individually by specifying a norm:

```
>>> for norm in (0, 1, 2, np.inf):
... value = calculate_mutual_coherence(Phi, Psi, norm=norm)
... print('{}-norm: {:.3f}'.format(norm, value))
0-norm: 0.889
1-norm: 0.298
2-norm: 0.333
inf-norm: 0.667
```

magni.cs.indicators. calculate_relative_energy (Phi, Psi, method=None)

Calculate the relative energy of Phi Psi matrix product atoms.

Parameters: • Phi

- Phi (magni.utils.matrices.Matrix or numpy.ndarray) The measurement matrix.
- **Psi** (*magni.utils.matrices.Matrix or numpy.ndarray*) The dictionary matrix.
- **method** (*str*) The method used for summarising the relative energies of the Phi Psi matrix product atoms.

Returns:

relative_energy (*float or dict*) – The relative_energy summary value(s).

Notes

The summary *method* used is either 'mean' for mean value, 'std' for standard deviation, 'min' for minimum value, or 'diff' for difference between the maximum and minimum values.

If *method* is None, the function returns a dict containing all of the above summaries. Otherwise, the function returns the specified summary.

The relative energies, which are summarised by the given method, are calculated as:

11/2/16 Page 101 of 209

$$\frac{|| \quad 1||_2}{|| \quad 1||_2}, \quad , \frac{|| \quad n||_2}{|| \quad n||_2}$$

where *n* is the number of columns in *Psi*.

Examples

For example,

```
>>> import numpy as np
>>> import magni
>>> from magni.cs.indicators import calculate_relative_energy
>>> Phi = np.zeros((5, 9))
>>> Phi[0, 0] = Phi[1, 2] = Phi[2, 4] = Phi[3, 6] = Phi[4, 8] = 1
>>> Psi = magni.imaging.dictionaries.get_DCT((3, 3))
>>> for item in sorted(calculate_relative_energy(Phi, Psi).items()):
... print('{}: {:.3f}'.format(*item))
diff: 0.423
mean: 0.735
min: 0.577
std: 0.126
```

The above values can be calculated individually by specifying a norm:

```
>>> for method in ('mean', 'std', 'min', 'diff'):
... value = calculate_relative_energy(Phi, Psi, method=method)
... print('{}: {:.3f}'.format(method, value))
mean: 0.735
std: 0.126
min: 0.577
diff: 0.423
```

magni.imaging package

Subpackage providing functionality for image manipulation.

Routine listings

dictionaries

Module providing fast linear operations wrapped in matrix emulators.

domains

Module providing a multi domain image class.

evaluation

Module providing functions for evaluation of image reconstruction quality.

measurements

Module providing functions for constructing scan patterns for measurements.

preprocessing

Module providing functionality to remove tilt in images.

visualisation

Module providing functionality for visualising images.

mat2vec(x)

Function to reshape a matrix into vector by stacking columns.

11/2/16 Page 102 of 209

```
vec2mat(x, mn tuple)
```

Function to reshape a vector into a matrix.

```
double_mirror(img, fftstyle=True)
```

Function to mirror image in both the vertical and horisontal axes.

```
get inscribed masks(img, as vec=False)
```

Function to get inscribed masks covering the image.

Notes

See _util for documentation of mat2vec, vec2mat, double_mirror, and get_inscribed_masks.

Subpackages

magni.imaging.dictionaries package

Subpackage providing functionality for image dictionary manipulations.

Routine listings

```
get_DCT(shape, overcomplete_shape=None)
```

Get the DCT fast operation dictionary for the given image shape.

get_DFT(shape, overcomplete_shape=None)

Get the DFT fast operation dictionary for the given image shape.

def get_DFT_transform_matrix(N)

Return the normalised N-by-N discrete fourier transform (DFT) matrix.

analysis

Module providing functionality to analyse dictionaries.

utils

Module providing utility functions for the dictionaries subpackage.

Submodules

magni.imaging.dictionaries. fastops module

Module providing functionality related to linear transformations.

Routine listings

```
dct2(x, mn_tuple, overcomplete_mn_tuple=None)
```

2D discrete cosine transform.

idct2(x, mn_tuple, overcomplete_mn_tuple=None)

2D inverse discrete cosine transform.

dft2(x, mn_tuple, overcomplete_mn_tuple=None)

2D discrete Fourier transform.

idft2(x, mn_tuple, overcomplete_mn_tuple=None)

2D inverse discrete Fourier transform.

magni.imaging.dictionaries._fastops. **dct2**(*x*, *mn_tuple*, *overcomplete_mn_tuple=None*)

Apply the 2D Discrete Cosine Transform (DCT).

11/2/16 Page 103 of 209

x is assumed to be the column vector resulting from stacking the columns of the associated matrix which the transform is to be taken on.

- **Parameters: x** (*ndarray*) The m*n x 1 vector representing the associated column stacked matrix.
 - mn_tuple (tuple of int) -(m, n) m number of rows in the matrix represented by x; n number of columns in the matrix.
 - overcomplete_mn_tuple (tuple of int, optional) (mo, no) mo number of rows in the matrix represented by the function's output; no number of columns in the matrix.

Returns:

alpha (ndarray) – When overcomplete mn tuple is omitted: an m*n x 1 vector of coefficients scaled such that x = idct2(dct2(x)). When overcomplete mn tuple is supplied: an mo*no x 1 vector of coefficients scaled such that x = idct2(dct2(x)).

See also:

scipy.fftpack.dct()

1D DCT

Notes

The overcomplete transform feature invoked by supplying overcomplete mn tuple is only available for SciPy >= 0.16.1.

DCT The 2D uses the seperability property 1D DCT. of the http://stackoverflow.com/questions/14325795/scipys-fftpack-dct-and-idct Including the reshape operation, the full transform is: dct(dct(x.reshape(n, m).T).T).T.T.reshape(m * n, 1)

magni.imaging.dictionaries._fastops. idct2(x, mn_tuple, overcomplete_mn_tuple=None) Apply the 2D Inverse Discrete Cosine Transform (iDCT).

x is assumed to be the column vector resulting from stacking the columns of the associated matrix which the transform is to be taken on.

- **Parameters:** x (ndarray) The vector representing the associated column-stacked matrix. When *overcomplete mn tuple* is omitted: the shape must be m*n x 1. When *overcomplete mn tuple* is supplied: the shape must be mo*no x 1.
 - mn_tuple (tuple of int) (m, n) m number of rows in the associated matrix, n number of columns in the associated matrix. When overcomplete mn tuple is supplied, this is the shape of the function's output.
 - overcomplete mn tuple (tuple of int, optional) (mo, no) mo number of rows in the associated matrix, no number of columns in the associated matrix. When supplied, this is the shape of the function's input.

Returns:

ndarray - An m*n x 1 vector of coefficients scaled such that x = dct2(idct2(x)).

11/2/16 Page 104 of 209

See also:

scipy.fftpack.idct()

1D inverse DCT

Notes

The overcomplete transform feature invoked by supplying overcomplete mn tuple is only available for SciPy >= 0.16.1.

magni.imaging.dictionaries._fastops. dft2(x, mn tuple, overcomplete mn tuple=None)

Apply the 2D Discrete Fourier Transform (DFT).

x is assumed to be the column vector resulting from stacking the columns of the associated matrix which the transform is to be taken on.

- **Parameters: x** (*ndarray*) The m*n x 1 vector representing the associated column stacked matrix.
 - mn tuple (tuple of int) -(m, n) m number of rows in the matrix represented by *x*; *n* number of columns in the matrix.
 - overcomplete mn tuple (tuple of int, optional) (mo, no) mo number of rows in the matrix represented by the function's output; no number of columns in the matrix.

Returns:

ndarray - When overcomplete mn tuple is omitted: an m*n x 1 vector of coefficients scaled such that x = idft2(dft2(x)).overcomplete mn tuple is supplied: an mo*no x 1 vector of coefficients scaled such that x = idft2(dft2(x)).

See also:

numpy.fft.fft2()

The underlying 2D FFT used to compute the 2D DFT.

Notes

This is a normalised DFT, i.e. a normalisation constant of \overline{m} is used.

magni.imaging.dictionaries._fastops. idft2(x, mn_tuple, overcomplete_mn_tuple=None)

Apply the 2D Inverse Discrete Fourier Transform (iDFT).

x is assumed to be the column vector resulting from stacking the columns of the associated matrix which the transform is to be taken on.

11/2/16 Page 105 of 209

- **Parameters: x** (*ndarray*) The m*n x 1 vector representing the associated column stacked matrix.
 - mn tuple (tuple of int) -(m, n) m number of rows in the associated matrix, n number of columns in the associated matrix. When overcomplete mn tuple is supplied, this is the shape of the function's output.
 - overcomplete mn tuple (tuple of int, optional) (mo, no) mo number of rows in the associated matrix, no number of columns in the associated matrix. When supplied, this is the shape of the function's input.

Returns:

ndarray - An m*n x 1 vector of coefficients scaled such that x = dft2(idft2(x)).

See also:

numpy.fft.ifft2()

The underlying 2D iFFT used to compute the 2D iDFT.

Notes

This is a normalised iDFT, i.e. a normalisation constant of \overline{m} is used.

magni.imaging.dictionaries_fastops. **validate transform_fwd**(*x*, *mn_tuple*, overcomplete_mn_tuple=None)

Validate a 2D transform.

- **Parameters: x** (*ndarray*) The m*n x 1 vector representing the associated column stacked matrix.
 - mn_tuple (tuple of int) -(m, n) m number of rows in the matrix represented by x; n number of columns in the matrix.
 - overcomplete mn tuple (tuple of int, optional) (mo, no) mo number of rows in the matrix represented by the function's output; no number of columns in the matrix.

magni.imaging.dictionaries._fastops._validate_transform_bwd(x, mn_tuple, overcomplete mn tuple=None)

Validate a 2D transform.

- **Parameters: x** (*ndarray*) The m*n x 1 vector representing the associated column stacked matrix.
 - mn_tuple (tuple of int) (m, n) m number of rows in the associated matrix, *n* number of columns in the associated matrix.
 - overcomplete mn tuple (tuple of int, optional) (mo, no) mo number of rows in the associated matrix, no number of columns in the associated matrix.

magni.imaging.dictionaries._matrices module

Module providing fast linear operations wrapped in matrix emulators.

11/2/16 Page 106 of 209 Routine listings

get_DCT(shape, overcomplete_shape=None)

Get the DCT fast operation dictionary for the given image shape.

get DFT(shape, overcomplete shape=None)

Get the DFT fast operation dictionary for the given image shape.

See also:

magni.imaging.dictionaries._fastops

Fast linear operations.

magni.utils.matrices

Matrix emulators.

magni.imaging.dictionaries. matrices. **get DCT**(shape, overcomplete shape=None)

Get the DCT fast operation dictionary for the given image shape.

- **Parameters:** shape (*list or tuple*) The shape of the image for which the dictionary is the DCT dictionary.
 - overcomplete_shape (list or tuple, optional) The shape of the (overcomplete) frequency domain for the DCT dictionary. The entries must be greater than or equal to the corresponding entries in *shape*.

Returns:

matrix (*magni.utils.matrices.Matrix*) – The specified DCT dictionary.

See also:

magni.utils.matrices.Matrix()

The matrix emulator class.

Examples

Create a dummy image:

```
>>> import numpy as np, magni
>>> img = np.random.randn(64, 64)
>>> vec = magni.imaging.mat2vec(img)
```

Perform DCT in the ordinary way:

```
>>> dct_normal = magni.imaging.dictionaries. fastops.dct2(vec, img.shape)
```

Perform DCT using the present function:

```
>>> from magni.imaging.dictionaries import get_DCT
>>> matrix = get_DCT(img.shape)
>>> dct_matrix = matrix.T.dot(vec)
```

Check that the two ways produce the same result:

```
>>> np.allclose(dct_matrix, dct_normal)
True
```

11/2/16 Page 107 of 209

Compute the overcomplete transform (and back again) and check that the resulting image is identical to the original. Notice how this example first ensures that the necessary version of SciPy is available:

```
>>> from pkg resources import parse version
>>> from scipy import version as scipy version
>>> if parse_version(_scipy_version) >= parse_version('0.16.0'):
... matrix = get_DCT(img.shape, img.shape)
... dct_matrix = matrix.T.dot(vec)
... vec_roundtrip = matrix.dot(dct_matrix)
... np.allclose(vec, vec_roundtrip)
... else:
... True
True
```

magni.imaging.dictionaries._matrices. **get_DFT**(*shape*, *overcomplete_shape=None*)

Get the DFT fast operation dictionary for the given image shape.

- **Parameters:** shape (*list or tuple*) The shape of the image for which the dictionary is the DFT dictionary.
 - overcomplete_shape (list or tuple, optional) The shape of the (overcomplete) frequency domain for the DFT dictionary. The entries must be greater than or equal to the corresponding entries in *shape*.

Returns:

matrix (*magni.utils.matrices.Matrix*) – The specified DFT dictionary.

See also:

magni.utils.matrices.Matrix()

The matrix emulator class.

Examples

Create a dummy image:

```
>>> import numpy as np, magni
>>> img = np.random.randn(64, 64)
>>> vec = magni.imaging.mat2vec(img)
```

Perform DFT in the ordinary way:

```
>>> dft_normal = magni.imaging.dictionaries._fastops.dft2(vec, img.shape)
```

Perform DFT using the present function:

```
>>> from magni.imaging.dictionaries import get_DFT
>>> matrix = get DFT(img.shape)
>>> dft matrix = matrix.conj().T.dot(vec)
```

Check that the two ways produce the same result:

```
>>> np.allclose(dft matrix, dft normal)
True
```

Compute the overcomplete transform (and back again):

11/2/16 Page 108 of 209

```
>>> matrix = get_DFT(img.shape, img.shape)
>>> dft_matrix = matrix.conj().T.dot(vec)
>>> vec_roundtrip = matrix.dot(dft_matrix)
```

Check that the twice transformed image is identical to the original:

```
>>> np.allclose(vec, vec_roundtrip)
True
```

magni.imaging.dictionaries. mtx1D module

Module providing 1D matrices for building 2D separable transforms.

Routine listings

```
get DCT transform matrix(N)
```

Return the normalised N-by-N discrete cosine transform (DCT) matrix.

```
get DFT transform matrix(N)
```

Return the normalised N-by-N discrete fourier transform (DFT) matrix.

```
magni.imaging.dictionaries._mtx1D. get DCT transform matrix(N)
```

Return the normalised N-by-N discrete cosine transform (DCT) matrix.

Applying the returned transform matrix to a vector x: D.dot(x) yields the DCT of x. Applying the returned transform matrix to a matrix A: D.dot(A) applies the DCT to the columns of A. Taking D.dot(A.dot(D.T)) applies the DCT to both columns and rows, i.e. a full 2D separable DCT transform. The inverse transform (the 1D IDCT) is D.T.

```
Parameters: N (int) – The size of the DCT transform matrix to
```

return.

Returns: D (*ndarray*) – The DCT transform matrix.

Notes

The returned DCT matrix normalised such that is consitutes a orthonormal transform as given by equations (2.119) and (2.120) in [1].

References

[1] A.N. Akansu, R.A. Haddad, and P.R. Haddad, *Multiresolution Signal Decomposition: Transforms, Subbands, and Wavelets*, Academic Press, 2000.

Examples

For example, get a 5-by-5 DCT matrix

11/2/16 Page 109 of 209

and apply the 2D DCT transform to a dummy image

```
>>> np.random.seed(6021)
>>> img = np.random.randn(5, 5)
>>> img_dct = D.dot(img.dot(D.T))
>>> np.round(img_dct, 4)
array([[-0.5247, -0.0225, 0.9098, 0.369, -0.477],
        [ 1.7309, -0.4142, 1.9455, -0.6726, -1.3676],
        [ 0.6987, 0.5355, 0.7213, -0.8498, -0.1023],
        [ 0.0078, -0.0545, 0.3649, -1.4694, 1.732 ],
        [-1.5864, 0.156, 0.8932, -0.8091, 0.5056]])
```

magni.imaging.dictionaries._mtx1D. get_DFT_transform_matrix(N)

Return the normalised N-by-N discrete fourier transform (DFT) matrix.

Applying the returned transform matrix to a vector x: D.dot(x) yields the DFT of x. Applying the returned transform matrix to a matrix A: D.dot(A) applies the DFT to the columns of A. Taking D.dot(A.dot(D.T)) applies the DFT to both columns and rows, i.e. a full 2D separable DFT transform. The inverse transform (the 1D IDFT) is D.T.

Parameters: N (*int*) – The size of the DFT transform matrix to

return.

Returns: D (*ndarray*) – The DFT transform matrix.

See also:

scipy.linalg.dft()

The function used to generate the DFT transform matrix.

Notes

The returned DFT matrix normalised such that is consitutes a orthonormal transform as given by equations (2.105) and (2.109) in [2].

References

[2] A.N. Akansu, R.A. Haddad, and P.R. Haddad, *Multiresolution Signal Decomposition: Transforms, Subbands, and Wavelets*, Academic Press, 2000.

Examples

For example, get a 5-by-5 DFT matrix

and apply the 2D DFT transform to a dummy image

11/2/16 Page 110 of 209

```
>>> np.random.seed(6021)
>>> img = np.random.randn(5, 5)
>>> img_dft = D.dot(img.dot(D.T))
>>> np.round(img_dft, 2)
array([[-0.52+0.j , 0.44+0.48j, 0.11-0.39j, 0.11+0.39j, 0.44-0.48j],
        [ 1.04-0.59j, 1.32+0.13j, -1.20-0.39j, 0.35+0.66j, 0.19-0.36j],
        [ 0.18-1.24j, 0.75+0.44j, -0.75+0.72j, -0.52-0.8j , 0.77+0.13j],
        [ 0.18+1.24j, 0.77-0.13j, -0.52+0.8j , -0.75-0.72j, 0.75-0.44j],
        [ 1.04+0.59j, 0.19+0.36j, 0.35-0.66j, -1.20+0.39j, 1.32-0.13j]])
```

which may be shifted to have the zero-frequency component at the center of the spectrum

```
>>> np.round(scipy.fftpack.fftshift(img_dft), 2)
array([[-0.75-0.72j, 0.75-0.44j, 0.18+1.24j, 0.77-0.13j, -0.52+0.8j],
        [-1.20+0.39j, 1.32-0.13j, 1.04+0.59j, 0.19+0.36j, 0.35-0.66j],
        [ 0.11+0.39j, 0.44-0.48j, -0.52+0.j , 0.44+0.48j, 0.11-0.39j],
        [ 0.35+0.66j, 0.19-0.36j, 1.04-0.59j, 1.32+0.13j, -1.20-0.39j],
        [-0.52-0.8j , 0.77+0.13j, 0.18-1.24j, 0.75+0.44j, -0.75+0.72j]])
```

magni.imaging.dictionaries. visualisations module

Module providing functionality for visualising dictionary coefficients.

Routine listings

visualise DCT(shape)

Function for visualising DCT coefficients.

visualise DFT(shape)

Function for visualising DFT coefficients.

magni.imaging.dictionaries. visualisations. **visualise DCT**(*shape*)

Return utilities for visualising DCT coefficients.

A handle to a function to transform the coefficients into a 'displayable' format is returned along with a tuple of ranges of the axes in the 2D coefficient plane.

Parameter Returns:

Parameters: shape (*tuple*) – The shape of the 2D DCT being visualised.

- **display_coefficients** (*Function*) The function used to transform coefficients into a 'displayable' format.
- axes_extent (tuple) The ranges of the axes in the 2D coefficient plane.

Notes

The display_coefficients function takes log10 to the absolute value of the transform cofficient vector given to it as an argument. The returned displayable coefficients is a matrix.

The axes_extent consists of (abcissa_min, abcissa_max, ordinate_min, ordinate_max).

magni.imaging.dictionaries. visualisations. **visualise DFT**(*shape*)

Return utilities for visualising DFT coefficients.

A handle to a function to transform the coefficients into a 'displayable' format is returned along with a tuple of ranges of the axes in the 2D coefficient plane.

11/2/16 Page 111 of 209

Returns:

Parameters: shape (*tuple*) – The shape of the 2D DFT being visualised.

- display coefficients (Function) The function used to transform coefficients into a 'displayable' format.
- axes_extent (tuple) The ranges of the axes in the 2D coefficient plane.

Notes

The display coefficients function takes log10 to the absolute value of the transform cofficient vector given to it as an argument. The returned displayable coefficients is a matrix that is flipped up/down and fftshifted.

The axes extent consists of (abcissa min, abcissa max, ordinate min, ordinate max).

magni.imaging.dictionaries.analysis module

Module providing functionality for the analysis of dictionaries.

Routine listings

get reconstructions(img, transform, fractions)

Function to get a set of transform reconstructions.

- show coefficient histogram(img, transforms, bins=None, range=None, output path=None, fig ext='pdf') Function to show a transform coefficient histogram.
- show psnr_energy_rolloff(img, reconstructions, fractions, return_vals=False, output path=None, fig ext='pdf') Function to show PSNR and retained energy rolloff of reconstructions.
- show reconstructions(coefficients, reconstructions, transform, fractions, output path=None, fig_ext='pdf') Function to show transforms reconstructions and coefficients.
- show sorted coefficients(img, transforms, output path=None, fig ext='pdf') Function to show a plot of transform coefficients in sorted order.
- show transform coefficients(img, transforms, output path=None, fig ext='pdf') Function to show transform coefficients.
- show transform quantiles(img, transform, fraction=1.0, area mask=None) Function to show quantiles of transform coefficients.

magni.imaging.dictionaries.analysis. **qet reconstructions** (*img*, *transform*, *fractions*)

Return transform reconstructions with different fractions of coefficents.

The image img is transform coded using the specified transform. Reconstructions for the fractions of transform coefficients kept are returned along with the coefficients used in the reconstructions.

- **Parameters:** img (ndarray) The image to get reconstructions of.
 - **transform** (*str*) The transform to use to obtain the reconstructions.
 - fractions (list or tuple) The fractions of coefficients to be used in the reconstructions.

11/2/16 Page 112 of 209 Returns:

- coefficients (list) The list of coefficients (each an ndarray) used in the reconstructions.
- **reconstructions** (*list*) The list of reconstructions.

Examples

Get reconstructions from DCT based on 20% and 40% of the coefficients:

```
>>> import numpy as np
>>> from magni.imaging.dictionaries.analysis import get_reconstructions
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transform = 'DCT'
>>> fractions = (0.2, 0.4)
>>> coefs, recons = get reconstructions(img, transform, fractions)
>>> len(recons)
>>> tuple(int(s) for s in coefs[0].shape)
>>> tuple(int(s) for s in recons[0].shape)
```

magni.imaging.dictionaries.analysis. **show_coefficient_histogram** (*img*, *transforms*, *bins=None*, range=None, output_path=None, fig_ext='pdf')

Show a histogram of coefficient values for different transforms.

A histogram of the transform coefficient values for img using the different transforms are shown. If output_path is not None, the resulting figure and data used in the figure are saved.

- **Parameters:** img (*ndarray*) The image to get transform coefficients for.
 - transforms (list or tuple) The names as strings of the transforms to use.
 - **bins** (*int*) The number of bins to use in the histogram (the default is None, which implies that the number of bins is determined based on the size of img).
 - range (tuple) The lower and upper range of the bins to use in the histogram (the default is None, which implies that the min and max values of img are used).
 - output_path (str) The output path (see notes below) to save the figure and data to (the default is None, which implies that the figure and data are not saved).
 - fig ext (str) The figure extension determining the format of the saved figure (the default is 'pdf' which implies that the figure is saved as a PDF).

See also:

matplotlib.pyplot.hist()

The underlying histogram plot function.

Notes

The *output path* is specified as a path to a folder + an optional prefix to the file name. The remaining file name is fixed. If e.g., the fixed part of the file name was 'plot', then:

11/2/16 Page 113 of 209

- output_path = '/home/user/' would save the figure under /home/user/ with the name plot.pdf.
- output path = '/home/user/best' would save the figure under /home/user with the name best plot.pdf.

In addition to the saved figures, an annotated and chased HDF database with the data used to create the figures are also saved. The name of the HDF database is the same as for the figure with the exception that the file extension is '.hdf5'.

Examples

Save a coefficient histogram using the DCT and the DFT as transforms

```
>>> import os, numpy as np
>>> from magni.imaging.dictionaries import analysis as _a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transforms = ('DCT', 'DFT')
>>> output_path = './histogram_test'
>>> _a.show_coefficient_histogram(img, transforms, output_path=output_path)
>>> current_dir = os.listdir('./')
>>> for file in sorted(current_dir):
... if 'histogram_test' in file:
       print(file)
histogram_test_coefficient_histogram.hdf5
histogram test coefficient histogram.pdf
```

magni.imaging.dictionaries.analysis. **show_psnr_energy_rolloff**(img, reconstructions, fractions, return vals=False, output path=None, fig ext='pdf')

Show the PSNR and energy rolloff for the reconstructions.

A plot of the Peak Signal to Noise Ratio (PSNR) and retained energy in the recontructions versus the fractions of coefficients used in the reconstructions is shown. If return vals is True, the data used in the plot is returned. If output path is not None, the resulting figure and data used in the figure are saved.

- **Parameters:** img (*ndarray*) The image which the reconstructions are based on.
 - reconstructions (list or tuple) The reconstructions (each an ndarray) to show rolloff for.
 - fractions (list or tuple) The fractions of coefficients used in the reconstructions.
 - return vals (bool) The flag indicating wheter or not to return the PSNR and energy values (the default is False, which indicate that the values are not returned).
 - output path (str) The output path (see notes below) to save the figure and data to (the default is None, which implies that the figure and data are not saved).
 - fig ext (str) The figure extension determining the format of the saved figure (the default is 'pdf' which implies that the figure is saved as a PDF).

Returns:

- psnrs (ndarray) The PSNR values shown in the figure (only returned if return vals=True).
- energy (ndarray) The retained energy values shown in the figure (only returned if return vals=True).

Notes

11/2/16 Page 114 of 209 The *output path* is specified as a path to a folder + an optional prefix to the file name. The remaining file name is fixed. If e.g., the fixed part of the file name was 'plot', then:

- output_path = '/home/user/' would save the figure under /home/user/ with the name plot.pdf.
- output path = '/home/user/best' would save the figure under /home/user with the name best plot.pdf.

In addition to the saved figures, an annotated and chased HDF database with the data used to create the figures are also saved. The name of the HDF database is the same as for the figure with the exception that the file extension is '.hdf5'.

Examples

Save a PSNR and energy rolloff plot for reconstructions bases on the DCT:

```
>>> import os, numpy as np
>>> from magni.imaging.dictionaries import analysis as _a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transform = 'DCT'
>>> fractions = (0.2, 0.4)
>>> coefs, recons = _a.get_reconstructions(img, transform, fractions)
>>> o p = './rolloff test'
>>> _a.show_psnr_energy_rolloff(img, recons, fractions, output_path=o_p)
>>> current dir = os.listdir('./')
>>> for file in sorted(current_dir):
... if 'rolloff_test' in file:
       print(file)
rolloff_test_psnr_energy_rolloff.hdf5
rolloff_test_psnr_energy_rolloff.pdf
```

magni.imaging.dictionaries.analysis. **show_reconstructions** (coefficients, reconstructions, transform, fractions, output_path=None, fig_ext='pdf')

Show reconstructions and corresponding coefficients.

- **Parameters:** coefficients (list or tuple) The coefficients (each an ndarray) used in the reconstructions.
 - reconstructions (list or tuple) The reconstructions (each an ndarray) to show.
 - **transform** (*str*) The transform used to obtain the reconstructions.
 - fractions (list or tuple) The fractions of coefficents used in the reconstructions.
 - **output path** (*str*) The output path (see notes below) to save the figure and data to (the default is None, which implies that the figure and data are not saved).
 - fig_ext (str) The figure extension determining the format of the saved figure (the default is 'pdf' which implies that the figure is saved as a PDF).

Notes

The *output path* is specified as a path to a folder + an optional prefix to the file name. The remaining file name is fixed. If e.g., the fixed part of the file name was 'plot', then:

• output_path = '/home/user/' would save the figure under /home/user/ with the name plot.pdf.

11/2/16 Page 115 of 209 output path = '/home/user/best' would save the figure under /home/user with the name best_plot.pdf.

In addition to the saved figures, an annotated and chased HDF database with the data used to create the figures are also saved. The name of the HDF database is the same as for the figure with the exception that the file extension is '.hdf5'.

Examples

Save images of coefficients and reconstructions based on the DCT:

```
>>> import os, numpy as np
>>> from magni.imaging.dictionaries import analysis as _a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> trans = 'DCT'
>>> fracs = (0.2, 0.4)
>>> coefs, recons = _a.get_reconstructions(img, trans, fracs)
>>> o_p = './reconstruction test'
>>> a.show reconstructions(coefs, recons, trans, fracs, output path=o p)
>>> current dir = os.listdir('./')
>>> for file in sorted(current_dir):
... if 'reconstruction_test' in file:
       print(file)
reconstruction test reconstruction coefficients.hdf5
reconstruction_test_reconstruction_coefficients.pdf
reconstruction_test_reconstructions.hdf5
reconstruction test reconstructions.pdf
```

magni.imaging.dictionaries.analysis. **show_sorted_coefficients**(*img*, *transforms*, *output_path=None*, fig ext='pdf')

Show the transform coefficient values in sorted order.

A plot of the sorted coefficient values vs array index number is shown. Ifoutput path is not None, the resulting figure and data used in the figure are saved.

- **Parameters:** img (ndarray) The image to show the sorted transform coefficients values for.
 - transforms (list or tuple) The names as strings of the transforms to
 - output_path (str) The output path (see notes below) to save the figure and data to (the default is None, which implies that the figure and data are not saved).
 - fig_ext (str) The figure extension determining the format of the saved figure (the default is 'pdf' which implies that the figure is saved as a PDF).

Notes

The *output_path* is specified as a path to a folder + an optional prefix to the file name. The remaining file name is fixed. If e.g, the fixed part of the file name was 'plot', then:

- output_path = '/home/user/' would save the figure under /home/user/ with the name plot.pdf.
- output_path = '/home/user/best' would save the figure under /home/user with the name best_plot.pdf.

In addition to the saved figures, an annotated and chased HDF database with the data

11/2/16 Page 116 of 209

used to create the figures are also saved. The name of the HDF database is the same as for the figure with the exception that the file extension is '.hdf5'.

Examples

Save a sorted transform coefficient plot for the DCT and DFT transforms:

```
>>> import os, numpy as np
>>> from magni.imaging.dictionaries import analysis as a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transforms = ('DCT', 'DFT')
>>> output_path = './sorted_test'
>>> _a.show_sorted_coefficients(img, transforms, output_path=output_path)
>>> current_dir = os.listdir('./')
>>> for file in sorted(current_dir):
... if 'sorted_test' in file:
       print(file)
sorted test sorted coefficients.hdf5
sorted test sorted coefficients.pdf
```

magni.imaging.dictionaries.analysis.show transform coefficients (img, transforms, output path=None, fig ext='pdf')

Show the transform coefficients.

The transform coefficient of img are shown for the transforms. If output path is not None, the resulting figure and data used in the figure are saved.

- **Parameters:** img (*ndarray*) The image to show the transform coefficients for.
 - transforms (list or tuple) The names as strings of the transforms to
 - **output path** (*str*) The output path (see notes below) to save the figure and data to (the default is None, which implies that the figure and data are not saved).
 - fig ext (str) The figure extension determining the format of the saved figure (the default is 'pdf' which implies that the figure is saved as a PDF).

Notes

The *output path* is specified as a path to a folder + an optional prefix to the file name. The remaining file name is fixed. If e.g., the fixed part of the file name was 'plot', then:

- output_path = '/home/user/' would save the figure under /home/user/ with the name plot.pdf.
- output_path = '/home/user/best' would save the figure under /home/user with the name best plot.pdf.

In addition to the saved figures, an annotated and chased HDF database with the data used to create the figures are also saved. The name of the HDF database is the same as for the figure with the exception that the file extension is '.hdf5'.

Examples

Save a figure showing the coefficients for the DCT and DFT transforms:

11/2/16 Page 117 of 209

```
>>> import os, numpy as np
>>> from magni.imaging.dictionaries import analysis as a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transforms = ('DCT', 'DFT')
>>> o_p = './coefficient_test'
>>> _a.show_transform_coefficients(img, transforms, output_path=o_p)
>>> current dir = os.listdir('./')
>>> for file in sorted(current_dir):
... if 'coefficient_test' in file:
       print(file)
coefficient_test_transform_coefficients.hdf5
coefficient_test_transform_coefficients.pdf
```

magni.imaging.dictionaries.analysis. show transform quantiles (img, transform, fraction=1.0, area mask=None, ax=None)

Show a plot of the quantiles of the transform coefficients.

The fraction of transform coefficients holding the most energy for the image img is considered. The four quantiles within this fraction of coefficients are illustrated in the transform domain by showing coefficients between the different quantiles in different colours. If an area mask is specified, only this area in the plot is highlighted whereas the rest is darkened.

- **Parameters:** img (ndarray) The image to show the transform quantiles for.
 - **transform** (*str*) The transform to use.
 - fraction (float) The fraction of coefficients used in the quantiles calculation (the default value is 1.0, which implies that all coefficients are used).
 - area_mask (ndarray) Bool array of the same shape asimg which indicates the area of the image to highlight (the default value is None, which implies that no particular part of the image is highlighted).
 - ax (matplotlib.axes.Axes) The axes on which the image is displayed (the default is None, which implies that a separate figure is created).

Returns:

coef count (*dict*) – Different counts of coeffcients within the *area mask* (only returned if area_mask is not None).

Notes

The ticks on the colorbar shown below the figure are the percentiles of the entire set of coefficients corresponding to the quantiles with fraction of coefficients. For instance, if fraction is 0.10, then the percentiles are 92.5, 95.0, 97.5, 100.0, corresponding to the four quantiles within the 10 percent coefficients holding the most energy.

The *coef count* dictionary holds the following keys:

- C total: Total number of considered coefficients.
- Q potential: Number of potential coefficients within mask area.
- P fraction: The fraction of Q potential to the pixel count in img.
- Q total: Total number of (considered) coefficients within mask area.
- Q fraction: The fraction of Q total to Q potential
- QC fraction: The fraction of Q total to C total.
- Q0-Q1 : Number of coefficients smaller than the first quantile.
- Q1-Q2: Number of coefficients between the first and second quantile.
- Q2-Q3: Number of coefficients between the second and third quantile.
- Q3-Q4: Number of coefficients between the third and fourth quantile.

11/2/16 Page 118 of 209 Each of the QX-QY holds a tuple containing two values:

- 1. The number of coefficients.
- 2. The fraction of the number of coefficients to Q total.

Examples

For example, show quantiles for a fraction of 0.2 of the DCT coefficients:

```
>>> import numpy as np
>>> from magni.imaging.dictionaries import analysis as _a
>>> img = np.arange(64).astype(np.float).reshape(8, 8)
>>> transforms = 'DCT'
>>> fraction = 0.2
>>> _a.show_transform_quantiles(img, transform, fraction=fraction)
```

magni.imaging.dictionaries.analysis. **_save_output**(output_path, name, fig, fig_ext, datasets)

Save figure and data output.

- **Parameters:** output_path (*str*) The output_path to save to.
 - **name** (*str*) The 'fixed' part of the file name saved to.
 - **fig** (*matplotlib.figure.Figure*) The figure instance to save.
 - **fig_ext** (*str*) The file extension to use for the saved figure.
 - datasets (dict) The dict of dicts for datasets to save in a HDF database.

magni.imaging.dictionaries.utils module

Module providing support functionality for the dictionaries subpackage.

Routine listings

```
get_function_handle(type_, transform)
```

Function to get a function handle to a transform method.

```
get transform names()
```

Function to get a tuple of names of the available transforms.

magni.imaging.dictionaries.utils. **get_function_handle**(*type_, transform*)

Return a function handle to a given transform method.

- **Parameters:** type_ ({'matrix', 'transform_matrix', 'visualisation'}) Identifier of the type of method to return a handle to.
 - **transform** (*str*) Identifier of the transform method to return a handle

Returns:

f_handle (*function*) – Handle to *transform*.

Notes

The available types are:

- matrix: 2D transform matrix object.
- transform matrix: 1D transform matrix.
- visualisation: Coefficient visualisation preparation function.

11/2/16 Page 119 of 209

Examples

For example, return a handle to the matrix method for a DCT:

```
>>> from magni.imaging.dictionaries.utils import get_function_handle
>>> get_function_handle('matrix', 'DCT').__name__
'get_DCT'
```

or return a handle to the 1D DCT transform matrix method:

```
>>> get_function_handle('transform_matrix', 'DCT').__name__
'get_DCT_transform_matrix'
```

or return a handle to the visualisation method for a DFT:

```
>>> get_function_handle('visualisation', 'DFT').__name__
'visualise_DFT'
```

magni.imaging.dictionaries.utils. get_transform_names()

Return a tuple of names of the available transforms.

```
Returns: names (tuple) – The tuple of names of available transforms.
```

Examples

For example, get transform names and extract 'DCT' and 'DFT'

```
>>> from magni.imaging.dictionaries.utils import get_transform_names
>>> names = get_transform_names()
>>> tuple(sorted(name for name in names if name in ('DCT', 'DFT')))
('DCT', 'DFT')
```

and a handle to corresponding visualisation method for the DCT

```
>>> from magni.imaging.dictionaries.utils import get_function_handle
>>> f_handles = tuple(get_function_handle('visualisation', name)
... for name in names)
>>> tuple(f_h.__name__ for f_h in f_handles if 'DCT' in f_h.__name__)
('visualise_DCT',)
```

magni.imaging.measurements package

Subpackage providing functionality for constructing and visualising scan patterns for measurements.

This subpackage provides several pairs of scan pattern functions. The first function, named *_sample_surface, is used for sampling a given surface. The second function, named *_sample_image, is a wrapper that provides a pixel-oriented interface to the first function. In addition to these pairs of scan pattern functions, the module provides auxillary functions that may be used to visualise the scan patterns.

Routine listings

```
construct_measurement_matrix(coords, h, w)
```

Function for constructing a measurement matrix.

11/2/16 Page 120 of 209

- construct_pixel_mask(h, w, pixels)
 Construct a binary pixel mask.
- lissajous_sample_image(h, w, scan_length, num_points, f_y=1., f_x=1., theta_y=0., theta_x=np.pi / 2) Function for lissajous sampling an image.
- lissajous_sample_surface(I, w, speed, sample_rate, time, f_y=1., f_x=1., theta_y=0., theta_x=np.pi / 2, speed_mode=0) Function for lissajous sampling a surface.
- plot_pattern(I, w, coords, mode, output_path=None) Function for visualising a scan pattern.
- plot_pixel_mask(h, w, pixels, output_path=None)
 Function for visualising a pixel mask obtained from a scan pattern.
- random_line_sample_image(h, w, scan_length, num_points, discrete=None, seed=None) Function for random line sampling an image.
- random_line_sample_surface(I, w, speed, sample_rate, time, discrete=None, seed=None) Function for random line sampling a surface.
- spiral_sample_image(h, w, scan_length, num_points, rect_area=False) Function for spiral sampling an image.
- spiral_sample_surface(I, w, speed, sample_rate, time, rect_area=False) Function for spiral sampling a surface.
- square_spiral_sample_image(h, w, scan_length, num_points) Function for square spiral sampling an image.
- square_spiral_sample_surface(I, w, speed, sample_rate, time)
 Function for square spiral sampling a surface.
- uniform_line_sample_image(h, w, scan_length, num_points)
 Function for uniform line sampling an image.
- uniform_line_sample_surface(I, w, speed, sample_rate, time) Function for uniform line sampling a surface.
- uniform_rotated_line_sample_image(h, w, scan_length, num_points, angle=0., follow edge=True) Function for uniform rotated line sampling an image.
- uniform_rotated_line_sample_surface(I, w, speed, sample_rate, time, angle=0., follow_edge=True) Function for uniform rotated line sampling a surface.
- unique_pixels(coords)

Function for determining unique pixels from a set of coordinates.

- zigzag_sample_image(h, w, scan_length, num_points, angle=np.pi / 20) Function for zigzag sampling an image.
- zigzag_sample_surface(I, w, speed, sample_rate, time, angle=np.pi / 20) Function for zigzag sampling a surface.

Notes

In principle, most of the scan pattern related parameters need only be positive. However, it is assumed that the following requirements are fulfilled:

11/2/16 Page 121 of 209

```
Minimum length of scan area:
```

1 nm

Minimum width of scan area:

1 nm

Minimum scan speed:

1 nm/s

Minimum sample rate:

1 Hz

Minimum scan time:

1 s

Minimum scan length:

1 nm

Minimum number of scan points:

1

Examples

Sample a surface using a lissajous pattern:

```
>>> from magni.imaging.measurements import lissajous_sample_surface
>>> I = 13.0; w = 13.0; speed = 4.0; time = 27.0; sample_rate = 3.0;
>>> coords = lissajous_sample_surface(I, w, speed, sample_rate, time,
... speed_mode=1)
```

Display the resulting pattern:

```
>>> from magni.imaging.measurements import plot_pattern
>>> plot_pattern(l, w, coords, 'surface')
```

Sample a 128x128 pixel image using a spiral pattern:

```
>>> from magni.imaging.measurements import spiral_sample_image
>>> h = 128; w = 128; scan_length = 1000.0; num_points = 200;
>>> coords = spiral_sample_image(h, w, scan_length, num_points)
```

Display the resulting pattern:

```
>>> plot_pattern(h, w, coords, 'image')
```

Find the corresponding unique pixels and plot the pixel mask:

```
>>> from magni.imaging.measurements import unique_pixels, plot_pixel_mask
>>> unique_pixels = unique_pixels(coords)
>>> plot_pixel_mask(h, w, unique_pixels)
```

Submodules

magni.imaging.measurements. lissajous module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
lissajous_sample_image(h, w, scan_length, num_points, f_y=1., f_x=1., theta_y=0., theta_x=np.pi / 2) Function for lissajous sampling an image.
```

11/2/16 Page 122 of 209

lissajous sample surface(l, w, speed, sample rate, time, f y=1., f x=1., theta y=0., theta x=np.pi / 2, speed mode=0) Function for lissajous sampling a surface.

magni.imaging.measurements._lissajous. **lissajous sample image**(*h*, *w*, *scan length*, *num points*, f y=1.0, f x=1.0, theta y=0.0, theta x=1.5707963267948966

Sample an image using a lissajous pattern.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a lissajous pattern are determined. The scan length determines the length of the path scanned whereas *num points* indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - scan length (*float*) The length of the path to scan in units of pixels.
 - **num_points** (*int*) The number of samples to take on the scanned path.
 - **f v** (*float*) The frequency of the y-sinusoid (the default value is 1.0).
 - **f x** (*float*) The frequency of the x-sinusoid (the default value is 1.0).
 - theta_y (float) The starting phase of the y-sinusoid (the default is
 - theta_x (float) The starting phase of the x-sinusoid (the default is pi / 2).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height *h* is measured along the y-axis.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import lissajous_sample_image
>>> h = 10
>>> W = 10
>>> scan length = 50.0
>>> num points = 12
>>> np.set printoptions(suppress=True)
>>> lissajous sample image(h, w, scan length, num points)
array([[ 5.
             , 9.5
    [1.40370042, 7.70492686],
    [ 0.67656563, 3.75183526],
    [3.39871123, 0.79454232],
    [7.39838148, 1.19240676],
    [ 9.48459832, 4.62800824],
    [7.99295651, 8.36038857],
    [ 4.11350322, 9.41181634],
    [\ 0.94130617,\ 6.94345168],
    [1.0071768, 2.92458128],
    [4.25856283, 0.56150128],
     8.10147506, 1.7395012],
    [ 9.4699986 , 5.51876059]])
```

magni.imaging.measurements. lissajous. **lissajous sample surface**(*I*, *w*, *speed*, *sample rate*, *time*,

11/2/16 Page 123 of 209 f y=1.0, f x=1.0, theta y=0.0, theta x=1.5707963267948966, speed mode=0)Sample a surface area using a lissajous pattern.

The coordinates (in units of meters) resulting from sampling an area of sizel times w using a lissajous pattern are determined. The scanned path is determined from the probe speed and the scan time.

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - w (float) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample_rate** (*float*) The sample rate in units of Hertz.
 - **time** (*float*) The scan time in units of seconds.
 - **f_y** (*float*) The frequency of the y-sinusoid (the default value is 1.0).
 - **f x** (*float*) The frequency of the x-sinusoid (the default value is 1.0).
 - theta y (float) The starting phase of the y-sinusoid (the default is 0.0).
 - theta_x (float) The starting phase of the x-sinusoid (the default is pi
 - **speed mode** (*int*) The speed mode used to select sampling points (the default is 0 which implies that the speed argument determines the speed, and f y and f x determine the ratio between the relative frequencies used).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length / is measured along the y-axis.

Generally, the lissajous sampling pattern does not provide constant speed, and this cannot be compensated for without violating f y, f x, or both. Therefore, speed mode allows the user to determine how this issue is handled: In speed_mode 0, constant speed equal to *speed* is ensured by non-uniform sampling of a lissajous curve, whereby f y and f x are not constant frequencies. In speed mode 1, average speed equal to speed is ensured by scaling f y and f x by the same constant. In speed mode 2, f y and f x are kept constant and the speed is only used to determine the path length in combination with time.

Examples

For example,

11/2/16 Page 124 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import lissajous_sample_surface
>>> w = 1e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> lissajous_sample_surface(I, w, speed, sample_rate, time)
array([[ 0.0000005, 0.000001 ],
    [\ 0.00000001,\ 0.00000058],
    [\ 0.00000033,\ 0.00000003],
    [ 0.00000094, 0.00000025],
    [ 0.00000082, 0.00000089], [ 0.00000017, 0.00000088], [ 0.00000007, 0.00000024],
    [ 0.00000068, 0.00000003],
    [ 0.00000099, 0.0000006 ],
    [0.00000048, 0.000001],
            , 0.0000057],
    [0.00000035, 0.00000002],
    [0.00000094, 0.00000027]])
```

magni.imaging.measurements. matrices module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
construct measurement matrix(coords, h, w)
```

Function for constructing a measurement matrix.

magni.imaging.measurements. matrices. **construct measurement matrix** (coords, h, w)

Construct a measurement matrix extracting the specified measurements.

- **Parameters:** coords (*ndarray*) The *k* floating point coordinates arranged into a 2D array where each row is a coordinate pair (x, y), such that *coords* has size $k \times 2$.
 - **h** (*int*) The height of the image measured in pixels.
 - **w** (*int*) The width of the image measured in pixels.

Returns:

Phi (magni.utils.matrices.Matrix) - The constructed measurement matrix.

See also:

magni.utils.matrices.Matrix()

The matrix emulator class.

Notes

The function constructs two functions: one for extracting pixels at the coordinates specified and one for the transposed operation. These functions are then wrapped by a matrix emulator which is returned.

Examples

Create a dummy 5 by 5 pixel image and an example sampling pattern:

11/2/16 Page 125 of 209

```
>>> import numpy as np, magni
>>> img = np.arange(25, dtype=np.float).reshape(5, 5)
>>> vec = magni.imaging.mat2vec(img)
>>> coords = magni.imaging.measurements.uniform line sample image(
          5, 5, 16., 17)
```

Sample the image in the ordinary way:

```
>>> unique = magni.imaging.measurements.unique_pixels(coords)
>>> samples_normal = img[unique[:, 1], unique[:, 0]]
>>> samples_normal = samples_normal.reshape((len(unique), 1))
```

Sample the image using the present function:

```
>>> from magni.imaging.measurements import construct measurement matrix
>>> matrix = construct measurement matrix(coords, *img.shape)
>>> samples matrix = matrix.dot(vec)
```

Check that the two ways produce the same result:

```
>>> np.allclose(samples matrix, samples normal)
```

magni.imaging.measurements. random line module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

random_line_sample_image(h, w, scan_length, num_points, discrete=None, seed=None) Function for random line sampling an image.

random line sample surface(I, w, speed, sample rate, time, discrete=None, seed=None) Function for random line sampling a surface.

magni.imaging.measurements._random_line.random_line_sample_image(h, w, scan_length, num points, discrete=None, seed=None)

Sample an image using a set of random straight lines.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a pattern based on a set of random straight lines are determined. The scan length determines the length of the path scanned whereas num points indicates the number of samples taken on that path. If discrete is set, it specifies the finite number of equally spaced lines from which the scan lines are be chosen at random. For reproducible results, the *seed* may be used to specify a fixed seed of the random number generator.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - **scan_length** (*float*) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned
 - **discrete** (*int or None, optional*) The number of equally spaced lines from which the scan lines are chosen (the default is None, which implies that no discritisation is used).
 - **seed** (*int or None, optional*) The seed used for the random number generator (the defaul is None, which implies that the random number generator is not seeded).

11/2/16 Page 126 of 209 **Returns:**

coords (ndarray) – The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height h is measured along the y-axis.

Each of the scanned lines span the entire width of the image with the exception of the last line that may only be partially scanned if the *scan_length* implies this. The top and bottom lines of the image are always included in the scan.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import random line sample image
>>> h = 10
>>> w = 10
>>> scan length = 50.0
>>> num_points = 12
>>> seed = 6021
>>> np.set printoptions(suppress=True)
>>> random_line_sample_image(h, w, scan_length, num_points, seed=seed)
array([[ 0.5 , 0.5
    [4.59090909, 0.5
    [8.68181818, 0.5
    [7.01473938, 1.28746666],
   [2.92383029, 1.28746666],
          , 2.95454545],
            , 7.04545455],
    [4.03665944, 7.59970419],
    [8.12756853, 7.59970419],
    [ 8.68181818, 9.5
    [4.59090909, 9.5]
   [ 0.5
          , 9.5 ]])
```

magni.imaging.measurements._random_line.random_line_sample_surface(*I*, *w*, *speed*, *sample_rate*, *time*, *discrete=None*, *seed=None*)

Sample a surface area using a set of random straight lines.

The coordinates (in units of meters) resulting from sampling an image of size/ times w using a pattern based on a set of random straight lines are determined. The scanned path is determined from the probe speed and the scan time. If discrete is set, it specifies the finite number of equally spaced lines from which the scan lines are be chosen at

11/2/16 Page 127 of 209

random. For reproducible results, the seed may be used to specify a fixed seed of the random number generator.

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - **w** (*float*) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample rate** (*float*) The sample rate in units of Hertz.
 - **time** (*float*) The scan time in units of seconds.
 - **discrete** (int or None, optional) The number of equally spaced lines from which the scan lines are chosen (the default is None, which implies that no discritisation is used).
 - **seed** (int or None, optional) The seed used for the random number generator (the defaul is None, which implies that the random number generator is not seeded).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length / is measured along the y-axis.

Each of the scanned lines span the entire width of the image with the exception of the last line that may only be partially scanned if the speed and time implies this. The top and bottom lines of the image are always included in the scan and are not included in the *discrete* number of lines.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import random line sample surface
>>> 1 = 2e-6
>>> w = 2e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> seed = 6021
>>> np.set_printoptions(suppress=True)
>>> random_line_sample_surface(I, w, speed, sample_rate, time, seed=seed)
array([[ 0.
          . 0.
                    ],
    [ 0.00000067, 0.
    [ 0.00000133, 0.
                       ],
    [ 0.000002 , 0.
    [ 0.000002 , 0.00000067]
    [0.000002, 0.00000133],
    [ 0.00000158, 0.00000158],
    [0.00000091, 0.00000158],
    [0.00000024, 0.00000158],
    [ 0.
          , 0.000002 ],
    [0.00000067, 0.000002],
    [0.00000133, 0.000002],
    [0.000002, 0.000002]
```

magni.imaging.measurements. spiral module

Module providing public functions for the magni.imaging.measurements subpackage.

11/2/16 Page 128 of 209 Routine listings

spiral_sample_image(h, w, scan_length, num_points, rect_area=False) Function for spiral sampling an image.

spiral sample surface(I, w, speed, sample rate, time, rect area=False) Function for spiral sampling a surface.

magni.imaging.measurements._spiral_sample_image(h, w, scan_length, num_points, rect_area=False)

Sample an image using an archimedean spiral pattern.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using an archimedean spiral pattern are determined. The scan length determines the length of the path scanned whereas *num_points* indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - scan_length (float) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned
 - rect area (bool) A flag indicating whether or not the full rectangular area is sampled (the default value is False which implies that the "corners" of the rectangular area are not sampled).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height h is measured along the y-axis. The width must equal the height for an archimedian spiral to make sense.

If the rect_area flag is True, then it is assumed that the sampling continues outside of the rectangular area specified by h and w such that the "corners" of the rectangular area are also sampled. The sample points outside of the rectangular area are discarded and, hence, not returned.

Examples

For example,

11/2/16 Page 129 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import spiral_sample_image
>>> w = 10
>>> scan_length = 50.0
>>> num_points = 12
>>> np.set_printoptions(suppress=True)
>>> spiral_sample_image(h, w, scan_length, num_points)
array([[ 6.28776846, 5.17074073],
    [3.13304898, 5.24133767],
    [6.07293751, 2.93873701],
    [6.99638041, 6.80851189],
    [2.89868434, 7.16724999],
    [2.35773914, 3.00320067],
    [6.41495385, 1.71018152],
    [8.82168896, 5.27557847],
    [6.34932919, 8.83624957],
    [ 2.04885699, 8.11199373],
    [0.6196052, 3.96939755]])
```

magni.imaging.measurements._spiral_sample_surface(/, w, speed, sample_rate, time, rect area=False)

Sample a surface area using an archimedean spiral pattern.

The coordinates (in units of meters) resulting from sampling an area of size times w using an archimedean spiral pattern are determined. The scanned path is determined from the probe speed and the scan time.

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - **w** (*float*) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample_rate** (*float*) The sample rate in units of Hertz.
 - **time** (*float*) The scan time in units of seconds.
 - rect area (bool) A flag indicating whether or not the full rectangular area is sampled (the default value is False which implies that the "corners" of the rectangular area are not sampled).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length I is measured along the y-axis. The width must equal the length for an archimedian sprial to make sense.

If the rect area flag is True, then it is assumed that the sampling continues outside of the rectangular area specified by I and w such that the "corners" of the rectangular area are also sampled. The sample points outside of the rectangular area are discarded and, hence, not returned.

Examples

For example,

11/2/16 Page 130 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import spiral_sample_surface
>>> 1 = 1e-6
>>> w = 1e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> spiral_sample_surface(I, w, speed, sample_rate, time)
array([[ 0.00000036, 0.00000046],
    [\ 0.00000052,\ 0.00000071],
    [0.00000052, 0.00000024],
    [0.00000059, 0.00000079],
    [ 0.00000021, 0.00000033], [ 0.00000084, 0.00000036],
    [ 0.00000049, 0.0000009 ],
    [0.0000001, 0.00000036],
    [ 0.00000072, 0.00000011],
    [0.00000089, 0.00000077],
    [ 0.00000021, 0.00000091]])
```

magni.imaging.measurements. square spiral module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
square_spiral_sample_image(h, w, scan_length, num_points)
    Function for square spiral sampling an image.
```

```
square spiral sample surface(I, w, speed, sample rate, time)
    Function for square spiral sampling a surface.
```

magni.imaging.measurements. square spiral. square spiral sample image(h, w, scan length, num points)

Sample an image using a square spiral pattern.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a square spiral pattern are determined. The scan length determines the length of the path scanned whereas *num_points* indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - **scan_length** (*float*) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned path.

Returns:

coords (*ndarray*) – The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height *h* is measured along the y-axis.

Examples

For example,

11/2/16 Page 131 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import square_spiral_sample_image
>>> w = 10
>>> scan_length = 50.0
>>> num_points = 12
>>> np.set_printoptions(suppress=True)
>>> square_spiral_sample_image(h, w, scan_length, num_points)
            , 5.
    [6.28571429, 5.97619048],
    [4.38095238, 3.71428571],
    [ 2.42857143, 5.92857143],
    [ 4.95238095, 7.57142857],
    [7.57142857, 6.02380952],
         , 2.42857143],
    [ 2.83333333, 2.42857143],
    [1.14285714, 4.9047619],
    [ 1.35714286, 8.85714286],
    [5.52380952, 8.85714286],
    [8.85714286, 8.02380952]])
```

magni.imaging.measurements._square_spiral.square_spiral_sample_surface(I, w, speed, sample rate, time)

Sample a surface area using a square spiral pattern.

The coordinates (in units of meters) resulting from sampling an area of size times w using a square spiral pattern are determined. The scanned path is determined from the probe speed and the scan time.

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - w (float) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample_rate** (*float*) The sample rate in units of Hertz.
 - time (float) The scan time in units of seconds.

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length / is measured along the y-axis.

Examples

For example,

11/2/16 Page 132 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import square_spiral_sample_surface
>>> I = 1e-6
>>> w = 1e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> square_spiral_sample_surface(I, w, speed, sample_rate, time)
array([[ 0.0000005, 0.0000005],
    [0.0000004, 0.0000004],
    [0.0000006, 0.0000007],
    [0.0000005, 0.0000003],
    [ 0.0000002, 0.0000007],
    [0.0000008, 0.0000007],
    [0.0000006, 0.0000002],
    [ 0.0000001, 0.0000004],
    0.0000003, 0.0000009]
    [0.0000009, 0.0000008],
    [ 0.0000009, 0.0000001],
    [0.0000002, 0.0000001]])
```

magni.imaging.measurements._uniform_line module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
uniform line sample image(h, w, scan length, num points)
    Function for uniform line sampling an image.
uniform line sample surface(I, w, speed, sample rate, time)
```

Function for uniform line sampling a surface.

magni.imaging.measurements. uniform line.uniform line sample image(h, w, scan length, num points)

Sample an image using a set of uniformly distributed straight lines.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a pattern based on a set of uniformly distributed straight lines are determined. The scan length determines the length of the path scanned whereas num points indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - **scan_length** (*float*) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned path.

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height h is measured along the y-axis.

Each of the scanned lines span the entire width of the image with the exception of the last line that may only be partially scanned if the scan length implies this. The top and bottom lines of the image are always included in the scan.

11/2/16 Page 133 of 209

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import uniform line sample image
>>> w = 10
>>> scan_length = 50.0
>>> num_points = 12
>>> np.set_printoptions(suppress=True)
>>> uniform_line_sample_image(h, w, scan_length, num_points)
array([[ 0.5
           , 0.5
    [4.59090909, 0.5
    [8.68181818, 0.5
    [9.22727273, 3.5]
    [5.13636364, 3.5]
    [1.04545455, 3.5]
    [ 1.04545455, 6.5
    [5.13636364, 6.5
    [ 9.22727273, 6.5
    [8.68181818, 9.5
    [4.59090909, 9.5
    [ 0.5
          , 9.5 ]])
```

magni.imaging.measurements._uniform_line.uniform_line_sample_surface(*I*, *w*, *speed*, *sample_rate*, *time*)

Sample as surface area using a set of uniformly distributed straight lines.

The coordinates (in units of meters) resulting from sampling an area of sizel times w using a pattern based on a set of uniformly distributed straight lines are determined. The scanned path is determined from the probe speed and the scan time.

I (float) – The length of the area to scan in units of meters.
w (float) – The width of the area to scan in units of meters.
speed (float) – The probe speed in units of meters/second.
sample_rate (float) – The sample rate in units of Hertz.
time (float) – The scan time in units of seconds.
coords (ndarray) – The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height I is measured along the y-axis.

Each of the scanned lines span the entire width of the image with the exception of the last line that may only be partially scanned if the *scan_length* implies this. The top and bottom lines of the image are always included in the scan.

Examples

For example,

11/2/16 Page 134 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import uniform_line_sample_surface
>>> 1 = 2e-6
>>> w = 2e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> uniform_line_sample_surface(I, w, speed, sample_rate, time)
             , 0.
array([[ 0.
                     ],
    [ 0.0000067, 0.
    [ 0.00000133, 0.
                        ],
    0.000002 , 0.
    [0.000002, 0.00000067],
    [0.00000167, 0.000001],
    [0.000001, 0.000001],
    [ 0.00000033, 0.000001 ],
           , 0.00000133],
           , 0.000002 ],
    [ 0.
    [ 0.00000067, 0.000002 ],
    [0.00000133, 0.000002],
    [0.000002, 0.000002]])
```

magni.imaging.measurements. uniform rotated line module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
uniform rotated line sample image(h, w, scan length, num points, angle=0.,
   follow edge=True) Function for uniform rotated line sampling an image.
```

uniform rotated line sample surface(I, w, speed, sample rate, time, angle=0., follow_edge=True) Function for uniform rotated line sampling a surface.

scan length, num points, angle=0.0, follow edge=True)

Sample an image using a uniform rotated line pattern.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a uniform rotated line pattern are determined. The scan length determines the length of the path scanned whereas *num points* indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - **scan_length** (*float*) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned
 - angle (float) The angle measured in radians by which the uniform lines are rotated (the default is 0.0 resulting in a pattern identical to that of uniform line sample image).
 - follow edge (bool) A flag indicating whether or not the pattern follows the edges of the rectangular area in-between lines (the default is True).

Returns:

coords (*ndarray*) – The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

11/2/16 Page 135 of 209

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height h is measured along the y-axis.

The angle is limited to the interval

-). An angle of 0 results in the same behaviour as that of uniform_line_sample_image. An increase in the angle rotates the overall direction counterclockwise, i.e., at
- $\frac{1}{2}$, the uniform_line_sample_image sampling pattern is rotated 90 degrees counterclockwise.

If the *follow_edge* flag is True, then the pattern follows the edges of the rectangular area when moving from one line to the next. If the flag is False, then the pattern follows a line perpendicular to the uniform lines when moving from one line to the next. In the latter case, some of the uniform lines are shortened to allow the described behaviour.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import \
... uniform rotated line sample image
>>> h = 10
>>> w = 10
>>> scan_length = 50.0
>>> num_points = 12
>>> np.set_printoptions(suppress=True)
>>> uniform_rotated_line_sample_image(h, w, scan_length, num_points)
array([[ 0.5
             , 0.5
    [4.59090909, 0.5
    [8.68181818, 0.5
    [ 9.22727273, 3.5
    [5.13636364, 3.5]
    [1.04545455, 3.5]
    [1.04545455, 6.5]
     5.13636364, 6.5
     9.22727273, 6.5
    [8.68181818, 9.5]
    [4.59090909, 9.5]
    [ 0.5
           , 9.5
```

magni.imaging.measurements._uniform_rotated_line.uniform_rotated_line_sample_surface(*l*, *w*, speed, sample_rate, time, angle=0.0, follow_edge=True)

Sample a surface area using a uniform rotated line pattern.

The coordinates (in units of meters) resulting from sampling an area of size l times l using uniform rotated line pattern are determined. The scanned path is determined from the probe l and the scan l time.

11/2/16 Page 136 of 209

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - w (float) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample_rate** (*float*) The sample rate in units of Hertz.
 - **time** (*float*) The scan time in units of seconds.
 - angle (float) The angle measured in radians by which the uniform lines are rotated (the default is 0.0 resulting in a pattern identical to that of uniform line sample image).
 - follow_edge (bool) A flag indicating whether or not the pattern foolows the edges of the rectangular area in-between lines (the default is True).

Returns:

coords (*ndarray*) – The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length / is measured along the y-axis.

The *angle* is limited to the interval

-). An angle of 0 results in the same behaviour that uniform line sample surface. An increase in the angle rotates the overall direction counterclockwise, i.e., at
- $\frac{1}{2}$, the uniform_line_sample_surface sampling pattern is rotated 90 degrees counterclockwise.

If the follow edge flag is True, then the pattern follows the edges of the rectangular area when moving from one line to the next. If the flag is False, then the pattern follows a line perpendicular to the uniform lines when moving from one line to the next. In the latter case, some of the uniform lines are shortened to allow the described behaviour.

Examples

For example,

11/2/16 Page 137 of 209

```
>>> import numpy as np
>>> from magni.imaging.measurements import \
... uniform rotated line sample surface
>>> 1 = 1e-6
>>> w = 1e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> uniform_rotated_line_sample_surface(I, w, speed, sample_rate, time)
              , 0.
    [ 0.0000067, 0.
    [ 0.00000083, 0.00000017], [ 0.00000017, 0.00000017],
    [0.00000033, 0.00000033],
    [ 0.000001 , 0.00000033], [ 0.0000005 , 0.0000005 ],
            , 0.0000067],
    [ 0.00000067, 0.00000067],
    [ 0.00000083, 0.00000083],
    [0.00000017, 0.00000083],
    [ 0.00000033, 0.000001 ],
    [ 0.000001 , 0.000001 ]])
```

magni.imaging.measurements._util module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
construct pixel mask(h, w, pixels)
    Construct a binary pixel mask.
```

unique pixels(coords)

Function for determining unique pixels from a set of coordinates.

```
magni.imaging.measurements. util. construct pixel mask(h, w, pixels)
```

Construct a binary pixel mask.

An image (2D array) of shape w x h is created where all pixels are marked True.

- **Parameters:** h (*int*) The height of the image in pixels.
 - **w** (*int*) The width of the image in pixels.
 - pixels (ndarray) The 2D array of pixels that make up the mask. Each row is a coordinate pair (x, y), such that coords has size len(pixels) x 2.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import construct_pixel_mask
>>> h = 3
>>> w = 3
>>> pixels = np.array([[0, 0], [1, 1], [2, 1]])
>>> construct_pixel_mask(h, w, pixels)
array([[ True, False, False],
    [False, True, True],
    [False, False, False]], dtype=bool)
```

11/2/16 Page 138 of 209 magni.imaging.measurements. util. **unique pixels** (coords)

Identify unique pixels from a set of coordinates.

The floating point coords are reduced to a unique set of integer pixels by flooring the floating point values.

Parameters: coords (*ndarray*) – The *k* floating point coordinates arranged into a 2D

array where each row is a coordinate pair (x, y), such that *coords* has

size $k \times 2$.

unique pixels (ndarray) – The $l \le k$ unique (integer) pixels, such that **Returns:**

unique pixels is a 2D array and has size / x 2.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import unique pixels
>>> coords = np.array([[1.7, 1.0], [1.0, 1.2], [3.3, 4.3]])
>>> np.int_(unique_pixels(coords))
array([[1, 1],
    [3, 4]])
```

magni.imaging.measurements. visualisation module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

plot_pattern(I, w, coords, mode, output_path=None)

Function for visualising a scan pattern.

plot pixel mask(h, w, pixels, output path=None)

Function for visualising a pixel mask obtained from a scan pattern.

magni.imaging.measurements._visualisation. plot_pattern(I, w, coords, mode, output_path=None) Display a plot that shows the pattern given by a set of coordinates.

The pattern given by the *coords* is displayed on an w x I area. If mode is 'surface', I and w are regarded as measured in meters. If mode is 'image', I and w are regarded as measured in pixels. The coords are marked by filled circles and connected by straight dashed lines.

- **Parameters:** I (float or int) The length/height of the area. If mode is 'surface', it must be a float. If *mode* is 'image', it must be an integer.
 - w (float or int) The width of the area. If mode is 'surface', it must be a float. If *mode* is 'image', it must be an integer.
 - coords (ndarray) The 2D array of pixels that make up the mask. Each row is a coordinate pair (x, y).
 - **mode** ({'surface', 'image'}) The display mode that dertermines the axis labeling and the type of I and w.
 - **output path** (*str. optional*) Path (including file type extension) under which the plot is saved (the default value is None which implies, that the plot is not saved).

11/2/16 Page 139 of 209

Notes

The resulting plot is displayed in a figure using matplotlib's pyplot.plot.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import plot pattern
>>> I = 3
>>> W = 3
>>> coords = np.array([[0, 0], [1, 1], [2, 1]])
>>> mode = 'image'
>>> plot_pattern(I, w, coords, mode)
```

magni.imaging.measurements. visualisation. **plot pixel mask**(h, w, pixels, output path=None)

Display a binary image that shows the given pixel mask.

A black image with w x h pixels is created and the pixels are marked with white.

- **Parameters:** h (*int*) The height of the image in pixels.
 - **w** (*int*) The width of the image in pixels.
 - pixels (ndarray) The 2D array of pixels that make up the mask. Each row is a coordinate pair (x, y), such that coords has size len(*pixels*) x 2.
 - **output_path** (*str, optional*) Path (including file type extension) under which the plot is saved (the default value is None which implies, that the plot is not saved).

Notes

The resulting image is displayed in a figure using magni.imaging.visualisation.imshow.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import plot_pixel_mask
>>> h = 3
>>> W = 3
>>> pixels = np.array([[0, 0], [1, 1], [2, 1]])
>>> plot_pixel_mask(h, w, pixels)
```

magni.imaging.measurements. zigzag module

Module providing public functions for the magni.imaging.measurements subpackage.

Routine listings

```
zigzag sample image(h, w, scan length, num points, angle=np.pi / 20)
    Function for zigzag sampling an image.
```

zigzag sample surface(I, w, speed, sample rate, time, angle=np.pi / 20) Function for zigzag sampling a surface.

11/2/16 Page 140 of 209 magni.imaging.measurements._zigzag. zigzag_sample_image(h, w, scan_length, num_points, angle=0.15707963267948966)

Sample an image using a zigzag pattern.

The coordinates (in units of pixels) resulting from sampling an image of sizeh times w using a zigzag pattern are determined. The scan length determines the length of the path scanned whereas *num points* indicates the number of samples taken on that path.

- **Parameters:** h (*int*) The height of the area to scan in units of pixels.
 - **w** (*int*) The width of the area to scan in units of pixels.
 - scan length (*float*) The length of the path to scan in units of pixels.
 - **num points** (*int*) The number of samples to take on the scanned path.
 - angle (float) The angle measured in radians by which the lines deviate from being horizontal (the default is pi / 20).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the height *h* is measured along the y-axis.

The angle is measured clockwise relative to horizontal and is limited to the interval

```
ct n \frac{h}{w}
```

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import zigzag_sample_image
>>> h = 10
>>> w = 10
>>> scan_length = 50.0
>>> num points = 12
>>> np.set_printoptions(suppress=True)
>>> zigzag_sample_image(h, w, scan_length, num_points)
array([[ 0.5
              , 0.5
    [ 4.98949246, 1.21106575],
    [9.47898491, 1.9221315],
    [5.03152263, 2.63319725],
    [ 0.54203017, 3.344263 ],
    [ 4.94746229, 4.05532875],
    [ 9.43695474, 4.7663945 ],
    [5.0735528, 5.47746025],
    [ 0.58406034, 6.188526 ],
    [ 4.90543212, 6.89959175],
[ 9.39492457, 7.6106575 ],
[ 5.11558297, 8.32172325]])
```

magni.imaging.measurements._zigzag. zigzag_sample_surface(/, w, speed, sample_rate, time, angle=0.15707963267948966)

Sample a surface area using a zigzag pattern.

The coordinates (in units of meters) resulting from sampling an area of size times w

11/2/16 Page 141 of 209 using a zigzag pattern are determined. The scanned path is determined from the probe speed and the scan time.

- **Parameters:** I (*float*) The length of the area to scan in units of meters.
 - w (float) The width of the area to scan in units of meters.
 - **speed** (*float*) The probe speed in units of meters/second.
 - **sample rate** (*float*) The sample rate in units of Hertz.
 - **time** (*float*) The scan time in units of seconds.
 - angle (float) The angle measured in radians by which the lines deviate from being horizontal (the default is pi / 20).

Returns:

coords (ndarray) - The coordinates of the samples arranged into a 2D array, such that each row is a coordinate pair (x, y).

Notes

The orientation of the coordinate system is such that the width w is measured along the x-axis whereas the length / is measured along the y-axis.

The *angle* is measured clockwise relative to horizontal and is limited to the interval

ct n
$$\frac{h}{w}$$

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.measurements import zigzag_sample_surface
>>> 1 = 1e-6
>>> w = 1e-6
>>> speed = 7e-7
>>> sample_rate = 1.0
>>> time = 12.0
>>> np.set_printoptions(suppress=True)
>>> zigzag_sample_surface(l, w, speed, sample_rate, time)
array([[ 0.
            , 0.
    [ 0.00000069, 0.00000011],
    [ 0.0000062, 0.00000022],
    [0.00000007, 0.00000033],
    [0.00000077, 0.00000044],
    [0.00000054, 0.00000055],
    [ 0.00000015, 0.00000066],
    [0.00000084, 0.00000077],
    [0.00000047, 0.00000088],
    [0.00000022, 0.00000099],
    [ 0.00000091, 0.0000009 ],
    [0.00000039, 0.0000008],
    [ 0.0000003 , 0.00000069]])
```

Submodules

magni.imaging. util module

Module providing the public functions of the magni.imaging subpackage.

```
magni.imaging. util. double mirror(img, fftstyle=False)
```

Mirror image in both the vertical and horisontal axes.

11/2/16 Page 142 of 209

The image is mirrored around its upper left corner first in the horizontal axis and then in the vertical axis such that an image of four times the size of the original is returned. If fftstyle is True, the image is constructed such it would represent a fftshifted version of the mirrored *img* such that entry (0, 0) is the DC component.

- **Parameters:** img (*ndarray*) The image to mirror.
 - **fftstyle** (bool) The flag that indicates if the fftstyle mirrored image is returned.

Returns:

mirrored img (*ndarray*) – The mirrored image.

Examples

For example, mirror a very simple 2-by-3 pixel image.

```
>>> import numpy as np
>>> from magni.imaging._util import double_mirror
>>> img = np.arange(6).reshape(2, 3)
>>> img
array([[0, 1, 2],
    [3, 4, 5]])
>>> double_mirror(img)
array([[5, 4, 3, 3, 4, 5],
    [2, 1, 0, 0, 1, 2],
    [2, 1, 0, 0, 1, 2],
    [5, 4, 3, 3, 4, 5]]
>>> double_mirror(img, fftstyle=True)
array([[0, 0, 0, 0, 0, 0],
    [0, 5, 4, 3, 4, 5],
    [0, 2, 1, 0, 1, 2],
    [0, 5, 4, 3, 4, 5]]
```

magni.imaging. util. **get inscribed masks** (*img*, *as vec=False*)

Return a set of inscribed masks covering the image.

Two masks are returned. One is the disc with radius equal to that of the inscribed circle for img. The other is the inscribed square of the first mask. If as vec is True, the img must be a vector representation of the (matrix) image. In this case, the masks are also returned in vector representation.

- **Parameters:** img (ndarray) The square image of even height/width which the masks should cover.
 - as_vec (bool) The indicator of whether or not to treating as a vector instead of an image (the default is False, which implies that img is treated as a matrix.

Returns:

- cicle_mask (ndarray) The inscribed cicle mask.
- **square_mask** (*ndarray*) The inscribed square mask.

Examples

For example, get the inscribed masks of an 8-by-8 image:

11/2/16 Page 143 of 209

```
>>> import numpy as np
>>> from magni.imaging._util import get_inscribed_masks, mat2vec
>>> img = np.arange(64).reshape(8, 8)
>>> circle_mask, square_mask = get_inscribed_masks(img)
>>> np_printoptions = np.get_printoptions()
>>> np.set_printoptions(formatter={'bool': lambda x: str(int(x))})
>>> circle_mask
array([[0, 1, 1, 1, 1, 1, 1, 0],
      [1, 1, 1, 1, 1, 1, 1, 1],
      [1, 1, 1, 1, 1, 1, 1, 1],
      [1, 1, 1, 1, 1, 1, 1, 1],
      [1, 1, 1, 1, 1, 1, 1, 1, 1],
[1, 1, 1, 1, 1, 1, 1, 1],
      [1, 1, 1, 1, 1, 1, 1, 1],
[0, 1, 1, 1, 1, 1, 1, 0]], dtype=bool)
>>> square_mask
array([[0, 0, 0, 0, 0, 0, 0, 0],
      [0, 1, 1, 1, 1, 1, 1, 0],

[0, 1, 1, 1, 1, 1, 1, 0],

[0, 1, 1, 1, 1, 1, 1, 0],

[0, 1, 1, 1, 1, 1, 1, 0],

[0, 1, 1, 1, 1, 1, 1, 0],

[0, 1, 1, 1, 1, 1, 1, 0],
      [0, 0, 0, 0, 0, 0, 0, 0]], dtype=bool)
```

Or get the same masks based on a vector:

11/2/16 Page 144 of 209

```
>>> img_vec = mat2vec(img)
>>> c_vec_mask, s_vec_mask = get_inscribed_masks(img_vec, as_vec=True)
>>> c_vec_mask
array([[0],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [0],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
[1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [0],
     [1],
     [1],
     [1],
     [1],
     [1],
     [1],
     [0]], dtype=bool)
>>> np.set_printoptions(**np_printoptions)
```

11/2/16 Page 145 of 209

```
magni.imaging. util. mat2vec(x)
```

Reshape x from matrix to vector by stacking columns.

Parameters: x (*ndarray*) – Matrix that should be reshaped to vector. Returns: ndarray – Column vector formed by stacking the columns of the matrix Χ.

See also:

vec2mat()

The inverse operation

Notes

The returned column vector is C contiguous.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging._util import mat2vec
>>> x = np.arange(4).reshape(2, 2)
>>> X
array([[0, 1],
    [2, 3]])
>>> mat2vec(x)
array([[0],
    [2],
    [1],
    [3]])
```

magni.imaging._util. **vec2mat**(*x*, *mn_tuple*)

Reshape *x* from column vector to matrix.

- **Parameters: x** (*ndarray*) Matrix that should be reshaped to vector.
 - mn_tuple (tuple) A tuple (m, n) containing the parameters m, n as listed below.
 - **m** (*int*) Number of rows in the resulting matrix.
 - **n** (*int*) Number of columns in the resulting matrix.

Returns:

ndarray - Matrix formed by taking n columns of lenght m from the column vector x.

See also:

mat2vec()

The inverse operation

Notes

The returned matrix is C contiguous.

11/2/16 Page 146 of 209

Examples

For example,

magni.imaging.domains module

Module providing a multi domain image class.

Routine listings

MultiDomainImage(object)

Provide access to an image in the domains of a compressed sensing context.

class magni.imaging.domains. MultiDomainImage(Phi, Psi)

Bases: object

Provide access to an image in the domains of a compressed sensing context.

Given a measurement matrix and a dictionary, an image can be supplied in either the measurement domain, the image domain, or the coefficient domain. This class then provides access to the image in all three domains.

```
    Parameters:
    Phi
        magni.utils.matrices.MatrixCollection,)
        measurement matrix.
    Psi
        magni.utils.matrices.MatrixCollection,)
        magni.utils.matrices.Matrix,
        magni.utils.matrices.Matrix
```

Notes

The image is only converted to other domains than the supplied when the the image is requested in another domain. The image is, however, stored in up to three versions internally in order to reduce computation overhead. This may introduce a memory overhead.

Examples

Define a measurement matrix which skips every other sample:

```
>>> import numpy as np, magni
>>> func = lambda vec: vec[::2]
>>> func_T = lambda vec: np.float64([vec[0], 0, vec[1]]).reshape(3, 1)
>>> Phi = magni.utils.matrices.Matrix(func, func_T, (), (2, 3))
```

11/2/16 Page 147 of 209

Define a dictionary which is simply a rotated identity matrix:

Instantiate the current class to handle domains:

```
>>> from magni.imaging.domains import MultiDomainImage
>>> domains = MultiDomainImage(Phi, Psi)
```

An image can the be supplied in any domain and likewise retrieved in any domain. For example, setting the measurements and getting the coefficients:

Or setting the coefficients and getting the image:

Or setting the image and getting the measurements:

```
>>> domains.image = np.ones(3).reshape(-1, 1)
>>> print(domains.measurements)
[[ 1.]
    [ 1.]]
```

coefficients

Get the image in the coefficient domain.

Returns: coefficients (numpy.ndarray) – The dictionary coefficients of the image.

image

Get the image in the image domain.

Returns: image (*numpy.ndarray*) – The image.

measurements

Get the image in the measurement domain.

Returns: measurements (numpy.ndarray) – The measurements of the image.

magni.imaging.evaluation module

Module providing functions for evaluation of image reconstruction quality.

11/2/16 Page 148 of 209

Routine listings

calculate_mse(x_org, x_recons)

Function to calcualte Mean Squared Error (MSE).

calculate psnr(x org, x recons, peak)

Function to calculate Peak Signal to Noise Ratio (PSNR).

calculate_retained_energy(x_org, x_recons)

Function to calculate the percentage of energy retained in reconstruction.

magni.imaging.evaluation. **calculate** $mse(x \ org, x \ recons)$

Calculate Mean Squared Error (MSE) between x recons and x org.

Parameters: • **x_org** (*ndarray*) – Array of original values.

• x_recons (ndarray) - Array of reconstruction

values.

Returns: mse (float) – Mean Squared Error (MSE).

Notes

The Mean Squared Error (MSE) is calculated as:

$$\frac{1}{-}\cdot\sum(x_o - x_{cons})^2$$

where *N* is the number of entries in *x_org*.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.evaluation import calculate_mse
>>> x_org = np.arange(4).reshape(2, 2)
>>> x_recons = np.ones((2,2))
>>> print('{:.2f}'.format(calculate_mse(x_org, x_recons)))
1.50
```

magni.imaging.evaluation. **calculate_psnr**(*x_org*, *x_recons*, *peak*)

Calculate Peak Signal to Noise Ratio (PSNR) between*x_recons* and *x_org*.

Parameters: • **x_org** (*ndarray*) – Array of original values.

• **x_recons** (*ndarray*) – Array of reconstruction values.

• peak (int or float) - Peak value.

Returns: psnr (float) – Peak Signal to Noise Ratio (PSNR) in

dB.

Notes

The PSNR is as calculated as

11/2/16 Page 149 of 209

$$1 \cdot \circ {}_1 \quad \frac{k^2}{1/ \quad \cdot \quad (x_o \quad -x \quad {}_{cons})^2}$$

where N is the number of entries in x org.

lf

$$|x_{o} - x_{cons}| = (1 + 1 |x_{cons}|)$$
 then `np.inf` is returned.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.evaluation import calculate psnr
>>> x_org = np.arange(4).reshape(2, 2)
>>> x_recons = np.ones((2,2))
>>> peak = 3
>>> print('{:.2f}'.format(calculate_psnr(x_org, x_recons, peak)))
7.78
```

magni.imaging.evaluation.calculate retained energy (x org, x recons)

Calculate percentage of energy retained in reconstruction.

- Parameters: x_org (ndarray) Array of original values (must not be all zeros).
 - x recons (ndarray) Array of reconstruction values.

Returns:

energy (*float*) – Percentage of retained energy in reconstruction.

Notes

The retained energy is as calculated as

$$rac{x^2_{cons}}{x_o^2} \cdot 1$$

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.evaluation import calculate_retained_energy
>>> x_org = np.arange(4).reshape(2, 2)
>>> x_recons = np.ones((2,2))
>>> print('{:.2f}'.format(calculate_retained_energy(x_org, x_recons)))
28.57
```

magni.imaging.preprocessing module

Module providing functionality to remove tilt in images.

Routine listings

detilt(img, mask=None, mode='plane_flatten', degree=1, return_tilt=False) Function to remove tilt from an image.

11/2/16 Page 150 of 209 magni.imaging.preprocessing. detilt(img, mask=None, mode='plane flatten', degree=1, return tilt=False)

Estimate the tilt in an image and return the detilted image.

- **Parameters:** img (*ndarray*) The image that is to be detilted.
 - mask (ndarray, optional) Bool array of the same size asimg indicating the pixels to use in detilt (the default is None, which implies, that the the entire image is used)
 - mode ({'line_flatten', 'plane_flatten'}, optional) The type of detilting applied (the default is plane flatten).
 - degree (int, optional) The degree of the polynomial used in line flattening (the default is 1).
 - return_tilt (bool, optional) If True, the detilted image and the estimated tilt is returned (the default is False).

Returns:

- **img_detilt** (*ndarray*) Detilted image.
- tilt (ndarray, optional) The estimated tilt (image). Only returned if return tilt is True.

Notes

If mode is line flatten, the tilt in each horizontal line of pixels in the image is estimated by a polynomial fit independently of all other lines. If mode is plane flatten, the tilt is estimated by fitting a plane to all pixels.

If a custom *mask* is specified, only the masked (True) pixels are used in the estimation of the tilt.

Examples

For example, line flatten an image using a degree 1 polynomial

```
>>> import numpy as np
>>> from magni.imaging.preprocessing import detilt
>>> img = np.array([[0, 2, 3], [1, 5, 7], [3, 6, 8]], dtype=np.float)
>>> np.set printoptions(suppress=True)
>>> detilt(img, mode='line_flatten', degree=1)
array([[-0.16666667, 0.33333333, -0.16666667],
    [-0.33333333, 0.66666667, -0.33333333],
    [-0.16666667, 0.33333333, -0.16666667]])
```

Or plane flatten the image based on a mask and return the tilt

```
>>> mask = np.array([[1, 0, 0], [1, 0, 1], [0, 1, 1]], dtype=np.bool)
>>> im, ti = detilt(img, mask=mask, mode='plane_flatten', return_tilt=True)
>>> np.set_printoptions(suppress=True)
array([[ 0.11111111, -0.66666667, -2.44444444],
    [-0.33333333, 0.88888889, 0.11111111]
    [ 0.2222222, 0.4444444, -0.33333333]])
array([[-0.11111111, 2.66666667, 5.44444444],
    [1.33333333, 4.11111111, 6.88888889],
    [ 2.77777778, 5.55555556, 8.333333333]])
```

magni.imaging.preprocessing. **line_flatten_tilt**(*img*, *mask*, *degree*)

Estimate tilt using the line flatten method.

11/2/16 Page 151 of 209

- **Parameters:** img (ndarray) The image from which the tilt is estimated.
 - mask (ndarray, or None) If not None, a bool ndarray of the the shape as img indicating which pixels should be used in estimate of tilt.
 - **degree** (*int*) The degree of the polynomial in the estimated line tilt.

Returns:

tilt (ndarray) - The estimated tilt.

magni.imaging.preprocessing. **plane_flatten_tilt**(img, mask)

Estimate tilt using the plane flatten method.

- **Parameters:** img (ndarray) The image from which the tilt is estimated.
 - mask (ndarray, or None) If not None, a bool ndarray of the the shape as img indicating which pixels should be used in estimate of tilt.

Returns:

tilt (*ndarray*) – The estimated tilt.

magni.imaging.visualisation module

Module providing functionality for visualising images.

The module provides functionality for adjusting the intensity of an image. It provides a wrapper of the matplotlib.pyplot.imshow function that may exploit the provided functions for adjusting the image intensity. Also it include a function may be used to display a set of related images using a common colormapping.

Routine listings

imshow(X, ax=None, intensity_func=None, intensity_args=(), **kwargs)

Function that may be used to display an image.

imsubplot(imgs, rows, titles=None, x_labels=None, y_labels=None,

y ticklabels=None, cbar label=None, normalise=True. x ticklabels=None, fixed_clim=None, **kwargs) Function that may be used to display a set of related images.

mask img from coords(img, coords)

Function for masking certain parts of an image based on coordinates.

shift mean(x mod, x org)

Function for shifting mean intensity of an image based on another image.

stretch image(img, max val, min val=0)

Function for stretching the intensity of an image.

magni.imaging.visualisation. **imshow** (X, ax=None, intensity func=None, intensity args=(), show axis='frame', **kwargs)

Display an image.

Wrap matplotlib.pyplot.imshow to display a possibly intensity manipulated version of the image X.

11/2/16 Page 152 of 209

- **Parameters: X** (*ndarray*) The image to be displayed.
 - ax (matplotlib.axes.Axes, optional) The axes on which the image is displayed (the default is None, which implies that the current axes is used).
 - **intensity func** (*FunctionType*, *optional*) The handle to the function used to manipulate the image intensity before the image is displayed (the default is None, which implies that no intensity manipulation is used).
 - intensity_args (list or tuple, optional) The arguments that are passed to the intensity func (the default is (), which implies that no arguments are passed).
 - show axis ({'none', 'top', 'inherit', 'frame'}) How the x- and y-axis are display. If 'none', no axis are displayed. If 'top', the x-axis is displayed at the top of the image. If 'inherit', the axis display is inherited from matplotlib.pyplot.imshow. If 'frame' only the frame is shown and not the ticks.

Returns:

im_out (*matplotlib.image.AxesImage*) – The AxesImage returned by matplotlibs imshow.

See also:

matplotlib.pyplot.imshow()

Matplotlib's imshow function.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.visualisation import imshow
>>> X = np.arange(4).reshape(2, 2)
>>> add_k = lambda X, k: X + k
>>> im_out = imshow(X, intensity_func=add_k, intensity_args=(2,))
```

magni.imaging.visualisation. imsubplot(imgs, rows, titles=None, x labels=None, y labels=None, x ticklabels=None, y ticklabels=None, cbar label=None, normalise=True, fixed clim=None, **kwarqs)

Display a set of related images as subplots in a figure.

The images images are shown in a figure with a subplot layout based on the number of rows. The titles, x_labels, y_labels, x_ticklabels, and y_ticklabels are shown in the If normalise is True, all the images will share the same normalised colorbar/colormapping i.e. a particular colour will correspond to the same value across all images.

11/2/16 Page 153 of 209

- **Parameters:** imgs (*list or tuple*) The images (as ndarrays) that is to be displayed.
 - rows (int) The number of rows to use in the subplot layout.
 - titles (list or tuple) The titles (as strings) to use for each of the subplots (the default is None, which implies that no titles are displayed).
 - x_labels (list or tuple) The x_labels (as strings) to use for each of the subplots (the default is None, which implies that no x labels are displayed).
 - y_labels (list or tuple) The y_labels (as strings) to use for each of the subplots (the default is None, which implies that no y labels are displayed).
 - x ticklabels (list or tuple) The x ticklabels (as strings or lists of strings) to use for the subplots (the default is None, which implies that no x ticklabels are displayed).
 - y_ticklabels (list or tuple) The y ticklabels (as strings or lists of strings) to use for the subplots (the default is None, which implies that no y ticklabels are displayed).
 - cbar label (str) The colorbar label to use with a normalised colormapping (the default is None, which implies that no colorbar label is displayed).
 - fixed clim (list or tuple) The colorbar limits as a (min, max) sequence (the default is None, which implies that the colorbar limits are inferred from the data).
 - normalise (bool) The flag that indicates whether to use a normalised colormapping.

Returns:

fig (*matplotlib.figure.Figure*) – The resulting figure instance.

See also:

matplotlib.pyplot.subplots()

Underlying subplot function.

Notes

The x_ticklabels and y_ticklabels may be either a collection of strings or a collections of collections of strings depending on wherther the labels should be shared across all subplots or different labels are to be used for each subplot.

Additional kwargs given to the function will be passed to the underlying suplot instantiation function matplotlib.pyplot.subplots.

If *normalise* is True, the common colorbar is shown below the subplots.

The implementation of the normalisation feature is based on the Pylab example: http://matplotlib.org/examples/pylab examples/multi image.html.

Examples

For example, show two images next to each other with a common colormapping:

11/2/16 Page 154 of 209

```
>>> import numpy as np
>>> from magni.imaging.visualisation import imsubplot
>>> img1 = np.arange(4).reshape(2, 2)
>>> img2 = np.ones((4, 4))
>>> fig = imsubplot([img1, img2], 1, titles=['arange', 'ones'],
... x_labels=['x1', 'x2'], y_labels=['y1', 'y2'], cbar_label='Example',
... normalise=True)
```

or show the same images with shared ticklabels and fixed colorbar limits:

```
>>> common_x_ticklabels = ['a', 'b']
>>> common_y_ticklabels = ['c', 'd']
>>> fig = imsubplot([img1, img2], 1, x_ticklabels=common_x_ticklabels,
... y_ticklabels=common_y_ticklabels, fixed_clim=(0, 2), normalise=True)
```

or show the same images with different colormappings and ticklabels:

```
>>> fig = imsubplot([img1, img2], 1, x_ticklabels=x_ticklabels,
... y ticklabels=y ticklabels, normalise=False)
```

magni.imaging.visualisation. mask_img_from_coords(img, coords)

Mask coordinates in an image.

The coordinates coords in the image img are masked such that only the cooordinates are shown.

- **Parameters:** img (ndarray) The image to mask.
 - coords (ndarray) The coordinates arranged into a 2D array, such that each row is a coordinate pair (x, y).

Returns:

masked_img (*numpy.ma.MaskedArray*) – The masked image.

See also:

magni.imaging.measurements()

Further description of the coordinate format.

Examples

For example, display only center pixel in a 3-by-3 image

```
>>> import numpy as np
>>> from magni.imaging.visualisation import mask_img_from_coords
>>> img = np.arange(9).reshape(3, 3)
>>> coords = np.array([[1, 1]])
>>> mask img from coords(img, coords)
masked_array(data =
[[-- -- --]
[-- 4 --]
[-- -- --]],
       mask =
[[ True True True]
[True False True]
[True True True]]
    fill value = 999999)
```

magni.imaging.visualisation. **shift_mean**(x_mod, x_org)

11/2/16 Page 155 of 209 Shift the mean value of x mod such that it equals the mean of x org.

Parameters: • **x_org** (*ndarray*) – The array which hold the "true" mean value. • x mod (ndarray) - The modified copy of x org which must have its mean value shifted. **Returns: shifted x mod** (ndarray) – A copy of x mod with the same mean value as x_org.

Examples

For example,

```
>>> import numpy as np
>>> from magni.imaging.visualisation import shift mean
>>> x_org = np.arange(4).reshape(2, 2)
>>> x \mod = np.ones((2, 2))
>>> print('{:.1f}'.format(x_org.mean()))
>>> print('{:.1f}'.format(x mod.mean()))
>>> shifted x mod = shift mean(x mod, x org)
>>> print('{:.1f}'.format(shifted_x_mod.mean()))
>>> np.set_printoptions(suppress=True)
>>> shifted x mod
array([[ 1.5, 1.5],
    [ 1.5, 1.5]])
```

magni.imaging.visualisation. **stretch_image**(*img*, *max_val*, *min_val=0*)

Stretch image such that pixels values are in [min val, max val].

- **Parameters:** img (*ndarray*) The (float) image that is to be stretched.
 - max val (int or float) The maximum value in the stretched image.
 - min_val (int or float) The minimum value in the stretched image.

Returns:

stretched_img (*ndarray*) – A stretched copy of the input image.

Notes

The pixel values in the input image are scaled to lie in the interval finin val, max val using a linear stretch.

Examples

For example, stretch an image between 0 and 1

```
>>> import numpy as np
>>> from magni.imaging.visualisation import stretch image
>>> img = np.arange(4, dtype=np.float).reshape(2, 2)
>>> stretched_img = stretch_image(img, 1)
>>> np.set_printoptions(suppress=True)
>>> stretched img
array([[ 0.
             , 0.33333333],
    [ 0.6666667, 1.
```

or stretch the image between -1 and 1

11/2/16 Page 156 of 209

```
>>> stretched_img = stretch_image(img, 1.0, min_val=-1.0)
>>> stretched_img
array([[-1. , -0.33333333],
    [ 0.33333333, 1.
```

or re-stretch the strecthed image between -3.0 and -1.5

```
>>> stretched_img = stretch_image(stretched_img, -1.5, min_val=-3.0)
>>> stretched img
array([[-3., -2.5],
    [-2., -1.5]])
```

or re-stretch that image between 1.25 and 8.00

```
>>> stretched_img = stretch_image(stretched_img, 8.00, min_val=1.25)
>>> stretched img
array([[ 1.25, 3.5],
    [5.75, 8.]])
```

magni.imaging.visualisation. **handle_ticklabels**(ax, k, x_ticklabels, y_ticklabels)

Handle and format ticks and ticklabels for use in imsubplot.

The imsubplot function creates a figure showing an abitrary number of subplots. The imsubplot function allows for custom ticklabels along the x- and y-axes. This function handles the formatting of the ticklabels for a given subplot.

- **Parameters:** ax (matplotlib.axes.Axes) The matplotlib axes (subplot) to format ticklabels for.
 - **k** (int) The axes index, i.e. the subplot number out of the total number of subplots.
 - x_ticklabels (list or tuple) The x_ticklabels (as strings or lists of strings) to use for the subplots (the default is None, which implies that no x ticklabels are displayed).
 - y ticklabels (list or tuple) The y ticklabels (as strings or lists of strings) to use for the subplots (the default is None, which implies that no y ticklabels are displayed).

class magni.imaging.visualisation. _ ImageColourTracker(tracker)

Bases: object

Track a common 'clim' in a set of matplotlib image subplots.

Parameters: tracker (*matplotlib.image.AxesImage*) – The image instance that must track a given 'clim'.

```
__call__(tracked)
```

magni.reproducibility package

Module providing functionality for aiding in quest for more reproducible research.

Routine listings

data

Module providing functions that return various data about the system, magni, files, etc.

11/2/16 Page 157 of 209 io

Module providing input/output functions to databases containing results from reproducible research.

Submodules

magni.reproducibility._annotation module

Module providing functions that may be used to annotate data.

Routine Listings

get conda info()

Function that returns information about an Anaconda install.

get_datetime()

Function that returns information about the current date and time.

get_git_revision(git_root_dir=None)

Function that returns information about the current git revision.

get_file_hashes(path, blocksize=2**30)

Function that returns the md5 and sha256 checksums of a file.

get_magni_config()

Function that returns information about the current configuration of Magni.

get magni info()

Function that returns genral information about Magni.

get_platform_info()

Function that returns information about the platform used to run the code.

Notes

The returned annotations are any nested level of dicts of dicts of strings.

magni.reproducibility._annotation.get_conda_info()

Return a dictionary contianing information from Conda.

Conda is the package manager for the Anaconda scientific Python distribution. This function will return various information about the Anaconda installation on the system by querying the Conda package database.

Warning: THIS IS HIGHLY EXPERIMENTAL AND MAY BREAK WITHOUT FURTHER NOTICE.

Note: This only works with the conda root environment.

Returns: conda_info (*dict*) – Various information from conda (see notes below for further details).

Notes

11/2/16 Page 158 of 209

If the Python interreter is unable to locate and import the conda package, an empty dicionary is returned.

The returned dictionary contains the same infomation that is returned by "conda info" in addition to an overview of the linked modules in the Anaconda installation as well as information about the "conda env" if it is available. Specifically, the returned dictionary has the following keys:

- platform
- conda version
- root_prefix
- default prefix
- envs dirs
- package_cache
- channels
- config_file
- linked_modules
- env_export

For conda < 4.2.0, it also includes the key:

is_foreign_system

For conda \geq 4.2.0, it also includes the keys:

- python_version
- · conda is private
- offline mode

Additionally, the returned dictionary has a key named *status*, which can have either of the following values:

- 'Succeeded' (Everything seems to be OK)
- 'Failed' (Something went wrong a few details are incluede in the key)

If "conda-env" is installed on the system, the *env_export* essentially holds the infomation from "conda env export -n root" as a dictionary. The information provided by this key partially overlaps with the infomation in the *linked modules* and *modules info* keys.

magni.reproducibility._annotation.get_datetime()

Return a dictionary holding the current date and time.

Returns: date_time (*dict*) – The dictionary holding the current date and time.

Notes

The returned dictionary has the following keys:

- today (date and time including timezone offset)
- utcnow (UTC date and time)
- pretty utc (UTC date and time formatted according to current locale)
- status

The status entry informs about the success of the pretty_utc formatting. It has one of the follwing values:

11/2/16 Page 159 of 209

- Succeeded (Everything seems OK)
- Failed (It was not possible to format the time)

magni.reproducibility._annotation. **get_git_revision**(*git_root_dir=None*)

Return a dictionary containing information about the current git revision.

Parameters: git_root_dir (*str*) – The path to the git root directory to get git revision

for (the default is None, which implies that the git revision of the magni

directory is returned).

Returns: git_revision (*dict*) – Information about the current git revision.

Notes

If the git revision extract succeeded, the returned dictionary has the following keys:

status (with value 'Succeeded')

- tag (output of "git describe")
- branch (output of "git describe -all")
- remote (output of "git remote -v")

If the git revision extract failed, the returned dictionary has the following keys:

- status (with value 'Failed')
- returncode (returncode from failing git command)
- output (output from failing git command)

The "git" commands are run in the git root directory.

magni.reproducibility._annotation. **get_file_hashes**(path, blocksize=1073741824)

Return a dictionary with md5 and sha256 checksums of a file.

Parameters: • path (*str*) – The path to the file to checksum.

- blocksize (int) - The chunksize (in bytes) to read from the file one at

a time.

Returns: file_hashes (dict) - The dictionary holding the md5 and sha256

hexdigests of the file.

magni.reproducibility._annotation.get_magni_config()

Return a dictionary holding the current configuration of Magni.

Returns: magni_config (*dict*) – The dictionary holding the current configuration of Magni.

Notes

The returned dictionary has a key for each of the *config* modules in Magni and its subpackages. The value of a given key is a dictionary with the current configuration of the corresponding *config* module. Furthermore, the returned dictionary has a status key, which can have either of the following values:

- Succeeded (The entire configuration was extracted)
- Failed (It was not possible to get information from one or more modules)

magni.reproducibility._annotation.get_magni_info()

11/2/16 Page 160 of 209

Return a string representation of the output of help(magni).

Returns: magni info (*dict*) – Information about magni.

Notes

The returned dictionary has a single key:

help magni (a string representation of help(magni))

magni.reproducibility._annotation.get_platform_info()

Return a dictionary containing information about the system platform.

Returns: platform_info (*dict*) — Various information about the system platform.

See also:

platform()

The Python module used to query information about the system.

Notes

The returned dictionary has the following keys:

- system
- node
- release
- version
- processor
- python
- libc
- linux
- mac os
- win32
- status

The linux/mac os/win32 entries are "empty" if they are not applicable.

If the processor information returned by platform is "empty", a query of *lscpu* is attempted in order to provide the necessary information.

The status entry informs about the success of the queries. It has one of the follwing values:

- 'All OK' (everything seems to be OK)
- 'Used Iscpu in processor query' (*Iscpu* was used)
- 'Processor guery failed' (failed to get processor information)

magni.reproducibility. chase module

Module providing functions that may be used to chase data.

11/2/16 Page 161 of 209

Routine listings

```
get_main_file_name()
```

Function that returns the name of the main file/script.

```
get main file source()
```

Function that returns the source code of the main file/script.

```
get main source()
```

Function that returns the 'local' source code of the main file/script.

```
get stack trace()
```

Function that returns the complete stack trace.

```
magni.reproducibility._chase.get_main_file_name()
```

Return the name of the main file/script which called this function.

This will (if possible) return the name of the main file/script invoked by the Python interpreter i.e., the returned name is the name of the file, in which the bottom call on the stack is defined. Thus, the main file is still returned if the call to this function is buried deep in the code.

Returns: source_file (*str*) – The name of the main file/script which called this function.

Notes

It may not be possible to associate the call to this function with a main file. This is for instance the case if the call is handled by an IPython kernel. If no file is found, the returned string will indicate so.

```
magni.reproducibility. chase. get main file source()
```

Return the source code of the main file/script which called this function.

This will (if possible) return the source code of the main file/script invoked by the Python interpreter i.e., the returned source code is the source code of the file, in which the bottom call on the stack is defined. Thus, the source code is still returned if the call to this function is buried deep in the code.

Returns: source (*str*) – The source code of the main file/script which called this function.

Notes

It may not be possible to associate the call to this function with a main file. This is for instance the case if the call is handled by an IPython kernel. If no file is found, the returned string will indicate so.

```
magni.reproducibility._chase. get_main_source()
```

Return the local source code of the main file which called this function.

This will (if possible) return the part of the source code of the main file/script surrounding (local to) the call to this function. The returned source code is a part of the source code of the file, in which the bottom call on the stack is defined. Thus, the source code is still

11/2/16 Page 162 of 209

returned if the call to this function is buried deep in the code.

Returns: source (*str*) – The local source code of the main file/script.

Notes

It may not be possible to associate the call to this function with a main file. This is for instance the case if the call is handled by an IPython kernel. If no file is found, the returned string will indicate so.

magni.reproducibility._chase. get_stack_trace()

Return the complete stack trace that led to the call to this function.

A (pretty) print version of the stack trace is returned.

Returns: printed_stack (str) – The stack trace that led to the call to this function.

magni.reproducibility.data module

Module providing functions that return data about the system, magni, files, etc.

Functions used to get various annotations and chases data are available from this module.

See also:

magni.reproducibility._annotation

Documentation of annotations

magni.reproducibility._chase

Documentation of chases

magni.reproducibility.io module

Module providing input/output functions to databases containing results from reproducible research.

Routine listings

annotate_database(h5file)

Function for annotating an existing HDF5 database.

chase database(h5file)

Function for chasing an existing HDF5 database.

create database(h5file)

Function for creating a new annotated and chased HDF5 database.

read annotations(h5file)

Function for reading annotations in an HDF5 database.

read chases(h5file)

Function for reading chases in an HDF5 database.

11/2/16 Page 163 of 209 remove annotations(h5file)

Function for removing annotations in an HDF5 database.

remove chases(h5file)

Function for removing chases in an HDF5 database.

write_custom_annotation(h5file, annotation_name, annotation_value, annotations_sub_group=None) Write a custom annotation to an HDF5 database.

See also:

magni.reproducibility. annotation.get conda info

Conda annotation

magni.reproducibility. annotation.get git revision

Git annotation

magni.reproducibility._annotation.get_platform_info

Platform annotation

magni.reproducibility._annotation.get_datetime

Date and time annotation

magni.reproducibility._annotation.get_magni_config

Magni config annotation

magni.reproducibility._annotation.get_magni_info

Magni info annotation

magni.reproducibility._chase.get_main_file_name

Magni main file name chase

magni.reproducibility._chase.get_main_file_source

Magni source code chase

magni.reproducibility._chase.get_main_source

Magni main source code chase

magni.reproducibility._chase.get_stack_trace

Magni stack trace chase

magni.reproducibility.io. annotate_database(h5file)

Annotate an HDF5 database with information about Magni and the platform.

The annotation consists of a group in the root of the h5file having nodes that each provide information about Magni or the platform on which this function is run.

Parameters: h5file (*tables.file*) – The handle to the HDF5 database that should be annotated.

See also:

magni.reproducibility._annotation.get_conda_info()

Conda annotation

magni.reproducibility._annotation.get_git_revision()

Git annotation

magni.reproducibility._annotation.get_platform_info()

Platform annotation

11/2/16 Page 164 of 209

```
magni.reproducibility._annotation.get_datetime()

Date and time annotation

magni.reproducibility._annotation.get_magni_config()

Magni config annotation

magni.reproducibility._annotation.get_magni_info()

Magni info annotation
```

Notes

The annotations of the database includes the following:

- conda info Information about Continuum Anacononda install
- git_revision Git revision and tag of Magni
- platform_info Information about the current platform (system)
- datetime The current date and time
- magni_config Infomation about the current configuration of Magni
- magni_info Information from help(magni)

Examples

Annotate the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import annotate_database
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='a') as h5file:
... annotate_database(h5file)
```

magni.reproducibility.io. **chase** database (*h5file*)

Chase an HDF5 database to track information about stack and source code.

The chase consist of a group in the root of the h5file having nodes that each profide information about the program execution that led to this chase of the database.

Parameters: h5file (*tables.file*) – The handle to the HDF5 database that should be chased.

See also:

```
magni.reproducibility._chase.get_main_file_name()
    Name of main file

magni.reproducibility._chase.get_main_file_source()
    Main file source code

magni.reproducibility._chase.get_main_source()
    Source code around main

magni.reproducibility. chase.get_stack_trace()
```

Complete stack trace

Notes

11/2/16 Page 165 of 209

The chase include the following information:

- main file name Name of the main file/script that called this function
- main file source Full source code of the main file/script
- main source Extract of main file source code that called this function
- stack trace Complete stack trace up until the call to this function

Examples

Chase the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import chase database
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='a') as h5file:
    chase_database(h5file)
```

magni.reproducibility.io. **create database**(*path*, *overwrite=True*)

Create a new HDF database that is annotated and chased.

Α new HDF database is created and it is annotated using magni.reproducibility.io.annotate_database and chased using magni.reproducibility.io.annotate_database. If the overwrite flag is true and existing database at path is overwritten.

- **Parameters:** path (*str*) The path to the HDF file that is to be created.
 - overwrite (bool) The flag that indicates if an existing database should be overwritten.

See also:

magni.reproducibility.io.annotate_database()

Database annotation

magni.reproducibility.io.chase database()

Database chase

Examples

Create a new database named 'new_db.hdf5':

```
>>> from magni.reproducibility.io import create database
>>> create database('new db.hdf5')
```

magni.reproducibility.io. read annotations (h5file)

Read the annotations to an HDF5 database.

Parameters: h5file (tables.file.File) – The handle to the HDF5 database from which

the annotations are read.

annotations (*dict*) – The annotations read from the HDF5 database. Returns: ValueError – If the annotations to the HDF5 database does not conform Raises:

to the Magni annotation standard.

11/2/16 Page 166 of 209

Notes

The returned dict holds a key for each annotation in the database. The value corresponding to a given key is in itself a dict. See *magni.reproducibility.annotate_database* for examples of such annotations.

Examples

Read annotations from the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import read_annotations
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='r') as h5file:
... annotations = read_annotations(h5file)
```

magni.reproducibility.io. **read_chases**(*h5file*)

Read the chases to an HDF5 database.

Parameters: h5file (*tables.file*. *File*) – The handle to the HDF5 database from which

the chases are read.

Returns: chasess (*dict*) – The chases read from the HDF5 database.

Raises: ValueError – If the chases to the HDF5 database does not conform to the

Magni chases standard.

Notes

The returned dict holds a key for each chase in the database. The value corresponding to a given key is a string. See *magni.reproducibility.chase_database* for examples of such chases.

Examples

Read chases from the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import read_chases
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='r') as h5file:
... chases = read_chases(h5file)
```

magni.reproducibility.io. **remove annotations** (*h5file*)

Remove the annotations from an HDF5 database.

Parameters: h5file (*tables.file*) – The handle to the HDF5 database from which the annotations are removed.

Examples

Remove annotations from the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import remove_annotations
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='a') as h5file:
... remove_annotations(h5file)
```

11/2/16 Page 167 of 209

magni.reproducibility.io. **remove chases** (*h5file*)

Remove the chases from an HDF5 database.

Parameters: h5file (*tables.file*.*File*) – The handle to the HDF5 database from which the chases are removed.

Examples

Remove chases from the database named 'db.hdf5':

```
>>> import magni
>>> from magni.reproducibility.io import remove chases
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='a') as h5file:
... remove chases(h5file)
```

magni.reproducibility.io. write custom annotation (h5file, annotation name, annotation value, annotations sub group=None)

Write a custom annotation to an HDF5 database.

The annotation is written to the *h5file* under the *annotation name* such that it holds the annotation value.

- **Parameters:** h5file (tables.file.File) The handle to the HDF5 database to which the annotation is written.
 - **annotation name** (*str*) The name of the annotation to write.
 - annotation_value (a JSON serialisable object) The annotation value to write.
 - annotations_sub_group (str) The group node under "/annotations" to which the custom annotation is written (the default is None which implies that the custom annotation is written directly under "/annotations").

Notes

The annotation value must be a JSON seriablisable object.

Examples

Write a custom annotation to an HDF5 database.

```
>>> import magni
>>> from magni.reproducibility.io import write custom annotation
>>> annotation name = 'custom annotation'
>>> annotation value = 'the value'
>>> with magni.utils.multiprocessing.File('db.hdf5', mode='a') as h5file:
... write_custom_annotation(h5file, annotation_name, annotation_value)
... annotations = magni.reproducibility.io.read_annotations(h5file)
>>> str(annotations['custom_annotation'])
'the value'
```

magni.reproducibility.io. **recursive annotation read** (h5 annotations, out annotations dict)

Recursively read annotations from an annotation group

11/2/16 Page 168 of 209

- **Parameters:** h5_annotations (tables.group.Group) The group to read annotations from.
 - out annotations dict (dict) The dictionary to store the read annotations in.

magni.utils package

Subpackage providing support functionality for the other subpackages.

Routine listings

multiprocessing

Subpackage providing intuitive and extensive multiprocessing functionality.

config

Module providing a robust configger class.

matrices

Module providing matrix emulators.

plotting

Module providing utilities for control of plotting using matplotlib.

validation

Subpackage providing validation capability.

types

Module providing custom data types.

split path(path)

Split a path into folder path, file name, and file extension.

Notes

See _util for documentation of split_path.

Subpackages

magni.utils.multiprocessing package

Subpackage providing intuitive and extensive multiprocessing functionality.

Routine listings

Configger providing configuration options for this subpackage.

File()

Control pytables access to hdf5 files when using multiprocessing.

process(func, namespace={}, args_list=None, kwargs_list=None, maxtasks=None)

Map multiple function calls to multiple processors.

Notes

11/2/16 Page 169 of 209 See _util for documentation of File. See _processing for documentation of process.

Submodules

magni.utils.multiprocessing. config module

Module providing configuration options for the multiprocessing subpackage.

See also:

magni.utils.config.Configger

The Configger class used

Notes

This module instantiates the *Configger* class provided by magni.utils.config. The configuration options are the following:

max_broken_pool_restarts

: int or None

The maximum number of attempts at restarting the process pool upon a BrokenPoolError (the default is 0, which implies that the process pool may not be restarted). If set to None, the process pool may restart indefinitely.

prefer futures

: bool

The indicator of whether or not to prefer the concurrent futures module over the multiprocessing module (the default is False, which implies that the multiprocessing module is used).

re raise exceptions

: bool

A flag indicating if exceptions should be re-raised (the default is False, which implies that the exception are not re-raised). It is useful to set this to True if the processing pool supports proper exception handling as e.g. when using "futures".

silence exceptions

: bool

A flag indicating if exceptions should be silenced (the default is False, which implies that exceptions are raised).

workers

: int

The number of workers to use for multiprocessing (the default is 0, which implies no multiprocessing).

See the notes for the magni.utils.multiprocessing.process function for more details about the *prefer futures* and *max broken pool restarts* configuration parameters.

magni.utils.multiprocessing._processing module

Module providing the process function.

Routine listings

11/2/16 Page 170 of 209

process(func, namespace={}, args_list=None, kwargs_list=None, maxtasks=None) Map multiple function calls to multiple processors.

See also:

magni.utils.multiprocessing.config Configuration options.

magni.utils.multiprocessing. processing. process(func, namespace={}, args list=None, kwargs_list=None, maxtasks=None)

Map multiple function calls to multiple processes.

For each entry in args_list and kwargs_list, a task is formed which is used for a function call of the type func(*args, **kwargs). Each task is executed in a seperate process using the concept of a processing pool.

- **Parameters:** func (function) A function handle to the function which the calls should be mapped to.
 - namespace (dict, optional) A dict whose keys and values should be globally available in func (the default is an empty dict).
 - args_list (list or tuple, optional) A sequence of argument lists for the function calls (the default is None, which implies that no arguments are used in the calls).
 - kwargs_list (list or tuple, optional) A sequence of keyword argument dicts for the function calls (the default is None, which implies that no keyword arguments are used in the calls).
 - maxtasks (int, optional) The maximum number of tasks of a process before it is replaced by a new process (the default is None, which implies that processes are not replaced).

Returns:

results (*list*) – A list with the results from the function calls.

Raises:

BrokenPoolError - If using the concurrent futures module and one or more workers terminate abrubtly with the automatic broken pool restart funtionality disabled.

See also:

magni.utils.multiprocessing.config()

Configuration options.

Notes

If the workers configuration option is equal to 0, map is used. Otherwise, the map functionality of a processing pool is used.

Reasons for using this function over map or standard multiprocessing:

- Simplicity of the code over standard multiprocessing.
- Simplicity in switching between single- and multiprocessing.
- The use of both arguments and keyword arguments in the function calls.
- The reporting of exceptions before termination.
- The possibility of terminating multiprocessing with a single interrupt.

11/2/16 Page 171 of 209 • The option of automatically restarting a broken process pool.

As of Python 3.2, two different, though quite similar, modules exist in the standard library for managing processing pools: multiprocessing and concurrent.futures. According to Python core contributor Jesse Noller, the plan is to eventually make concurrent.futures the only interface to the high level processing pools (futures), whereas multiprocessing is supposed to serve more low level needs for individually handling processes, locks, queues, etc. (see https://bugs.python.org/issue9205#msg132661). As of Python 3.5, both the multiprocessing.Pool and concurrent.futures.ProcessPoolExecutor serve almost the same purpose and provide very similar interfaces. The main differences between the two are:

- The option of using a worker initialiser is only available inmultiprocessing.
- The option of specifing a maximum number of tasks for a worker to execute before being replaced to free up ressources (the maxtasksperchild option) is only available in multiprocessing.
- The option of specifying a context is only available in multiprocessing
- "Reasonable" handling of abrubt worker termination and exceptions is only available in concurrent futures.

Particularly, the "reasonable" handling of a broken process pool may be a strong argument to prefer concurrent.futures over multiprocessing. The matter of handling a broken process pool has been extensively discussed in https://bugs.python.org/issue9205 which led to the fix for concurrent.futures. A similar fix for multiprocessing has been proposed in https://bugs.python.org/issue22393.

Both the multiprocessing and concurrent.futures interfaces are available for use with this function. If the configuration parameter *prefer_futures* is set to True and the concurrent.futures module is available, this is used. Otherwise, the multiprocessing module is used. A Python 2 backport of concurrent.futures is available at https://pythonhosted.org/futures/.

When using concurrent.futures, the *namespace* and *init_args* are ignored since these are not supported by that module. Support for *maxtasks* is emulated in a way such that each process on average does *maxtasks* before it is restarted. The *init_args* functionality may be added later on if an initialiser is added to concurrent.futures - see http://bugs.python.org/issue21423. If the *max_broken_pool_restarts* configuration parameter is set to a value different from 0, the Pool is automatically restarted and the tasks are re-run should a broken pool be encountered. If *max_broken_pool_restarts* is set to 0, a BrokenPoolError is raised should a broken pool be encountered.

When using multiprocessing, the max_broken_pool_restarts is ignored since the BrokenPoolError handling has not yet been implemented for the multiprocessing.Pool - see https://bugs.python.org/issue22393 as well as https://bugs.python.org/issue9205.

Examples

An example of how to use args list and kwargs list:

11/2/16 Page 172 of 209

```
>>> import magni
>>> from magni.utils.multiprocessing._processing import process
>>> def calculate(a, b, op='+'):
... if op == '+':
       return a + b
... elif op == '-':
       return a - b
>>> args_list = [[5, 7], [9, 3]]
>>> kwargs_list = [{'op': '+'}, {'op': '-'}]
>>> process(calculate, args_list=args_list, kwargs_list=kwargs_list)
[12, 6]
```

or the same example preferring concurrent.futures over multiprocessing:

```
>>> magni.utils.multiprocessing.config['prefer futures'] = True
>>> process(calculate, args list=args list, kwargs list=kwargs list)
[12, 6]
```

magni.utils.multiprocessing._processing._get_tasks(func, args_list, kwargs_list)

Prepare a list of tasks.

- **Parameters:** func (function) A function handle to the function which the calls should be mapped to.
 - args_list (list or tuple, optional) A sequence of argument lists for the function calls (the default is None, which implies that no arguments are used in the calls).
 - kwargs list (list or tuple, optional) A sequence of keyword argument dicts for the function calls (the default is None, which implies that no keyword arguments are used in the calls).

Returns:

tasks (*list*) – The list of tasks.

magni.utils.multiprocessing._processing._map_using_futures(func, tasks, init_args, maxtasks, max_broken_pool_restarts)

Map a set of *tasks* to *func* and run them in parallel using a futures.

If max broken pool restarts is different from 0, the tasks must be an explicit collection, e.g. a list or tuple, for the restart to work. If an exception occurs in one of the function calls, the process pool terminates ASAP and re-raises the first exception that occurred. All exceptions that may have occurred in the workers are available as the last element in the exceptions args.

- **Parameters:** func (function) A function handle to the function which the calls should be mapped to.
 - tasks (iterable) The list of tasks to use as arguments in the function
 - init args (tuple) The (func, namespace) tuple for the process init initialisation function.
 - maxtasks (int) The maximum number of tasks of a process before it is replaced by a new process. If set to None, the process is never replaced.
 - max broken pool restarts (int or None) The maximum number of attempts at restarting the process pool upon a BrokenPoolError. If set to None, the process pool may restart indefinitely.

Returns:

results (*list*) – The list of results from the map operation.

11/2/16 Page 173 of 209

Notes

The namespace and init args are ignored since these are not supported by concurrent.futures. Support for maxtasks is emulated in a way such that each process on average does maxtasks before it is restarted. The init args functionality may be added if an initialiser is added to concurrent.futures - see http://bugs.python.org/issue21423.

magni.utils.multiprocessing. processing. map using mppool (func, tasks, init args, maxtasks, max broken pool restarts)

Map a set of *tasks* to *func* and run them in parallel via multiprocessing

- **Parameters:** func (function) A function handle to the function which the calls should be mapped to.
 - tasks (iterable) The list of tasks to use as arguments in the function
 - init args (tuple) The (func, namespace) tuple for the process init initialisation function.
 - maxtasks (int) The maximum number of tasks of a process before it is replaced by a new process. If set to None, the process is never replaced.
 - max broken pool restarts (int or None) The maximum number of attempts at restarting the process pool upon a BrokenPoolError. If set to None, the process pool may restart indefinitely.

Returns:

results (*list*) – The list of results from the map operation.

Notes

The max broken pool restarts is ignored since the BrokenPoolError handling has not been implemented multiprocessing.Pool in the see yet https://bugs.python.org/issue22393 and https://bugs.python.org/issue9205.

magni.utils.multiprocessing._processing._process_init(func, namespace)

Initialise the process by making global variables available to it.

- **Parameters:** func (function) A function handle to the function which the calls should be mapped to.
 - namespace (dict) A dict whose keys and values should be globally available in func.

magni.utils.multiprocessing._processing._process_worker(fak_tuple)

Unpack and map a task to the function.

- Parameters: fak_tuple (tuple) A tuple (func, args, kwargs) containing the parameters listed below.
 - func (function) A function handle to the function which the calls should be mapped to.
 - args (list or tuple) The sequence of arguments that should be unpacked and passed.
 - kwargs (list or tuple) The sequence of keyword arguments that should be unpacked and passed.

Notes

11/2/16 Page 174 of 209

If an exception is raised in func, the stacktrace of that exception is printed since the exception is otherwise silenced until every task has been executed when using multiple workers.

Also, a workaround has been implemented to allow KeyboardInterrupts to interrupt the current tasks and all remaining tasks. This is done by setting a global variable, when catching a KeyboardInterrupt, which is checked for every call.

magni.utils.multiprocessing. util module

Module providing the public class of the magni.utils.multiprocessing subpackage.

class magni.utils.multiprocessing._util. File (*args, **kwargs)

Bases: object

Control pytables access to hdf5 files when using multiprocessing.

File retains the interface of tables.open file and should only be used in 'with' statements (see Examples).

- **Parameters:** args (tuple) The arguments that are passed to 'tables.open file'.
 - kwargs (dict) The keyword arguments that are passed to 'tables.open file'.

See also:

tables.open_file

The wrapped function.

Notes

Internally the module uses a global lock which is shared amongst all files. This solution is simple and does not entail significant overhead. However, the wait time introduced when using multiple files at the same time can be significant.

Examples

The class is used in the following way:

```
>>> from magni.utils.multiprocessing._util import File
>>> with File('database.hdf5', 'a') as f:
    pass # execute something involving the opened file
enter__()
    Acquire the global lock before opening and returning the file.
    Returns: file (tables.File) – The file specified in the call to
                  init .
 exit (type, value, traceback)
```

Release the global lock after closing the file.

11/2/16 Page 175 of 209

- **Parameters:** type (*type*) The type of the exception raised, if any.
 - **value** (*Exception*) The exception rasied, if any.
 - traceback (traceback) The traceback of the exception raised, if

magni.utils.validation package

Subpackage providing validation capability.

The intention is to validate all public functions of the package such that erroneous arguments in calls are reported in an informative fashion rather than causing arbitrary exceptions or unexpected results. To avoid performance impairments, the validation can be disabled globally.

Routine listings

types

Module providing abstract superclasses for validation.

decorate validation(func)

Decorate a validation function (see Notes).

disable_validation()

Disable validation globally (see Notes).

enable validate once()

Enable validating inputs only once (see Notes).

validate_generic(name, type, value_in=None, len_=None, keys_in=None,

has keys=None, superclass=None, ignore none=False, var=None) Validate nonnumeric objects.

validate levels(name, levels)

Validate containers and mappings as well as contained objects.

validate numeric(name, type, range ='[-inf;inf]', shape=(), precision=None,

ignore_none=False, var=None) Validate numeric objects.

validate once(func)

Decorate a function to allow for a one-time input validation (see Notes).

Notes

To be able to disable validation (and to ensure consistency), every public function or method should define a nested validation function with the name 'validate input' which takes no arguments. This function should be decorated by decorate_validation, be placed in the beginning of the parent function or method, and be called as the first thing after its definition.

Functions in magni may be decorated by validate once. If the validate once functionality is enabled, these functions only validate their input arguments on the first call to the function.

Examples

If, for example, the following function is defined:

11/2/16 Page 176 of 209

```
>>> def greet(person, greeting):
... print('{}, {} {} {} {} .'.format(greeting, person['title'], person['name']))
```

This function expects its argument, 'person' to be a dictionary with keys 'title' and 'name' and its argument, 'greeting' to be a string. If, for example, a list is passed as the first argument, a TypeError is raised with the description 'list indices must be integers, not str'. While obviously correct, this message is not excessively informative to the user of the function. Instead, this module can be used to redefine the function as follows:

```
>>> from magni.utils.validation import decorate_validation, validate_generic
>>> def greet(person, greeting):
... @decorate_validation
... def validate_input():
... validate_generic('person', 'mapping', has_keys=('title', 'name'))
... validate_generic('greeting', 'string')
... validate_input()
... validate_input()
... print('{}, {} {}.'.format(greeting, person['title'], person['name']))
```

If, again, a list is passed as the first argument, a TypeError with the description "The value(s) of >>type(person)<<, <type 'list'>, must be in ('mapping',)." is raised. Now, the user of the function can easily identify the mistake and correct the call to read:

```
>>> greet({'title': 'Mr.', 'name': 'Anderson'}, 'You look surprised to see me')
You look surprised to see me, Mr. Anderson.
```

Submodules

magni.utils.validation._generic module

Module providing the validate_generic function.

Routine listings

validate_generic(name, type, value_in=None, len_=None, keys_in=None, has_keys=None, superclass=None, ignore_none=False, var=None) Validate non-numeric objects.

```
magni.utils.validation._generic. validate_generic(name, type_, value_in=None, len_=None, keys_in=None, has_keys=None, superclass=None, ignore_none=False, var=None) Validate non-numeric objects.
```

The present function is meant to validate the type or class of an object. Furthermore, if the object may only take on a limited number of values, the object can be validated against this list. In the case of collections (for example lists and tuples) and mappings (for example dictionaries), a specific length can be required. Furthermore, in the case of mappings, the keys can be validated by requiring and/or only allowing certain keys.

If the present function is called with *name* set to None, an iterable with the value 'generic' followed by the remaining arguments passed in the call is returned. This is useful in combination with the validation function magni.utils.validation.validate_levels.

11/2/16 Page 177 of 209

- **Parameters:** name (*None*) The name of the variable to be validated.
 - type (None) One or more references to groups of types, specific types, and/or specific classes.
 - value_in (set-like) The list of values accepted. (the default is None, which implies that all values are accepted)
 - **len**_ (*int*) The length required. (the default is None, which implies that all lengths are accepted)
 - keys_in (set-like) The list of accepted keys. (the default is None, which implies that all keys are accepted)
 - has_keys (set-like) The list of required keys. (the default is None, which implies that no keys are required)
 - **superclass** (*class*) The required superclass. (the default is None, which implies that no superclass is required)
 - **ignore none** (bool) A flag indicating if the variable is allowed to be none. (the default is False)
 - var (None) The value of the variable to be validated.

See also:

magni.utils.validation.validate_levels()

Validate contained objects.

magni.utils.validation.validate_numeric()

Validate numeric objects.

Notes

name must refer to a variable in the parent scope of the function or method decorated by magni.utils.validation.decorate_validation which is closest to the top of the call stack. If name is a string then there must be a variable of that name in that scope. If name is a set-like object then there must be a variable having the first value in that set-like object as name. The remaining values are used as keys/indices on the variable to obtain the variable to be validated. For example, the name ('name', 0, 'key') refers to the variable "name[0] ['key']".

type is either a single value treated as a list with one value or a set-like object containing at least one value. Each value is either a specific type or class, or it refers to one or more types by having one of the string values 'string', 'explicit collection', 'impliciet collection', 'collection', 'mapping', 'function', 'class'.

- 'string' tests if the variable is a str.
- 'explicit collection' tests if the variable is a list or tuple.
- 'implicit collection' tests if the variable is iterable.
- 'collection' is a combination of the two above.
- 'mapping' tests if the variable is a dict.
- 'function' tests if the variable is a function.
- 'class' tests if the variable is a type.

var can be used to pass the value of the variable to be validated. This is useful either when the variable cannot be looked up by *name* (for example, if the variable is a property of the argument of a function) or to remove the overhead of looking up the value.

11/2/16 Page 178 of 209

Examples

Every public function and method of the present package (with the exception of the functions of this subpackage itself) validates every argument and keyword argument using the functionality of this subpackage. Thus, for examples of how to use the present function, browse through the code.

magni.utils.validation._generic._check_inheritance(name, var, superclass)

Check the superclasses of a variable.

- Parameters: name (None) The name of the variable to be validated.
 - var (None) The value of the variable to be validated.
 - **superclass** (*class*) The required superclass.

magni.utils.validation. generic. **check keys**(*name*, *var*, *keys in*, *has keys*)

Check the keys of a variable.

- Parameters: name (None) The name of the variable to be validated.
 - var (None) The value of the variable to be validated.
 - **keys_in** (*set-like*) The allowed keys.
 - has keys (set-like) The required keys.

magni.utils.validation. generic. **check len**(name, var, len)

Check the length of a variable.

- Parameters: name (None) The name of the variable to be validated.
 - var (None) The value of the variable to be validated.
 - **len** (*int*) The required length.

magni.utils.validation. generic. **check type**(name, var, types)

Check the type of a variable.

- Parameters: name (None) The name of the variable to be validated.
 - **var** (*None*) The value of the variable to be validated.
 - **types** (*set-like*) The allowed types.

magni.utils.validation._generic._check_value(name, var, value_in)

Check the value of a variable.

- Parameters: name (None) The name of the variable to be validated.
 - var (None) The value of the variable to be validated.
 - value in (set-like) The allowed values.

magni.utils.validation. levels module

Module providing the validate_levels function.

Routine listings

11/2/16 Page 179 of 209 validate levels(name, levels)

Validate containers and mappings as well as contained objects.

magni.utils.validation._levels. **validate_levels**(name, levels)

Validate containers and mappings as well as contained objects

The present function is meant to valdiate the 'levels' of a variable. That is, the value of the variable itself, the values of the second level (in case the value is a list, tuple, or dict), the values of the third level, and so on.

- Parameters: name (None) The name of the variable to be validated.
 - **levels** (*list or tuple*) The list of levels.

See also:

magni.utils.validation.validate_generic()

Validate non-numeric objects.

magni.utils.validation.validate_numeric()

Validate numeric objects.

Notes

name must refer to a variable in the parent scope of the function or method decorated by magni.utils.validation.decorate_validation which is closest to the top of the call stack. If name is a string then there must be a variable of that name in that scope. If name is a set-like object then there must be a variable having the first value in that set-like object as name. The remaining values are used as keys/indices on the variable to obtain the variable to be validated. For example, the *name* ('name', 0, 'key') refers to the variable "name[0] ['key']".

levels is a list containing the levels. The value of the variable is validated against the first level. In case the value is a list, tuple, or dict, the values contained in this are validated against the second level and so on. Each level is itself a list with the first value being either 'generic' or 'numeric' followed by the arguments that should be passed to the respective function (with the exception of *name* which is automatically prepended by the present function).

Examples

Every public function and method of the present package (with the exception of the functions of this subpackage itself) validates every argument and keyword argument using the functionality of this subpackage. Thus, for examples of how to use the present function, browse through the code.

magni.utils.validation. levels. **validate level** (name, var, levels, index=0) Validate a level.

11/2/16 Page 180 of 209

- **Parameters:** name (*None*) The name of the variable.
 - var (None) The value of the variable.
 - **levels** (set-like) The levels.
 - index (int) The index of the current level. (the default is 0)

magni.utils.validation. numeric module

Module providing the validate_numeric function.

Routine listings

```
validate_numeric(name, type, range_='[-inf;inf]', shape=(), precision=None,
    ignore none=False, var=None) Validate numeric objects.
```

```
magni.utils.validation._numeric.validate_numeric(name, type_, range_='[-inf;inf]', shape=(),
precision=None, ignore none=False, var=None)
```

Validate numeric objects.

The present function is meant to valdiate the type or class of an object. Furthermore, if the object may only take on a connected range of values, the object can be validated against this range. Also, the shape of the object can be validated. Finally, the precision used to represent the object can be validated.

If the present function is called with name set to None, an iterable with the value 'numeric' followed by the remaining arguments passed in the call is returned. This is useful in combination with the validation function magni.utils.validation.validate_levels.

- **Parameters:** name (*None*) The name of the variable to be validated.
 - **type** (*None*) One or more references to groups of types.
 - range_ (None) The range of accepted values. (the default is '[inf;inf]', which implies that all values are accepted)
 - shape (list or tuple) The accepted shape. (the default is (), which implies that only scalar values are accepted)
 - **precision** (*None*) One or more precisions.
 - ignore none (bool) A flag indicating if the variable is allowed to be none. (the default is False)
 - var (None) The value of the variable to be validated.

See also:

magni.utils.validation.validate_generic()

Validate non-numeric objects.

magni.utils.validation.validate_levels()

Validate contained objects.

Notes

name must refer to a variable in the parent scope of the function or method decorated by magni.utils.validation.decorate_validation which is closest to the top of the call stack. If name is a string then there must be a variable of that name in that scope. If name is a set-like object then there must be a variable having the first value in that set-like object as name.

11/2/16 Page 181 of 209 The remaining values are used as keys/indices on the variable to obtain the variable to be validated. For example, the *name* ('name', 0, 'key') refers to the variable "name[0] ['key']".

type_ is either a single value treated as a list with one value or a set-like object containing at least one value. Each value refers to a number of data types depending if the string value is 'boolean', 'integer', 'floating', or 'complex'.

- 'boolean' tests if the variable is a bool or has the data typenumpy.bool8.
- 'integer' tests if the variable is an int or has the data typenumpy.int8, numpy.int16, numpy.int32, Or numpy.int64.
- 'floating' tests if the variable is a float or has the data typenumpy.float16, numpy.float32, numpy.float64, Or numpy.float128.
- 'complex' tests if the variable is a complex or has the data type*numpy.complex32*, numpy.complex64, Or numpy.complex128.

Because *bool* is a subclass of *int*, a *bool* will pass validation as an 'integer'. This, however, is not the case for numpy.bool8.

range_ is either a list with two strings or a single string. In the latter case, the default value of the argument is used as the second string. The first value represents the accepted range of real values whereas the second value represents the accepted range of imaginary values. Each string consists of the following parts:

- One of the following delimiters: '[', '(', ']'.
- A numeric value (or '-inf').
- · A semi-colon.
- A numeric value (or 'inf').
- One of the following delimiters: ']', ')', '['.

shape is either None meaning that any shape is accepted or a list of integers. In the latter case, the integer -1 may be used to indicate that the given axis may have any length.

precision is either an integer treated as a list with one value or a list or tuple containing at least one integer. Each value refers to an accepted number of bits used to store each value of the variable.

var can be used to pass the value of the variable to be validated. This is useful either when the variable cannot be looked up by *name* (for example, if the variable is a property of the argument of a function) or to remove the overhead of looking up the value.

Examples

Every public function and method of the present package (with the exception of the functions of this subpackage itself) validates every argument and keyword argument using the functionality of this subpackage. Thus, for examples of how to use the present function, browse through the code.

magni.utils.validation._numeric._check_precision(name, dtype, types_, precision)

Check the precision of a variable.

11/2/16 Page 182 of 209

- **Parameters:** name (*None*) The name of the variable to be validated.
 - dtype (type) The data type of the variable to be validated.
 - **types_** (*set-like*) The list of accepted types.
 - **precision** (*None*) The accepted precision(s).

magni.utils.validation._numeric._check_range(name, bounds, range_)

Check the range of a variable.

- **Parameters:** name (*None*) The name of the variable to be validated.
 - bounds (list or tuple) The bounds of the variable to be validated.
 - range_ (*None*) The accepted range(s).

magni.utils.validation._numeric._check_shape(name, dshape, shape)

Check the shape of a variable.

- **Parameters:** name (*None*) The name of the variable to be validated.
 - dshape (list or tuple) The shape of the variable to be validated.
 - **shape** (*list or tuple*) The accepted shape.

magni.utils.validation._numeric._check_type(name, dtype, types_)

Check the type of a variable.

- **Parameters:** name (*None*) The name of the variable to be validated.
 - dtype (type) The data type of the variable to be validated.
 - **types**_ (*set-like*) The accepted types.

magni.utils.validation. numeric. **examine var**(*name*, *var*)

Examine a variable.

The present function examines the data type, the bounds, and the shape of a variable. The variable can be of a built-in type, a numpy type, or a subclass of the magni.utils.validation.types.MatrixBase Class.

- Parameters: name (None) The name of the variable to be
 - var (None) The value of the variable to be examined.

Returns:

- **dtype** (*type*) The data type of the examined variable.
- **bounds** (*tuple*) The bounds of the examined variable.
- **dshape** (*tuple*) The shape of the examined variable.

magni.utils.validation. util module

Module providing functionality for disabling validation.

The disabling functionality is provided through a decorator for validation functions and a function for disabling validation. Furthermore, the present module provides functionality which is internal to magnitudis.validation.

Routine listings

11/2/16 Page 183 of 209 decorate_validation(func)

Decorate a validation function to allow disabling of validation checks.

disable validation()

Disable validation checks in magni.

enable validate once()

Enable validating inputs only once in certain functions in magni.

get var(name)

Retrieve the value of a variable through call stack inspection.

report(type, description, format_args=(), var_name=None, var_value=None, expr='{}', prepend='') Raise an exception.

validate_once(func)

Decorate a function to allow for a one-time input validation only.

magni.utils.validation._util. **decorate_validation**(func)

Decorate a validation function to allow disabling of validation checks.

Parameters: func (function) – The validation function to be

decorated.

Returns: func (*function*) – The decorated validation function.

See also:

disable validation()

Disabling of validation checks.

Notes

This decorater wraps the validation function in another function which checks if validation has been disabled. If validation has been disabled, the validation function is not called. Otherwise, the validation function is called.

Examples

See disable validation for an example.

magni.utils.validation._util. disable_validation()

Disable validation checks in magni.

See also:

decorate_validation()

Decoration of validation functions.

Notes

This function merely sets a global flag and relies on decorate_validation to perform the actual disabling.

11/2/16 Page 184 of 209

Examples

An example of a function which accepts only an integer as argument:

```
>>> import magni
>>> def test(arg):
... @magni.utils.validation.decorate_validation
... def validate_input():
... magni.utils.validation.validate_numeric('arg', 'integer')
... validate_input()
```

If the function is called with anything but an integer, it fails:

```
>>> try:
... test('string')
... except BaseException:
... print('An exception occured')
... else:
... print('No exception occured')
An exception occured
```

However, if validation is disabled, the same call does not fail:

```
>>> from magni.utils.validation import disable_validation
>>> disable_validation()
>>> try:
... test('string')
... except BaseException:
... print('An exception occured')
... else:
... print('No exception occured')
No exception occured
```

Although strongly discouraged, validation can be re-enabled after being disabled in the following way:

```
>>> import magni.utils.validation._util
>>> magni.utils.validation._util._disabled = False
```

As this interface is not public it may be subject to changes in future releases without further notice.

```
magni.utils.validation._util.enable_validate_once()
```

Enable validating inputs only once in certain functions in magni.

```
See also:

validate_once()

Decoration of functions.
```

Notes

This function merely sets a global flag and relies on validate_once to implement the actual "validate once" functionality.

11/2/16 Page 185 of 209

Examples

An example of a function which accepts only an integer as argument:

```
>>> import magni
>>> @magni.utils.validation.validate_once
... def test(arg):
... @magni.utils.validation.decorate_validation
... def validate_input():
... magni.utils.validation.validate_numeric('arg', 'integer')
... validate_input()
```

If "validate once" is enabled, the function validation is only called on its first run:

```
>>> magni.utils.validation.enable_validate_once()
>>> try:
... test('string')
... except BaseException:
... print('An exception occured')
... else:
... print('No exception occured')
An exception occured

>>> try:
... test('string')
... except BaseException:
... print('An exception occured')
... else:
... print('No exception occured')
No exception occured
```

magni.utils.validation._util. get_var(name)

Retrieve the value of a variable through call stack inspection.

name must refer to a variable in the parent scope of the function or method decorated by magni.utils.validation.decorate_validation which is closest to the top of the call stack. Ifname is a string then there must be a variable of that name in that scope. If name is a set-like object then there must be a variable having the first value in that set-like object as name. The remaining values are used as keys/indices on the variable to obtain the variable to be validated. For example, the name ('name', 0, 'key') refers to the variable "name[0] ['key']".

```
Parameters: name (None) – The name of the variable to be retrieved.
```

Returns: var (*None*) – The value of the retrieved variable.

Notes

The present function searches the call stack from top to bottom until it finds a function named 'wrapper' defined in this file. That is, until it finds a decorated validation function. The present function then looks up the variable indicated by *name* in the parent scope of that decorated validation function.

```
magni.utils.validation._util. report(type_, description, format_args=(), var_name=None, var_value=None, expr='{}', prepend='')

Raise an exception.
```

The type of the exception is given by type, and the message can take on many forms.

11/2/16 Page 186 of 209

This ranges from a single description to a description formatted using a number of arguments, prepended by an expression using a variable name followed by the variable value, prepended by another description.

- **Parameters:** type (type) The exception type.
 - **description** (*str*) The core description.
 - format_args (list or tuple) The arguments which description is formatted with. (the default is (), which implies no arguments)
 - var name (*None*) The name of the variable which the description concerns. (the default is None, which implies that no variable name and value is prepended to the description)
 - var value (None) The value of the variable which the description concerns. (the default is None, which implies that the value is looked up from the *var name* in the call stack if *var name* is not None)
 - expr (str) The expression to evaluate with the variable name. (the default is '{}', which implies that the variable value is used directly)
 - **prepend** (*str*) The text to prepend to the message generated by the previous arguments. (the default is ")

Notes

name must refer to a variable in the parent scope of the function or method decorated by magni.utils.validation.decorate_validation which is closest to the top of the call stack. If name is a string then there must be a variable of that name in that scope. If name is a set-like object then there must be a variable having the first value in that set-like object as name. The remaining values are used as keys/indices on the variable to obtain the variable to be validated. For example, the name ('name', 0, 'key') refers to the variable "name[0] ['key']".

magni.utils.validation._util. **validate_once**(func)

Decorate a function to allow for a one-time input validation only.

Parameters: func (*function*) – The validation function to be

decorated.

Returns: **func** (*function*) – The decorated validation function.

See also:

enable_validate_once()

Enabling allow validate once functionality.

Notes

This decorater wraps any function. It checks if the "validate once" has been enabled. If it has been enabled, the validation function is only called on the first run. Otherwise, the validation function is always called if not otherwise disabled.

magni.utils.validation.types module

Module providing abstract superclasses for validation.

Routine listings

11/2/16 Page 187 of 209 MatrixBase(object)

Abstract base class of custom matrix classes.

MMSEInputChannel(object)

Abstract base class of a Minimum Mean Squaread Error (MMSE) input channel.

StopCriterion(object)

Abstract base class of a stop criterion.

ThresholdOperator(object)

Abstract base class of a threshold operator.

class magni.utils.validation.types. **MatrixBase**(dtype, bounds, shape)

Bases: object

Abstract base class of custom matrix classes.

The magnitutils.validation.validate_numeric function accepts built-in numeric types, numpy built-in numeric types, and subclasses of the present class. In order to perform validation checks, the validation function needs to know the data type, the bounds, and the shape of the variable. Thus, subclasses must call the init function of the present class with these arguments.

- **Parameters: dtype** (*type*) The data type of the values of the instance.
 - bounds (list or tuple) The bounds of the values of the instance.
 - **shape** (*list or tuple*) The shape of the instance.

bounds

list or tuple – The bounds of the values of the instance.

dtype

type – The data type of the values of the instance.

shape

list or tuple – The shape of the instance.

Notes

dtype is either a built-in numeric type or a numpy built-in numeric type.

If the matrix has complex values, bounds is a list with two values; The bounds on the real values and the bounds on the imaginary values. If, on the other hand, the matrix has real values, bounds has one value; The bounds on the real values. Each such bounds value is a list with two real, numeric values; The lower bound (that is, the minimum value) and the upper bound (that is, the maximum value).

bounds

dtype

shape

class magni.utils.validation.types. MMSEInputChannel(var)

Bases: object

11/2/16 Page 188 of 209 Abstract base class of a Minimum Mean Squaread Error (MMSE) input channel.

The magni.cs.reconstruction algorithms may make use of input channels to define the prior knowledge on the sought solution. In order for the validation of such an input channel to work, it must be based on this class.

Parameters: var (*dict*) – The input channel state variables.

compute (var)

Compute the input channel value.

Parameters: var (*dict*) – The variables used in computing of the input channel value.

Returns: • mean (ndarray) – The computed input channel mean.

• **variance** (*ndarray*) – The computed input channel variance.

class magni.utils.validation.types. MMSEOutputChannel(var)

Bases: object

Abstract base class of a Minimum Mean Squaread Error (MMSE) output channel.

The magni.cs.reconstruction algorithms may make use of output channels to define the observation model. In order for the validation of such an output channel to work, it must be based on this class.

Parameters: var (*dict*) – The output channel state variables.

compute (var)

Compute the output channel value.

Parameters: var (*dict*) – The variables used in computing of the output channel value.

Returns: • mean (ndarray) – The computed output channel mean.

• **variance** (*ndarray*) – The computed output channel variance.

class magni.utils.validation.types. **StopCriterion**(*var*)

Bases: object

Abstract base class of a stop criterion.

The magni.cs.reconstruction algorithms are typically iterative algorithms that make use of some stop criterion to determine if it has converged. In order for the validation of such a stop criterion to work, it must be based on this class.

Parameters: var (*dict*) – The stop criterion state variables.

compute (var)

Compute the stop criterion value.

11/2/16 Page 189 of 209

Parameters: var (*dict*) – The variables used in computing of the stop criterion

value.

• **stop** (*bool*) – The indicator of whether or not the stop criterion is satisfied.

• value (float) – The stop criterion value.

class magni.utils.validation.types. **ThresholdOperator** (var)

Bases: object

Abstract base class of a threshold operator.

The magni.cs.reconstruction algorithms may make use of threshold operators for "denoising". In order for the validation of such a threshold operator to work, it must be based on this class.

Parameters: var (*dict*) – The threshold operator state

variables.

compute_deriv_threshold(var)

Compute the entrywise derivative threshold.

Parameters: var (*dict*) – The variables used in computing the derivative

threshold.

Returns: eta deriv (*ndarray*) – The computed entrywise derivative threshold.

compute_threshold(var)

Compute the entrywise threshold.

Parameters: var (*dict*) – The variables used in computing the

threshold.

Returns: eta (*ndarray*) – The computed entrywise threshold.

update_threshold_level(var)

Update the threshold level state.

Parameters: var (*dict*) – The variables used in computing the threshold level

update.

Submodules

magni.utils._util module

Module providing the public function of the magni.utils subpackage.

magni.utils._util. **split_path**(*path*)

Split a path into folder path, file name, and file extension.

The returned folder path ends with a folder separation character while the returned file extension starts with an extension separation character. The function is independent of the operating system and thus of the use of folder separation character and extension separation character.

11/2/16 Page 190 of 209

Parameters: path (str) – The path of the file either absolute or relative to the current working directory.

Returns:

- path (str) The path of the containing folder of the input path.
- **name** (*str*) The name of the object which the input path points to.
- ext (str) The extension of the object which the input path points to (if any).

Examples

Concatenate a dummy path and split it using the present function:

```
>>> import os
>>> from magni.utils._util import split_path
>>> path = 'folder' + os.sep + 'file' + os.path.extsep + 'extension'
>>> parts = split_path(path)
>>> print(tuple((parts[0][-7:-1], parts[1], parts[2][1:])))
('folder', 'file', 'extension')
```

magni.utils.config module

Module providing a robust configger class.

Routine listings

Configger(object)

Provide functionality to access a set of configuration options.

Notes

This module does not itself contain any configuration options and thus has no access to any configuration options unlike the other config modules of magni.

class magni.utils.config. Configger (params, valids)

Bases: object

Provide functionality to access a set of configuration options.

The set of configuration options, their default values, and their validation schemes are specified upon initialisation.

- **Parameters:** params (*dict*) The configuration options and their default values.
 - valids (dict) The validation schemes of the configuration options.

See also:

magni.utils.validation

Validation.

Notes

valids must contain the same keys as params. For each key in 'valids', the first value is the validation function ('generic', 'levels', or 'numeric'), whereas the remaining values

11/2/16 Page 191 of 209 are passed to that validation function.

Examples

Instantiate Configger with the parameter 'key' with default value 'default' which can only assume string values.

```
>>> import magni
>>> from magni.utils.config import Configger
>>> valid = magni.utils.validation.validate_generic(None, 'string')
>>> config = Configger({'key': 'default'}, {'key': valid})
```

The number of parameters can be retrieved as the length:

```
>>> len(config)
1
```

That parameter can be retrieved in a number of ways:

```
>>> config['key']
'default'

>>> for key, value in config.items():
... print('key: {!r}, value: {!r}'.format(key, value))
key: 'key', value: 'default'

>>> for key in config.keys():
... print('key: {!r}'.format(key))
key: 'key'

>>> for value in config.values():
... print('value: {!r}'.format(value))
value: 'default'
```

Likewise, the parameter can be changed in a number of ways:

```
>>> config['key'] = 'value'
>>> config.update({'key': 'value changed by dict'})
>>> config['key']
'value changed by dict'

>>> config.update(key='value changed by keyword')
>>> config['key']
'value changed by keyword'
```

Finally, the parameter can be reset to the default value at any point:

```
>>> config.reset()
>>> config['key']
'default'
```

_funcs = {'numeric': <function validate_numeric at 0x7f7cf46e52f0>, 'levels': <function validate_levels at 0x7f7cf45daea0>, 'generic': <function validate_generic at 0x7f7cf4756bf8>}

11/2/16 Page 192 of 209

```
__getitem__(name)
```

Get the value of a configuration parameter.

Parameters: name (*str*) – The name of the parameter.

Returns: value (*None*) – The value of the

parameter.

__len__()

Get the number of configuration parameters.

Returns: length (int) – The number of

parameters.

__setitem__(name, value)

Set the value of a configuration parameter.

The value is validated according to the validation scheme of that parameter.

Parameters: • name (*str*) – The name of the parameter.

• value (None) — The new value of the parameter.

get(key=None)

Deprecated method.

See also:

Configger.__getitem__()

Replacing method.

Configger.items()

Replacing method.

Configger.keys()

Replacing method.

Configger.values()

Replacing method.

items()

Get the configuration parameters as key, value pairs.

Returns: items (set-like) – The list of parameters.

keys()

Get the configuration parameter keys.

Returns: keys (set-like) – The keys.

reset()

Reset the parameter values to the default values.

11/2/16 Page 193 of 209

```
set(dictionary={}, **kwargs)
```

Deprecated method.

See also:

Configger.__setitem__()

Replacing function.

update(params={}, **kwargs)

Update the value of one or more configuration parameters.

Each value is validated according to the validation scheme of that parameter.

- **Parameters:** params (dict, optional) A dictionary containing the key and values to update. (the default value is an empty dictionary)
 - **kwargs** (*dict*) Keyword arguments being the key and values to update.

values()

Get the configuration parameter values.

Returns: values (set-like) – The values.

magni.utils.matrices module

Module providing matrix emulators.

The matrix emulators of this module are wrappers of fast linear operations giving the fast linear operations the same basic interface as a numpy ndarray. Thereby allowing fast linear operations and numpy ndarrays to be used interchangably in other parts of the package.

Routine listings

Matrix(magni.utils.validation.types.MatrixBase)

Wrap fast linear operations in a matrix emulator.

MatrixCollection(magni.utils.validation.types.MatrixBase)

Wrap multiple matrix emulators in a single matrix emulator.

Separable2DTransform(Matrix)

Wrap a linear 2D separable transform in a matrix emulator.

SumApproximationMatrix(object)

Wrap a sum approximation in a matrix emulator.

norm(A, ord=None)

Compute a norm of a matrix.

See also:

magni.imaging._fastops

Fast linear operations.

magni.imaging._mtx1D

1D transforms matrices for use 2D separable transforms.

11/2/16 Page 194 of 209 class magni.utils.matrices. Matrix (func, conj_trans, args, shape, is_complex=False, is valid=True)

Bases: magni.utils.validation.types.MatrixBase

Wrap fast linear operations in a matrix emulator.

Matrix defines a few attributes and internal methods which ensures that instances have the same basic interface as a numpy ndarray instance without explicitly forming the array. This basic interface allows instances to be multiplied with vectors, to be transposed, to be complex conjuagted, and to assume a shape. Also, instances have an attribute which explicitly forms the matrix as an ndarray.

- **Parameters:** func (function) The fast linear operation applied to the vector when multiplying the matrix with a vector.
 - **conj trans** (*function*) The fast linear operation applied to the vector when multiplying the complex conjuated transposed matrix with a vector.
 - args (list or tuple) The arguments which should be passed to func and *conj trans* in addition to the vector.
 - **shape** (*list or tuple*) The shape of the emulated matrix.
 - is complex (bool) The indicator of whether or not the emulated matrix is defined for the complex numbers (the default is False, which implies that the emulated matrix is defined for the real numbers only).
 - is valid (bool) The indicator of whether or not the fast linear operation corresponds to a valid matrix (se discussion below) (the default is True, which implies that the matrix is considered valid).

See also:

magni.utils.validation.types.MatrixBase

Superclass of the present class.

Notes

The is valid indicator is an implementation detail used to control possibly missing operations in the fast linear operation. For instance, consider the Discrete Fourier Transform (DFT) and its inverse transform. The forward DFT transform is a matrix-vector product involving the corresponding DFT matrix. The inverse transform is a matrix-vetor product involving the complex conjugate transpose DFT matrix. That is, it involves complex conjugation and transposition, both of which are (individually) valid transformations of the DFT matrix. However, when implementing the DFT (and its inverse) using a Fast Fourier Transform (FFT), only the combined (complex conjugate, transpose) operation is available. Thus, when using an magni.util.matrices.Matrix, one can get, e.g., a Matrix.T object corresponding to its transpose. However, the operation of computing a matrix-vector product involving the tranpose DFT matrix is not available and the Matrix.T is, consequently, considered an invalid matrix. Only the combined Matrix. T.conj() or Matrix.conj(). T is considered a valid matrix.

Warning: The transpose operation changed in magni 1.5.0.

11/2/16 Page 195 of 209 For complex matrices, the the T tranpose operation now yields an invalid matrix as described above. Prior to magni 1.5.0, the T would yield the "inverse" which would usually be the complex conjugated transpose for complex matrices.

Examples

For example, the negative identity matrix could be emulated as

```
>>> import numpy as np, magni
>>> from magni.utils.matrices import Matrix
>>> func = lambda vec: -vec
>>> matrix = Matrix(func, func, (), (3, 3))
```

The example matrix will have the desired shape:

```
>>> matrix.shape (3, 3)
```

The example matrix will behave just like an explicit matrix:

```
>>> vec = np.float64([1, 2, 3]).reshape(3, 1)
>>> np.set_printoptions(suppress=True)
>>> matrix.dot(vec)
array([[-1.],
       [-2.],
       [-3.]])
```

If, at some point, an explicit representation of the matrix is required, this can easily be obtained:

Likewise, the transpose of the matrix can be obtained:

```
__array__()
```

Return ndarray representation of the matrix.

matrix state

Return a copy of the internal matrix state.

The internal matrix state consists of:

- func: The forward fast linear operator.
- conj_trans: The backward fast linear operator.
- args: The tuple of extra arguments passed to *func* and *conj trans*.
- is_complex: The indicator of whether or not the matrix is complex.
- is valid: The indicator of whether or not the matrix is valid.

11/2/16 Page 196 of 209

Returns: matrix_state (*dict*) – A copy of the internal matrix state

Α

Explicitly form the matrix.

The fast linear operations implicitly define a matrix which is usually not explicitly formed. However, some functionality might require a more advanced matrix interface than that provided by this class.

Returns: matrix (numpy.ndarray) – The explicit matrix.

Notes

The explicit matrix is formed by multiplying the matrix with the columns of an identity matrix and stacking the resulting vectors as columns in a matrix.

Т

Get the transpose of the matrix.

Returns: matrix (*Matrix*) – The transpose of the matrix.

Notes

The fast linear operation and the fast linear transposed operation of the resulting matrix are the same as those of the current matrix except swapped. The shape is modified accordingly. This returns an *invalid* matrix if the entries are complex numbers as only the complex conjugate transpose is considered valid.

conj()

Get the complex conjugate of the matrix.

Returns: matrix (*Matrix*) – The complex conjugate of the matrix.

Notes

The fast linear operation and the fast linear transposed operation of the resulting matrix are the same as those of the current matrix. This returns an *invalid* matrix if the entries are complex numbers as only the complex conjugate transpose is considered valid.

dot(vec)

Multiply the matrix with a vector.

Parameters: vec (numpy.ndarray) – The vector which the matrix is multiplied with.

Returns: vec (*numpy.ndarray*) – The result of the multiplication.

Notes

11/2/16 Page 197 of 209

This method honors magni.utils.validation.enable_validate_once.

class magni.utils.matrices. MatrixCollection(matrices)

Bases: magni.utils.validation.types.MatrixBase

Wrap multiple matrix emulators in a single matrix emulator.

MatrixCollection defines a few attributes and internal methods which ensures that instances have the same basic interface as a numpy ndarray instance without explicitly forming the ndarray. This basic interface allows instances to be multiplied with vectors, to be transposed, to be complex conjugated, and to assume a shape. Also, instances have an attribute which explicitly forms the matrix.

Parameters: matrices (list or tuple) – The collection of matrices (e.g. ndarrays or Matrix instances).

See also:

magni.utils.validation.types.MatrixBase

Superclass of the present class.

Matrix

Matrix emulator.

Examples

For example, two matrix emulators can be combined into one. That is, the matrix:

And the matrix:

Can be combined into one matrix emulator using the present class:

```
>>> from magni.utils.matrices import MatrixCollection
>>> matrix = MatrixCollection((negate, reverse))
```

The example matrix will have the desired shape:

```
>>> matrix.shape (3, 3)
```

The example matrix will behave just like an explicit matrix:

11/2/16 Page 198 of 209

If, at some point, an explicit representation of the matrix is required, this can easily be obtained:

Likewise, the transpose of the matrix can be obtained:

```
__array__()
```

Return ndarray representation of the matrix.

matrix_state

Return a copy of the internal matrix state.

The internal matrix state consists of:

• matrices: The tuple of matrices in the matrix collection

```
Returns: matrix_state (dict) – A copy of the internal matrix state
```

A

Explicitly form the matrix.

The collection of matrices implicitly defines a matrix which is usually not explicitly formed. However, some functionality might require a more advanced matrix interface than that provided by this class.

Returns: matrix (numpy.ndarray) – The explicit matrix.

Notes

The explicit matrix is formed by multiplying the matrix with the columns of an identity matrix and stacking the resulting vectors as columns in a matrix.

shape

Get the shape of the matrix.

The shape of the product of a number of matrices is the number of rows of the first matrix times the number of columns of the last matrix.

Returns: shape (*tuple*) – The shape of the matrix.

11/2/16 Page 199 of 209

Т

Get the transpose of the matrix.

The transpose of the product of the number of matrices is the product of the transpose of the matrices in reverse order.

Returns: matrix (*MatrixCollection*) – The transpose of the matrix.

conj()

Get the complex conjugate of the matrix.

The complex conjugate of the product of the number of matrices is the product of the complex conjugates of the matrices.

Returns: matrix (*MatrixCollection*) – The complex conjugate of the matrix.

dot(vec)

Multiply the matrix with a vector.

Parameters: vec (numpy.ndarray) – The vector which the matrix is multiplied

with.

Returns: **vec** (*numpy.ndarray*) – The result of the multiplication.

Notes

This method honors magni.utils.validation.enable_validate_once.

class magni.utils.matrices. **Separable2DTransform** (mtx 1, mtx r)

Bases: magni.utils.matrices.Matrix

Wrap a linear 2D separable transform in a matrix emulator.

A linear 2D separable transform is defined by the two matrices mtx / and mtx r in the sense that the Kronecker product kron(mtx I, mtx r) yields the full matrix of the linear 2D separable transform when the linear transform is implemented as a matrix-vector product. See e.g. [1] for the details.

Separable2DTransform defines a few attributes and internal methods which ensures that instances have the same basic interface as a numpy matrix instance without explicitly forming the matrix. This basic interface allows instances to be multiplied with vectors, to be transposed, and to assume a shape. Also, instances have an attribute which explicitly forms the matrix.

- **Parameters:** mtx I (ndarray) The "left" matrix in the defining Kronecker product.
 - mtx r (ndarray) The "right" matrix in the defining Kronecker product.

See also:

magni.utils.matrices.Matrix

Superclass of the present class.

11/2/16 Page 200 of 209

References

[1] A.N. Akansu, R.A. Haddad, and P.R. Haddad, *Multiresolution Signal Decomposition: Transforms, Subbands, and Wavelets*, Academic Press, 2000.

Examples

For example, the transform based on a 2-by-3 matrix and a 5-by-4 matrix

```
>>> import numpy as np
>>> from magni.imaging import vec2mat, mat2vec
>>> from magni.utils.matrices import Separable2DTransform
>>> mtx_l = np.arange(6).reshape(2, 3)
>>> mtx_r = np.arange(40)[::2].reshape(5, 4)
>>> sep_matrix = Separable2DTransform(mtx_l, mtx_r)
```

The full transform matrix shape is:

```
>>> print(tuple(int(s) for s in sep_matrix.shape)) (10, 12)
```

The matrix behaves like an explicit matrix

which, due to the separability of the transform, may also be computed as

```
>>> vec_as_mat = vec2mat(vec, (3, 4))
>>> mat2vec(mtx_l.dot(vec_as_mat.dot(mtx_r.T)))
array([[ 936],
        [ 3528],
        [ 2712],
        [10056],
        [ 4488],
        [16584],
        [ 6264],
        [ 8040],
        [ 29640]])
```

The explicit matrix is given by the kronecker product of the two matrices

```
>>> np.all(sep_matrix.A == np.kron(mtx_l, mtx_r))
True
```

The transpose of the matrix is also easily obtainable

11/2/16 Page 201 of 209

```
>>> sep_matrix_transpose = sep_matrix.T
>>> sep_matrix_transpose.A
array([[ 0, 0, 0, 0, 0, 0, 24, 48, 72, 96],
        [ 0, 0, 0, 0, 0, 6, 30, 54, 78, 102],
        [ 0, 0, 0, 0, 0, 12, 36, 60, 84, 108],
        [ 0, 0, 0, 0, 0, 18, 42, 66, 90, 114],
        [ 0, 8, 16, 24, 32, 0, 32, 64, 96, 128],
        [ 2, 10, 18, 26, 34, 8, 40, 72, 104, 136],
        [ 4, 12, 20, 28, 36, 16, 48, 80, 112, 144],
        [ 6, 14, 22, 30, 38, 24, 56, 88, 120, 152],
        [ 0, 16, 32, 48, 64, 0, 40, 80, 120, 160],
        [ 4, 20, 36, 52, 68, 10, 50, 90, 130, 170],
        [ 8, 24, 40, 56, 72, 20, 60, 100, 140, 180],
        [ 12, 28, 44, 60, 76, 30, 70, 110, 150, 190]])
```

matrix state

Return a copy of the internal matrix state.

The internal matrix state consists of:

- mtx_l: The "left" matrix in the defining Kronecker product.
- mtx_r: The "right" matrix in the defining Kronecker product.

```
Returns: matrix_state (dict) – A copy of the internal matrix state
```

Α

Explicitly form the matrix.

For a linear separable 2D transform, the full explicit transform matrix is given by the Kronecker product of the two matrices that define the separable transform.

Returns: matrix (ndarray) – The explicit matrix.

Т

Get the transpose of the matrix.

The transpose of a Kronekcer product of two matrices is the Kronecker product of the transposes of the two matrices.

```
Returns: matrix (Separable2DTransform) – The transpose of the matrix.
```

conj()

Get the complex conjugate of the matrix.

The complex conjugate of a Kronekcer product of two matrices is the Kronecker product of the complex conjugates of the two matrices.

```
Returns: matrix (Separable2DTransform) – The complex conjugate of the matrix.
```

dot(vec)

Multiply the matrix with a vector.

The matrix-vector product is efficiently computed by mtx_l.dot(V.dot(mtx_r.T)) where V is the matrix from which *vec* was created by stacking its columns.

11/2/16 Page 202 of 209

Parameters: vec (ndarray) – The vector corresponding to the stacked columns

of a matrix which this matrix is multiplied with.

Returns: vec (*ndarray*) – The result of the matrix-vector multiplication.

Notes

This method honors magni.utils.validation.enable_allow_validate_once.

class magni.utils.matrices. SumApproximationMatrix (scaling)

Bases: object

Wrap a sum approximation in a matrix emulator.

This simply emulates computing a scaled sum of the entries of a vector as a matrix vector product.

Parameters: scaling (int or float) – The scaling applied to the sum

approximation.

Examples

For example,

```
>>> import numpy as np
>>> from magni.utils.matrices import SumApproximationMatrix
>>> np.random.seed(seed=6021)
>>> n = 10
>>> vec = np.random.randn(n, 1)
>>> np_printoptions = np.get_printoptions()
>>> np.set_printoptions(precision=5)
>>> np.array([vec.mean()])
array([ 0.43294])
>>> matrix = SumApproximationMatrix(1.0/n)
>>> np.array([matrix.dot(vec)])
array([ 0.43294])
>>> np.set_printoptions(**np_printoptions)
```

dot(vec)

Form the matrix-vector product sum approximation.

Parameters: vec (*ndarray*) – The vector in the matrix-vector product.

Returns: scaled sum (float) – The scaled sum of the input

vector.

Notes

This method honors magni.utils.validation.enable_allow_validate_once.

magni.utils.matrices. **norm** (A, ord=None)

Compute a norm of a matrix.

Efficient norm computations for the magni.utils.matrices matrices.

11/2/16 Page 203 of 209

Parameters: • A (magni.util.matrices.Matrix or magni.utils.matrices.MatrixCollection)

- The matrix which norm is to be computed.

• **ord** (*str*) – The order of the norm to compute.

Returns:

norm (*float*) – The computed matrix norm.

See also:

numpy.linalg.norm()

Numpy's function for computing norms of ndarrays

Notes

This function is a simple wrapper around numpy.linalg.norm. It exploits some properties of the matrix classes in magni.utils.matrices to optimize the norm computation in terms of speed and/or memory usage.

Warning: This function builds the ndarray corresponding to the *A* matrix and passes that on to numpy.linalg.norm if no smarter way of computing the norm has been implemented. Beware of possible "out of memory" issues.

Examples

Compute the Frobenius norm of a Seperable2DTransform

```
>>> import numpy as np
>>> from magni.utils.matrices import Separable2DTransform, norm
>>> mtx_l = np.arange(6).reshape(2, 3)
>>> mtx_r = np.arange(40)[::2].reshape(5, 4)
>>> sep_matrix = Separable2DTransform(mtx_l, mtx_r)
>>> round(float(norm(sep_matrix, 'fro')), 2)
737.16
```

for comparison, the Frobenius computed using np.linalg.norm is

```
>>> round(float(np.linalg.norm(sep_matrix.A)), 2) 737.16
```

magni.utils.plotting module

Module providing utilities for control of plotting using matplotlib.

The module has a number of public attributes which provide settings for colormap cycles, linestyle cycles, and marker cycles that may be used in combination with matplotlib.

Routine listings

```
setup matplotlib(settings={}, cmap=None)
```

Function that set the Magni default matplotlib configuration.

```
colour collections
```

: dict

Collections of colours that may be used in e.g., amatplotlib color cycle / prop cycle.

11/2/16 Page 204 of 209

seq_cmaps

: list

Names of matplotlib.cm colormaps optimized for sequential data.

div cmaps

: list

Names of matplotlib.cm colormaps optimized for diverging data.

linestyles

: list

A subset of linestyles from matplotlib.lines

markers

: list

A subset of markers from matplotlib.markers

Examples

Use the default Magni matplotlib settings.

```
>>> import magni
>>> magni.utils.plotting.setup_matplotlib()
```

Get the normalised 'Blue' colour brew from the psp colour map:

```
>>> magni.utils.plotting.colour_collections['psp']['Blue'] ((0.1255, 0.502, 0.8745),)
```

class magni.utils.plotting. _ColourCollection(brews)

Bases: object

A container for colour maps.

A single colour is stored as an RGB 3-tuple of integers in the interval [0,255]. A set of related colours is termed a colour brew and is stored as a list of colours. A set of related colour brews is termed a colour collection and is stored as a dictionary. The dictionary key identifies the name of the colour collection whereas the value is the list of colour brews.

The default colour collections named "cb*" are colorblind safe, print friendly, and photocopy-able. They have been created using the online ColorBrewer 2.0 tool [1].

Parameters: brews (*dict*) – The dictionary of colour brews from which the colour collection is created.

Notes

Each colour brew is a list (or tuple) of length 3 lists (or tuples) of RGB values.

References

[1] M. Harrower and C. A. Brewer, "ColorBrewer.org: An Online Tool for Selecting Colour Schemes for Maps", *The Cartographic Journal*, vol. 40, pp. 27-37, 2003 (See also: http://colorbrewer2.org/)

11/2/16 Page 205 of 209

```
getitem (name)
```

Return a single colour brew.

The returned colour brew is normalised in the sense of matplotlib normalised rgb values, i.e., colours are 3-tuples of floats in the interval [0, 1].

Parameters: name (*str*) – Name of the colour brew to

return.

Returns: **brew** (tuple) – A colour brew list.

magni.utils.plotting. **setup_matplotlib**(*settings={}*, *cmap=None*)

Adjust the configuration of matplotlib.

Sets the default configuration of matplotlib to optimize for producing high quality plots of the data produced by the functionality provided in the Magni.

- **Parameters:** settings (dict, optional) A dictionary of custom matplotlibrc settings. See examples for details about the structure of the dictionary.
 - cmap (str or tuple, optional) Colormap to be used by matplotlib (the default is None, which implices that the 'coolwarm' colormap is used). supplied it must be a ('colormap name', tuple is matplotlib.colors.Colormap()) tuple.

UserWarning – If the supplied custom settings are invalid. Raises:

Examples

For example, set lines.linewidth=2 and lines.color='r'.

```
>>> from magni.utils.plotting import setup matplotlib
>>> custom_settings = {'lines': {'linewidth': 2, 'color': 'r'}}
>>> setup matplotlib(custom settings)
```

magni.utils.types module

Module providing custom data types.

Routine listings

ClassProperty(property)

Class property.

ReadOnlyDict(collections.OrderedDict)

Read-only ordered dict.

class magni.utils.types. ClassProperty (fget=None, fset=None, fdel=None, doc=None)

Bases: property

Class property.

The present class is a combination of the built-in property type and the built-in classmethod type. That is, it is a property which is invoked on a class rather than an instance but is available to both the class and its instances.

11/2/16 Page 206 of 209

- **Parameters:** fget (function, optional) The getter function of the property. (the default is None, which implies that the property cannot be read)
 - **fset** (function, optional) The setter function of the property. (the default is None, which implies that the property cannot be written)
 - **fdel** (function, optional) The deleter function of the property. (the default is None, which implies that the property cannot be deleted)
 - doc (string, optional) The docstring of the property. (the default is None, which implies that the docstring of the getter function, if any, is used)

See also:

property

Superclass from which all behaviour is inherited or extended.

Examples

The following example illustrates the difference between regular properties and class properties. First, the example class is defined and instantiated:

```
>>> from magni.utils.types import ClassProperty
>>> class Example(object):
   _x_class = 0
  def __init__(self, x):
     self.\_x\_instance = x
... @ClassProperty
  def x class(class ):
     return class_._x_class
... @property
... def x instance(self):
     return self._x_instance
>>> example = Example(1)
```

The regular read-only property works on the instance:

```
>>> print('{/r}'.format(example.x_instance))
>>> print('ls property? {/r}'.format(isinstance(Example.x_instance,
                               property)))
Is property? True
```

The class property, on the other hand, works on the class:

```
>>> print('{!r}'.format(example.x_class))
>>> print('ls property? {!r}'.format(isinstance(Example.x_class,
                               property)))
Is property? False
>>> print('{!r}'.format(Example.x_class))
```

```
get (obj, class =None)
```

The get method of the property.

11/2/16 Page 207 of 209

- **Parameters: obj** (*object*) The instance on which the property is requested.
 - class (type, optional) The class on which the property is requested. (the default is None, which implies that the class is retrieved from obi)

Returns:

value (*None*) – The value of the property.

See also:

```
property.__get__()
```

The method being extended.

```
set__(obj, value)
```

The set method of the property.

- **Parameters: obj** (*object*) The instance on which the property is requested.
 - value (None) The value which should be assigned to the property.

See also:

```
property.__set__()
```

The method being extended.

```
delete (obj)
```

The delete method of the property.

Parameters: obj (*object*) – The instance on which the property is

requested.

See also:

```
property.__delete__()
```

The method being extended.

class magni.utils.types. **ReadOnlyDict**(*args, **kwargs)

Bases: collections.OrderedDict

Read-only ordered dict.

The present ordered dict subclass has its non-read-only methods disabled.

See also:

collections.OrderedDict

Superclass from which all read-only behaviour is inherited.

Examples

This ordered dict subclass exposes the same interface as the OrderedDict class when

11/2/16 Page 208 of 209 not using methods that alter the dict.

```
>>> from magni.utils.types import ReadOnlyDict
>>> d = ReadOnlyDict(key='value')
>>> for item in d.items():
... print('{/r}'.format(item))
('key', 'value')
```

However, any attempt to assign another value to the property raises an exception:

```
>>> try:
... d['key'] = 'new value'
... except Exception:
... print('An exception occured.')
An exception occured.
```

```
__delitem__(name)
```

Prevent deletion of items.

```
Parameters: name (str) – The name of the item to be deleted.
```

```
__getattribute__(name)
```

Return the requested attribute unless it is a non-read-only method.

The purpose of this overwrite is to disable access to the following non-read-only dict methods: clear, pop, popitem, setdefault, and update. The first two methods are disabled otherwise.

```
Parameters: name (str) – The name of the requested
```

attribute.

Returns: attribute (*None*) – The requested attribute.

Notes

<u>__getattribute__</u> is implicitly called, when the attribute of an object is accessed as object.attribute.

```
__setitem__(name, value)
```

Prevent overwrite of items.

```
    name (str) – The name of the item to be overwritten.
    value (None) – The value to be written.
```

Magni 1.6.0 documentation »

© Copyright 2014-2016, Magni developers. Last updated on Nov 02, 2016. Created using Sphinx 1.4.1.

11/2/16 Page 209 of 209